



Equipment Type: 300 mm (~12-inch)
Automated Abrasive Cut-off Machine

Model: **MEGA-T300A**

Electrical Requirements: 208Y/380 Volts (3-phase)

Frequency: 50/60 Hz

Motor Horsepower: 3.1 KW @ 208Y 3-phase (4.16 HP)
5.7 KW @ 380V 3-phase (7.64 HP)

Manual Revision Date: January 31, 2023

This instruction manual is
provided with each piece of
delivered equipment.

WARRANTY

Terms and Conditions applying to all PACE Technologies Products

1. LIMITED WARRANTY AND DISCLAIMER:

PACE Technologies Products are warranted for two years from the purchase date to be free from defects in material and workmanship under correct use, normal operating conditions, and proper application.

PACE Technologies obligation under this warranty shall be limited to the repair or exchange, at PACE Technologies option, of any PACE Technologies Product or part which proves to be defective as provided herein. PACE Technologies reserves the right to either inspect the product at Buyer's location or require it to be returned to the factory for inspection. Buyer is responsible for freight to and from factory on all warranty claims. The above warranty does not extend to goods damaged or subjected to accident, abuse or misuse after release from PACE Technologies warehouse, nor goods altered or repaired by anyone other than specifically authorized PACE Technologies representatives. PACE Technologies shall not in any way be responsible for the consequences of any alteration, modification or misuse unless previously approved in writing by an officer of PACE Technologies. **Note: Corrosion is considered a maintenance issue and not a warranty issue.**

PACE TECHNOLOGIES MAKES NO EXPRESS WARRANTIES OTHER THAN THOSE WHICH ARE SPECIFICALLY DESCRIBED HEREIN. Any description of the goods sold hereunder, including any reference to Buyer's specifications and any description in catalogs, circulars and other written material published by PACE Technologies, is the sole purpose of identifying such goods and shall not create an express warranty that the goods shall conform to such description.

THIS WARRANTY IS EXPRESSLY IN LIEU OF ALL OTHER WARRANTIES, EXPRESSED OR IMPLIED. THERE ARE NO IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR PARTICULAR PURPOSE. THIS WARRANTY STATES PACE TECHNOLOGIES ENTIRE AND EXCLUSIVE LIABILITY AND BUYER'S EXCLUSIVE REMEDY FOR ANY CLAIM FOR DAMAGES IN CONNECTIONS WITH PACE TECHNOLOGIES PRODUCTS. PACE TECHNOLOGIES WILL IN NO EVENT BE LIABLE FOR INCIDENTAL OR CONSEQUENTIAL DAMAGES WHATSOEVER, NOR FOR ANY SUM IN EXCESS OF THE PURCHASE PRICE.

2. LIABILITY CAP:

PACE Technologies maximum aggregate liability for loss and damage arising under, resulting from or in connection with the supply or use of the Equipment and Consumables provided under this purchase, or from the performance or breach of any obligation (s) imposed hereunder, whether such liability arises from any one or more claims or actions for breach of contract, tort, (including negligence), delayed completion, warranty, indemnity, strict liability or otherwise, unless otherwise limited by the terms hereof, shall be limited to one hundred percent (100%) of the purchase price.

3. DELIVERY:

Customer assumes and shall bear the risk of all loss or damage to the Products from every cause whatsoever, whether or not insured, and title to such Products shall pass to Customer upon PACE Technologies delivery of the Products to the common carrier of Pace Technologies choice, or the carrier specified in writing by Customer, for shipment to Customer. Any claims for breakage, loss, delay, or damage shall be made to the carrier by the Customer and Pace Technologies will render customer reasonable assistance in prosecuting such claims.

4. ACCEPTANCE:

Customer shall inspect the Products promptly upon receipt of delivery. Unless customer objects in writing within thirty (30) business days thereafter, customer shall be deemed to have accepted the Products. All claims for damages, errors, or shortage in Products delivered shall be made by Customer in writing within such five (5) business day period. Failure to make any claim timely shall constitute acceptance of the Products.

5. PAYMENT:

Customer agrees to provide timely payment for the Products in accordance with the terms of payment set forth on the reverse side hereof or in any proposal submitted herewith. If any payment is not paid on or before its due date, Customer shall pay interest on such late payment from the due date until paid at the lesser of 12% per annum or the maximum rate allowed by law.

6. DEFAULT:

If Buyer is in default (including, but not limited to, the failure by Buyer to pay all amounts due and payable to Seller) under the work or purchase order or any other agreement between Buyer and Seller, Buyer's rights under the warranty shall be suspended during any period of such default and the original warranty period will not be extended beyond its original expiration date despite such suspension of warranty rights.

7. MISCELLANEOUS PROVISIONS:

This agreement has been made in and shall be governed by the laws of the State of Arizona. All disputes arising under or relating to the purchase of the equipment shall be brought and resolved solely and exclusively in the State of Arizona, Pima County. These terms and conditions and the description of the Products on the reverse side hereof or in any proposal submitted herewith constitute the entire agreement and understanding of the parties with respect to this sale and supersede all prior and contemporaneous agreements or understandings, inducements or representations, expressed or implied, written or oral, between the parties with respect hereto. Any term or provision of this Agreement may be amended, and any observance of any term of this Agreement may be waived, only by a writing signed by the party to be bounds. The waiver by a party of any breach shall not be deemed to constitute a waiver of any other breach. Should suit be brought on this Agreement, the prevailing party shall be entitled to recover its reasonable attorneys' fees and other costs of suit including costs and attorneys' fees incurred on appeal or in collection of any judgment., errors, or shortage in Products delivered shall be made by Customer in writing within such five (5) business day period. Failure to make any claim timely shall constitute acceptance of the Products.

8. RESTOCKING FEE:

All Returns are subject to a restocking charge equal to 15% (fifteen percent) of the Invoice, unless the Goods are proved to be non-conformed by PACE Technologies.



Contents

	PAGE
Nameplate	i
Warranty	ii
1.0 Product Description	1
2.0 Shipping, Unpacking and Installation	6
3.0 Safety Guidelines	9
4.0 Start-Up and Operation	12
5.0 Maintenance	23
6.0 Trouble Shooting	24
7.0 Cutting Basics	25
8.0 Electrical and Mechanical Drawings	28

1.0 Product Description

The **MEGA-T300A** is an automatic / manual table feed or manual wheel feed metallographic abrasive cut-off machine for cutting materials ranging from soft aluminum metals to hardened tool steels. It is ideal for the metallographic laboratory, as well as for small industrial or production applications.

The **MEGA-T300A** is very robust and durable with its cast aluminum alloy and stainless steel construction. Featuring a corrosion-free T-slot table, the **MEGA-T300A** is a very versatile wheel / table feed bench top model metallographic cutter.

Powerful motor

The **MEGA-T300A** is equipped with a powerful belt driven servo motor. In addition, the motor has an inductive brake for a quicker response of halting the blade (approx. 2 seconds).

Automated table feed cutting


The **MEGA-T300A** is an easy-to-use automated/manual cutter with simple controls. Cutting can be accomplished with either wheel feed or table feed.

The **MEGA-T300A** also has removable side port panels for sectioning longer samples.



INSTRUCTION MANUAL

3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

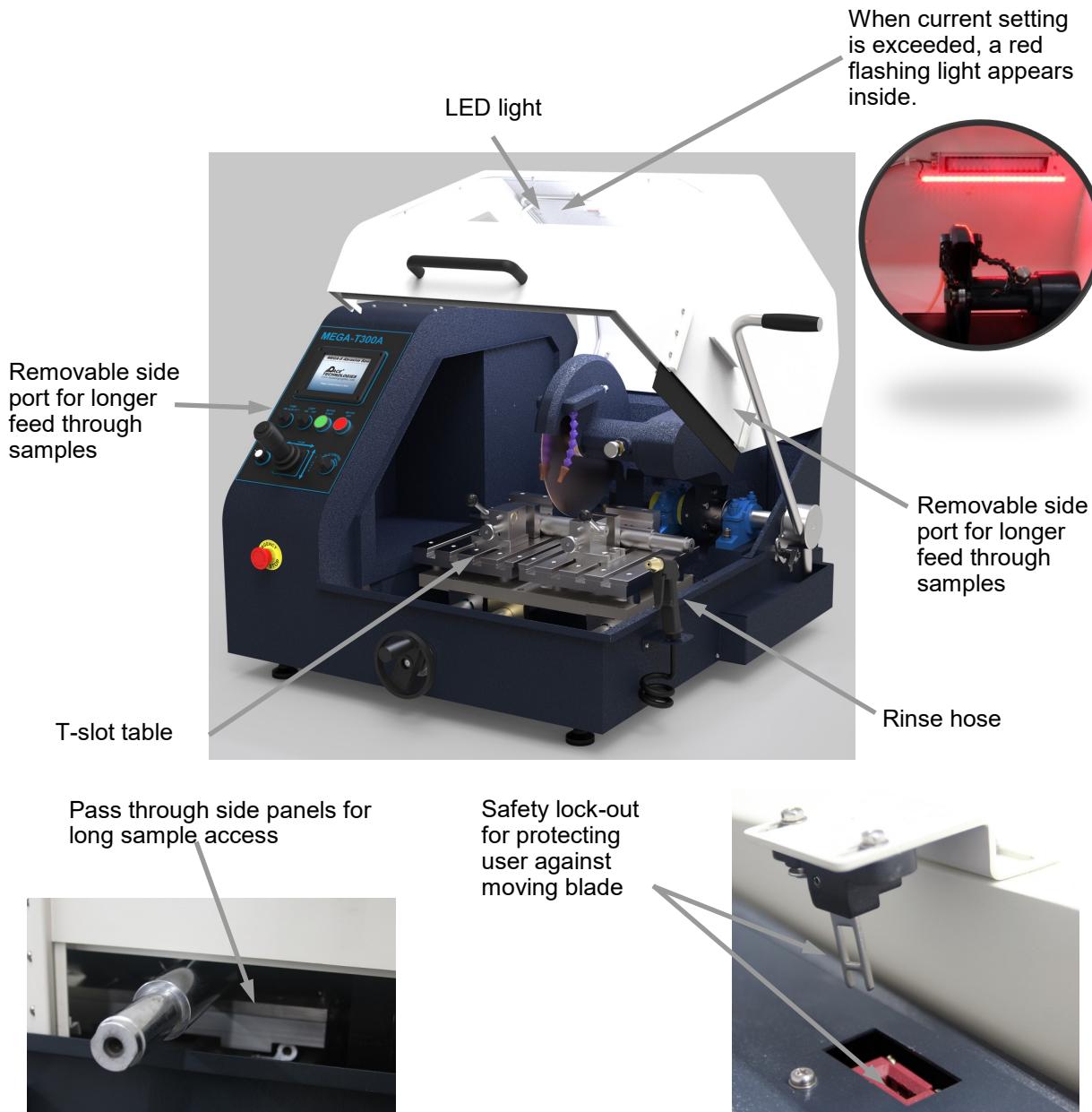

1. When closing the hood, press down on the hood to ensure that the safety switch is engaged.
2. After the cut has finished, if the hood lock does not release **DO NOT FORCE OPEN**, try pushing down on the hood to release the safety pin.




3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>



1.2 Technical Specifications


Electrical specifications:	208Y (3-phase or 120/120/120 to ground) or 380V 3-phase (Optional 480V to 380V 3-phase transformer)
Motor power:	Constant torque variable speed Servo Motor: -4.2 KW/(5.8 hp) 208Y -8 KW (10.6 hp) 380V -12 KW (16 hp) 480V
Speed (rpm):	Variable 1500-3500 rpm
Adjustable cutting load:	2-15 amps
Electronic motor brake	Stops in ~2 seconds
Cut-off wheels:	12-inch mm (~300 mm) diameter
Wheel arbor:	32 mm (~1.25 inch) diameter
Maximum sample diameter:	4-inch (100 mm) diameter 3 x 8-inches (75 x 200 mm) (HxD)
Wheel feed	Manual
Table feed	Automatic/Manual
Vertical movement	7-inch (180 mm)
Table feed movement (y-axis)	7.5-inch (190 mm)
Weight:	Approx. 500 lbs (120 kg)
Dimensions (WxDxH):	Hood Closed: 42 x 42 x 30 inch (1070 cm x 1070 cm x 760 cm) Hood Open: 42 x 42 x 36.5 inch (1070 cm x 1070 cm x 930 cm)
Table dimensions (WxD):	11 x 18.25-inch t-slot (280 x 470 mm)
Cutting chamber:	Cast aluminum block construction
Hood and body	Fabricated steel with shock resistant plastic window
Working temperature:	32° - 100°F (0 - 40°C)
Shipping temperature:	32° - 100°F (0 - 40°C)
Storage temperature:	32° - 100°F (0 - 40°C)
Recirculation system (included):	14 gallons (53 liters) EU Directives: CE Machinery directive 2006/42/EC RoHS Directive 2011/65/EU

1.3 Mechanical Schematic

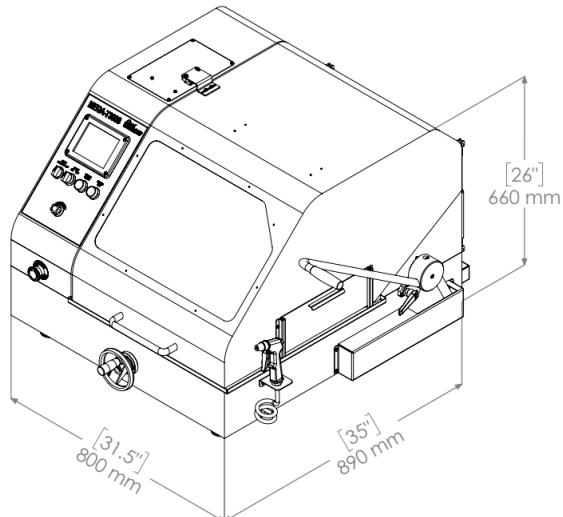

MEGA-T300A close hood dimensions

MEGA-T300A open hood dimensions

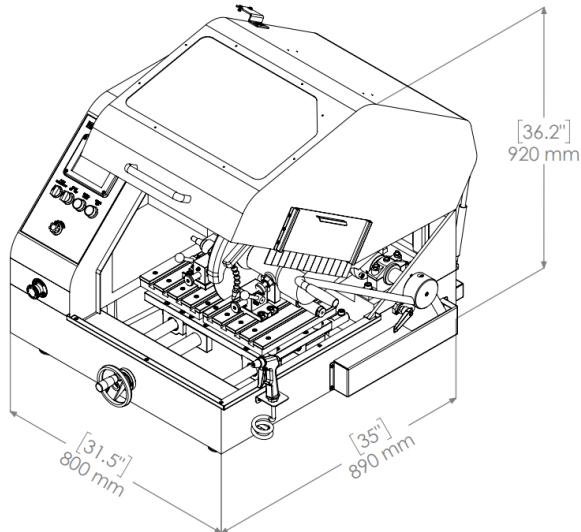

Note: Installation of the **MEGA-T300A** should be on a flat sturdy surface, with easy access to electrical connections.



3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

2.0 Unpacking, Shipping and Installation

2.1 Unpacking

Unit is delivered in a crate. Unpack and check for completeness of parts.

Measures WxHxD: Approx. 38 x 48 x 44-inch

Weight: Approximately 550 lbs

2.2 Shipping

When moving crate, lift from bottom. To move **MEGA-T300A** use lifting handles provided with the machine.

The **MEGA-T300A** is constructed of sensitive electronic and mechanical components. **DO NOT DROP.**

! **Caution:** Heavy equipment. Take care to avoid bodily injury.



(Installation continued on next page)

2.3 Installation

! Install unit carefully! Improper installation may void warranty.

The **MEGA-T300A** should be placed on a flat stable surface.

!


Connect coolant tank supply, drain and electrical connections.

Rinse Supply

Inlet Coolant

Drain Hose

Main Power
Pump Electrical

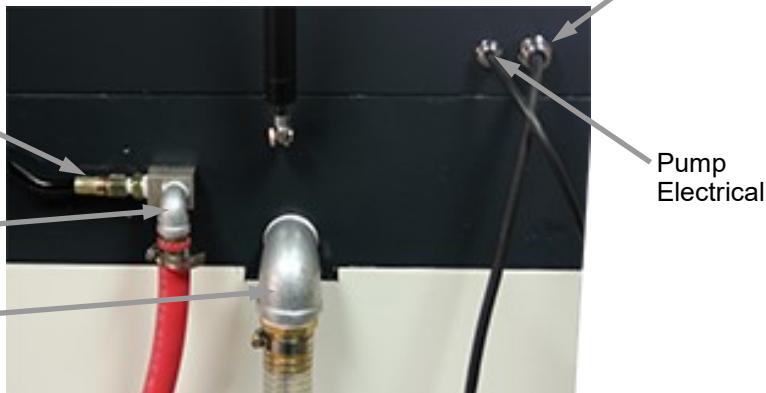

! **Electrical connections**—The standard saw is designed to operate at 208-220V 3-phase (120V-120V-120V-ground) (optional 380V). **Warning:** if not connected properly to 3-phase power damage may occur to the motor and void the warranty.

Verify voltage on nameplate



Main Power switch

! Verify the direction of rotation of the cut-off wheel. The wheel should turn from top to bottom as viewed from the front of the machine. If not reverse any two of the 3-phase power wires.

(Installation continued on next page)

Recirculating tank connections



External coolant supply: Attach 5/8-inch (17 mm) tube between pump and cutter.

Drain Hose: 1.5-inch (38 mm) hose.

Electrical connections: Connect electrical power cable to source.

Note: Inspect the operating voltage on the name plate.

Electrical connection for external coolant supply: Power for recirculation system comes from the **MEGA-T300S**,

Recommended Coolant: To minimize corrosion a coolant with a anti-corrosion additive or an oil emulsion coolant is recommended. **Note:** It is recommended that the coolant be changed every 50 cuts or once a week , whichever occurs first.

(Installation continued on next page)

3.0 Safety Guidelines

3.1 Warning Sign

! This sign points to special safety features on the machine.

3.2 Safety Precautions

Careful attention to this instruction manual and the recommended safety guidelines is essential for the safe operation of the **MEGA-T300A**.

Proper operator training is required for the safe operation of the **MEGA-T300A**. Any unauthorized mechanical and electrical change, as well as improper operation, voids all warranty claims. All service issues need to be reported to the manufacturer / supplier.

- ! Before operating, cutting chamber hood must be closed. After cutting, the safety latch will not open for approximately 2-8 seconds after pressing the emergency or red stop (red) button.
- ! Use only certified cut-off wheels from a professional supplier. Improper blades selection voids warranty. (For appropriate blade selection, refer to the Abrasive Blade Selection Guidelines Chart in Section 4.4)
- ! Disconnect power before opening the main unit.
- ! Replacement parts should be installed only by qualified personnel.
- ! Securely clamp the part / sample to the working table. During cutting, consider that the part may pinch and cause jamming of the cut-off wheel. Use the appropriate clamping devices to avoid this occurrence.
- ! Never start a cut under load.
- ! Make sure that the cut-off blade is rotating down and into the sample.


3.3 Emergency Statement

The **MEGA-T300A** abrasive cutter has been designed for cutting metallographic specimens. **DO NOT CUT** oversize samples. Always follow proper operational guidelines and avoid contact with moving parts, lubricants and abrasives. Seek appropriate medical care for cutting injuries.

Recirculation tank

The **MEGA-T300A** is equipped with a coarse filter to remove large particulate and broken blades

Note: Dispose of the old coolant in ordnance with federal, state and local regulations.

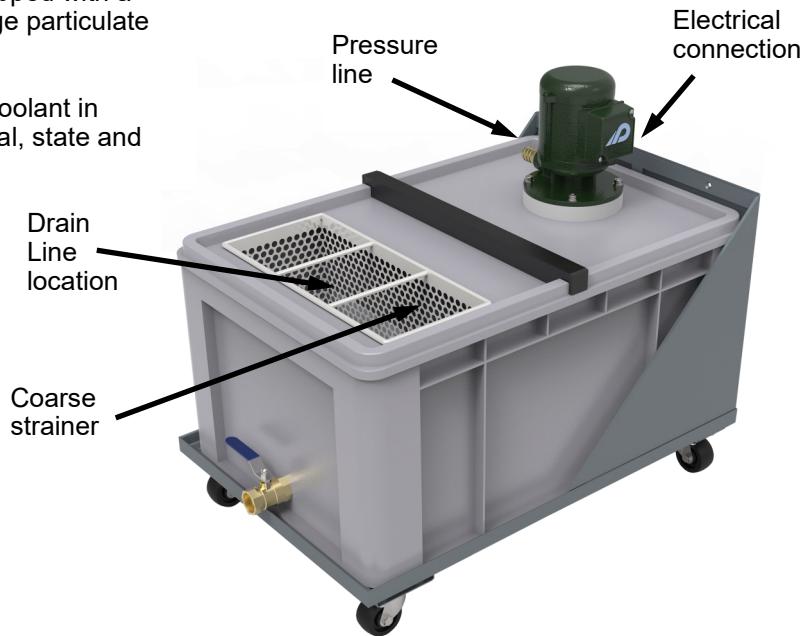

Rinse line

After use, water off the cutting chamber to prevent debris and corrosion build-up

! Leave the hood open to dry



3.4 Safety Tests



Emergency stop switch



Examine and verify that the **MEGA-T300A** safety devices and operating performance are in good working condition prior to use. The following safety checks are considered important:

Emergency stop switch

Test:	Activate main switch and close hood; depress emergency stop switch.
Proper Response:	Machine powers down and motor stops in 2-8 seconds
Malfunction:	Machine does not lose power.
Corrective measure:	If system does not power down, disconnect power supply cord and call service technician.



Magnetic Safety switch

Magnetic safety lock out switch

Test:	Activate main switch and close hood; turn cut-off wheel ON then OFF; Try to open the hood.
Proper Response:	Hood does not open for approximately 2-8 seconds.
Malfunction:	Hood opens immediately.
Corrective measure:	If cut-off wheel does not power down, disconnect power supply cord and call service technician.

HOOD (Closing / Opening)

1. When closing the hood, press down on the hood to ensure that the safety switch is engaged.
2. If the hood lock does not unlock **DO NOT FORCE OPEN**, or try pushing down on the hood to release the safety pin.

4.0 Start-up and Operation


4.1 General

The **MEGA-T300A** is a manual wheel and automatic / manual table feed cutter.

Note: When first powering up the **MEGA-T300A** the initialization time for the saw is approximately 15 seconds.

The servo motor also designed to increase speed slowly and typical takes 4-5 seconds to reach full speed.

4.2 Control Panel

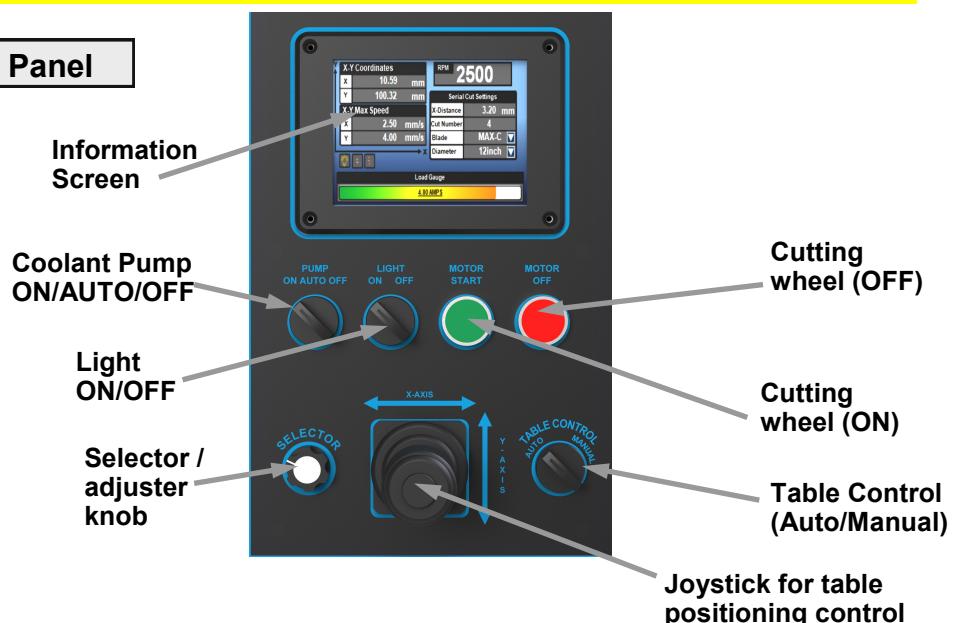

Cutting wheel ON/:	Starts the cutting procedure.
Cutting wheel OFF/:	Stops the cutting wheel (inductive brake engages to slow wheel). Hood cannot be opened until safety lock releases in approximately ~2-8 seconds after pressing the stop button.
Coolant ON/AUTO/OFF:	Operates the coolant pump in Auto or Manual mode. Must be selected before starting the motor.
Light ON/OFF:	Turns inside chamber light On or Off.
Selector / adjuster knob	Rotate to select feature, push to accept
TABLE CONTROL:	Select Auto for automated table cutting and Manual for manual table feed cutting
Joystick:	Controls the position of the table, X-Y axis

4.2.1 Control Panel Screen

START UP SCREEN



After turning on the main power switch the **START UP SCREEN** will be displayed.


Press the **SELECTOR** knob to progress to the **OPERATING SCREEN**



4.2.2 Control Panel Screen - Parameter Adjustments



To select the following parameters rotate the **SELECTOR** knob to highlight the parameter and then press the knob. Adjust the parameter by rotating the **SELECTOR** knob and press to accept

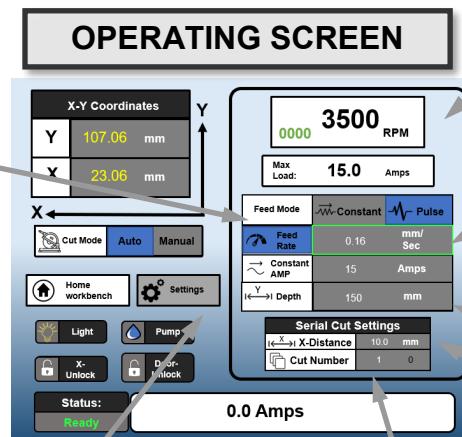

Select FEED MODE as constant or by pulsing based on parameters set in

SETTINGS adjustments for joystick Feed Rate, PULSE cutting force (amps), Overload protection (amps)

Set Operating

Set Y-axis FEED RATE

Sets the number of serial cuts

Set DEPTH of cut in y-
Sets the width of a serial cut adds the width of the blade as chosen in SETTINGS

INSTRUCTION MANUAL

3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

4.2.1 Control Panel Screen

START UP SCREEN



After turning on the main power switch the **START UP SCREEN** will be displayed.

Press the **SELECTOR** knob to progress to the **OPERATING SCREEN**




4.2.2 Control Panel Screen - Parameter Adjustments



To select the following parameters rotate the **SELECTOR** knob to highlight the parameter and then press the knob. Adjust the parameter by rotating the **SELECTOR** knob and press to accept

OPERATING SCREEN

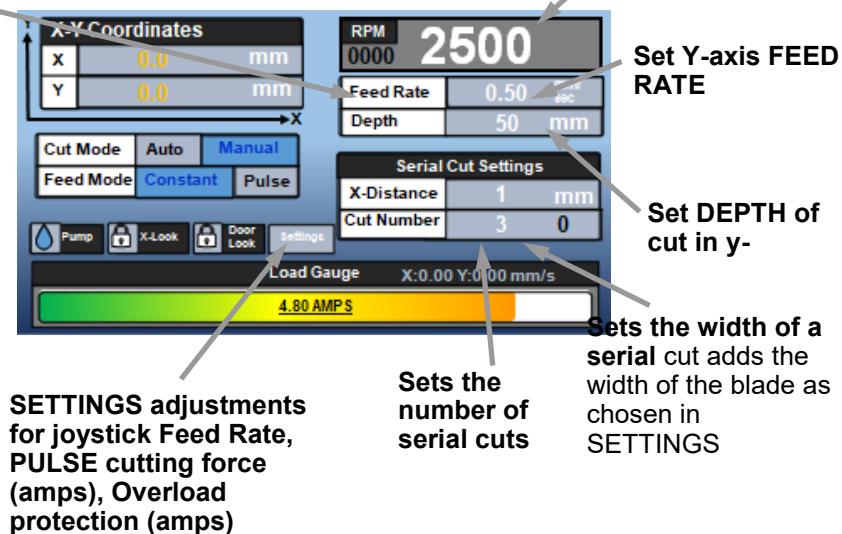

Select FEED MODE as constant or by pulsing based on parameters set in


4.2.3 Control Panel Screen - Manual vs. Auto Cutting

To change the following parameters rotate the **TABLE CONTROL** switch




Select CUT MODE to AUTO or MANUAL Cutting

OPERATING SCREEN

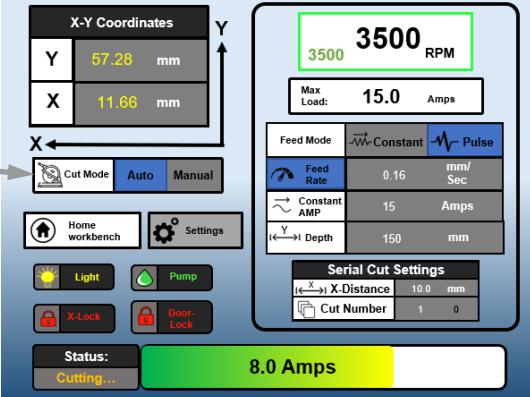


4.2.3 Control Panel Screen - Joystick Control

To change the position of the table use the Joystick

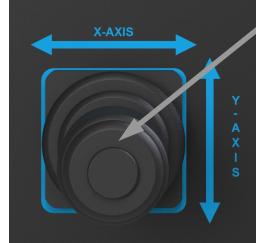

Adjust Joystick speed in the SETTINGS menu

OPERATING SCREEN

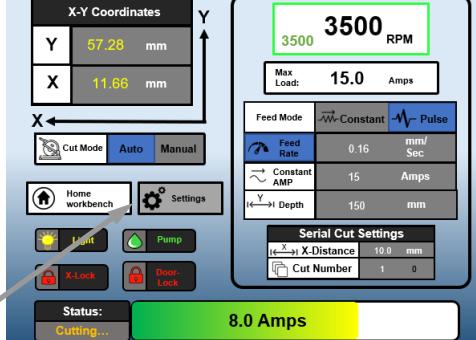


4.2.3 Control Panel Screen - Manual vs. Auto Cutting

To change the following parameters rotate the **TABLE CONTROL** switch




Select CUT MODE to AUTO or MANUAL Cutting

OPERATING SCREEN

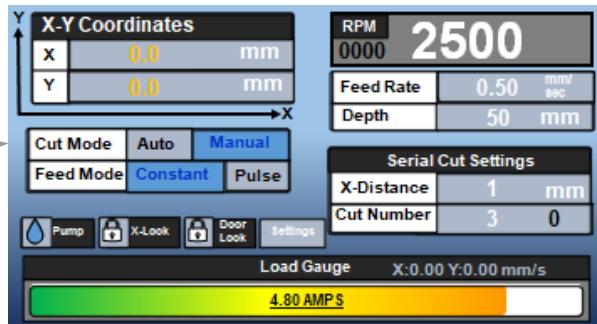


4.2.3 Control Panel Screen - Joystick Control

To change the position of the table use the Joystick

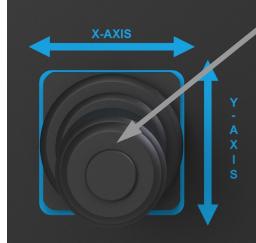

Adjust Joystick speed in the SETTINGS menu

OPERATING SCREEN

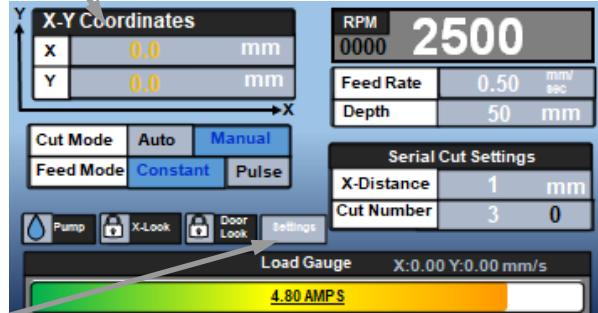


4.2.4 SETTINGS Screen - Joystick Control

Set Joystick Speed
 Maximum Values
 X-axis - 1.7 mm/sec
 Y-axis - 11.2 mm/sec

For Serial Cutting select present
 blade size or set CUSTOM
 thickness of blade


PULSE MODE
 Cutting force
 (amps)
 Cutting force
 feed rate and
 motor
 OVERLOAD
 (amps)

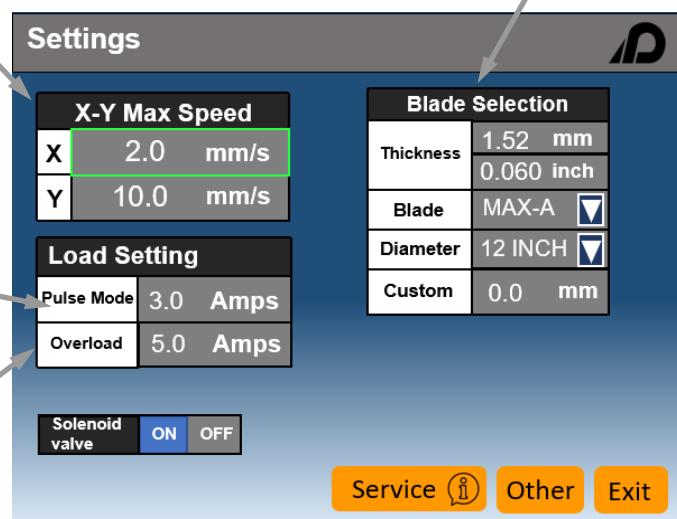

4.2.5 SETTINGS Screen - OTHER

For monitoring run time
 on cutter - useful for
 maintenance scheduling

Limit switch
 settings

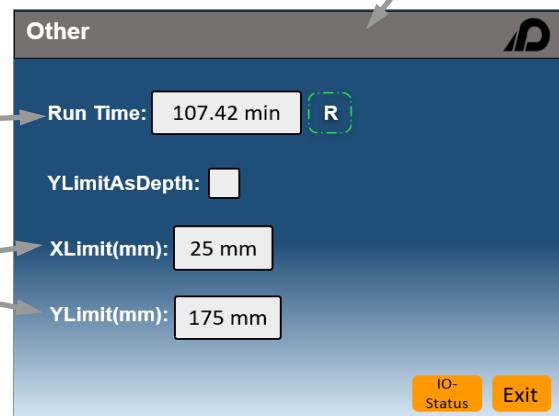

4.3 Wheel and Table Feed

4.3.1 Manual feed cutting




4.3.2 Table feed cutting




Open wheel spinner

4.3.2 Manual feed cutting

1. Move table to the front of the machine to install the appropriate blade
2. Properly fixture the specimen with the vises
3. Set the TABLE CONTROL switch to MANUAL
4. If cutting with the table feed, lower the cutter blade and lock into position. If manual cutting with the wheel, advance the table so the blade can advance all the way through the sample.
5. Set the cutting parameters on the SCREEN
 - CUTTING SPEED (rpm)
 - OVERLOAD Amps for controlling cutting force
6. Select the proper pump setting for the coolant
7. Close the hood
8. Press the GREEN Start Button
9. Either feed the table or cutting wheel into the specimen
10. If the OVERLOAD RED LIGHT flashes, reduce the cutting force
11. After completing the cut press the RED Stop Button
12. Allow blade and magnetic lock out to stop and disengage, respectively and then open the hood
13. Clean the cutting chamber and leave the

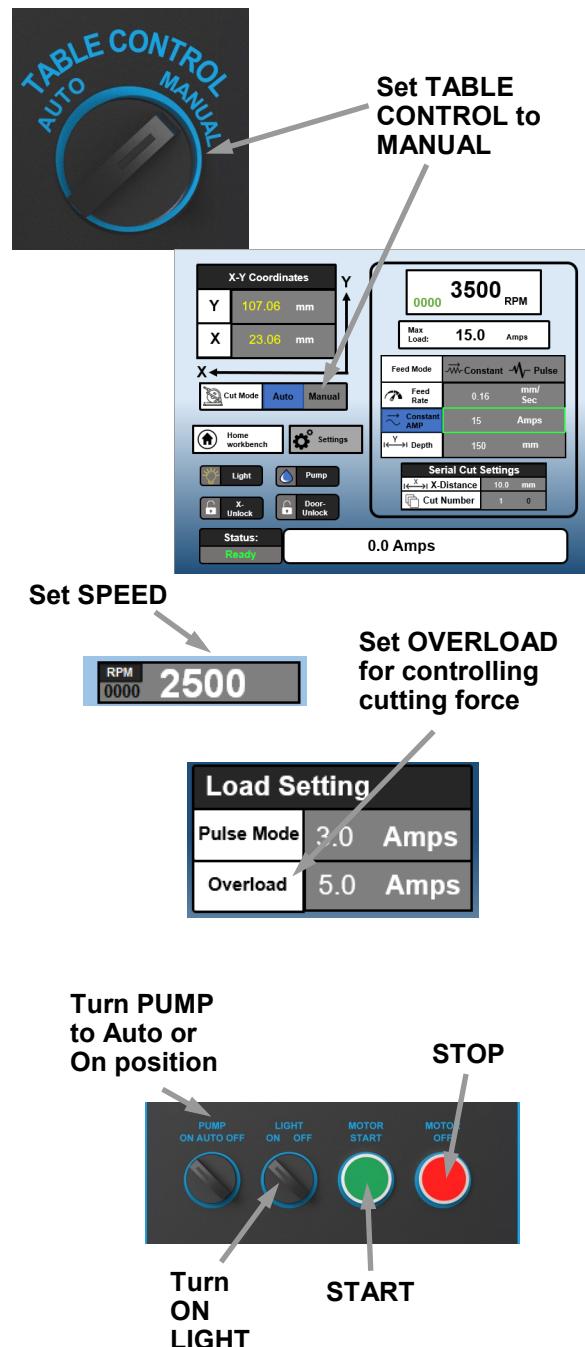

4.3.2 Automated Table feed cutting

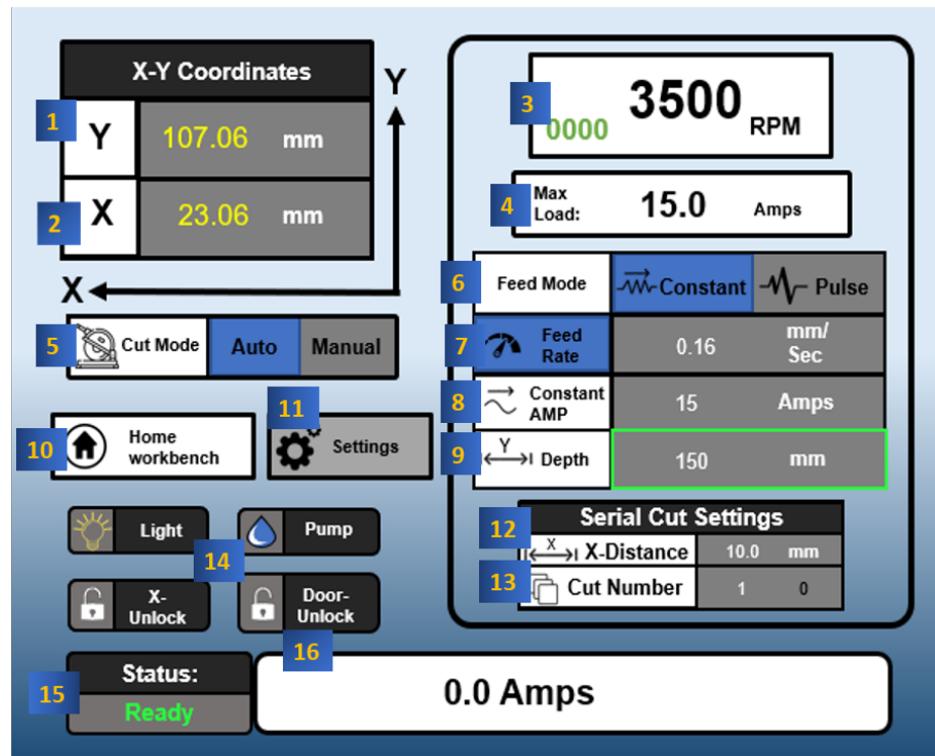

1. Y-axis coordinates of the table's current position
2. X-axis coordinates of the table's current position
3. RPM readout, Green displays is real-time rpm, the Black display is your current set RPM
4. Max AMP load setting
5. Cut mode, Auto or Manual, blue highlighted is the current selection
6. Feed Mode, Constant, or Pulse. Constant being a consistent feed rate, Pulse adheres to the amperage setting.
7. Feed Rate, Displays current feed rate of Y-Axis
8. Constant AMP, when Constant Mode is selected for Feed Mode this is the amperage setting for that mode.
9. Depth, Depth of cut you want to make in the Y axis
10. Home Workbench, returns the table to the home position (towards the front of the saw, closest to the user)
11. Settings, accesses the setting menu
12. X-axis serial cut width
13. Number of serial cuts to be performed
14. Saw's current status readout
15. Amount of amps the saw is currently drawing
16. RED BOX- Light on/off, Pump on/off, X-axis lock, Door lock/unlock status.



3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

4.3.2 Automated Table feed cutting (continued)

- 1. Move table to the front of the machine to install the appropriate blade**
- 2. Properly fixture the specimen with the vises**
- 3. Set the TABLE CONTROL switch to AUTO**
- 4. Lower the cutting blade and lock into position**
- 5. Set the cutting parameters on the SCREEN**
 - Wheel / blade SPEED (rpm)
 - FEED RATE (mm/sec)
 - DEPTH of cut (mm)
 - FEED MODE (constant or pulse mode) - Set PULSE MODE force (amps) in SETTINGS menu
 - For SERIAL CUTTING set the X-Distance for the width of the sectioned samples and adjust the number of cuts that you would like to make
 - Set OVERLOAD (amps) for controlling cutting force in SETTINGS menu. TIP: use left vise as a guide do not lock down so the sample can advance
- 6. Select the proper pump setting for the coolant**
- 7. Close the hood**
- 8. Position the sample by using the Joystick**
- 9. Press the GREEN Start Button to start the AUTO cutting operation**
- 10. The table will advance at the preset Feed Rate unless the sample cannot be cut at that rate. The amperage pull on the motor will increase as more force is put on the sample during cutting. If the maximum amperage exceeds the present value in SETTINGS the red overload light will flash and the table will backoff and then try again (similar to the PULSE Mode)**
- 11. The table will return to the position set with the Joystick and the motor and pump will turn off**
- 12. Allow blade and magnetic lock out to stop and disengage, respectively and then open the hood**
- 13. Clean the cutting chamber and leave the hood open when not in use**


4.4 Changing abrasive cut-off wheels

1. Remove blade locking bolt (reverse threaded).
 2. Position new 300 mm abrasive cut-off wheel into position.
- !** Use only certified abrasive cut-off wheels.
3. Gently tighten blade locking bolt (Note: locking bolt is reverse threaded).

Changing Blade

**Remove: Loosen by
turning CLOCKWISE**

**Tighten: Turn COUNTER-
CLOCKWISE**

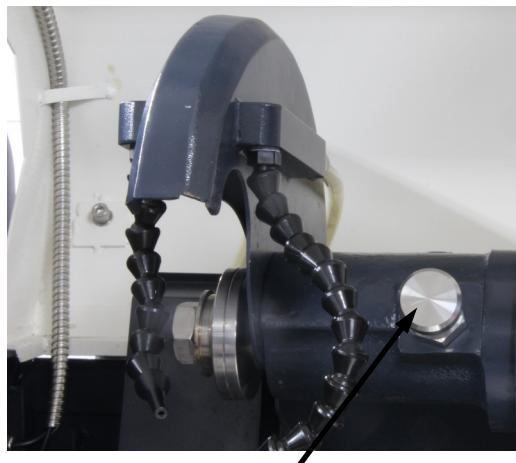

ABRASIVE Cutting CONSUMABLES SELECTION GUIDELINES CHART

Abrasive Cutting Fluids

<u>Pace Product Name</u>	<u>Catalog Number</u>	<u>Packaging</u>
MAXCUT™ Abrasive Cutting Fluid (32 oz)	MAXCUT 1000-32	32 oz
MAXCUT™ Abrasive Cutting Fluid (128 oz)	MAXCUT 1000-128	1 gallon
MAXCUT™ 2 Cutting Fluid with Anticorrosion Additive (32 oz)	MAXCUT2 1000-32	32 oz
MAXCUT™ 2 Cutting Fluid with Anticorrosion Additive (128 oz)	MAXCUT2 1000-128	1 gallon

300 mm (~12") MAXCUT™ Abrasive Blades 32mm (~1¼") Arbor

<u>Pace Product Name</u>	<u>Catalog Number</u>	<u>Packaging</u>
Soft Non-ferrous Materials (Aluminum, Brass, Zinc, Etc.)	MAX-E300	10
Hard Non-ferrous Materials	MAX-C300	10
Soft Steels	MAX-E30	10
Case Hardened Steels	MAX-VHS300	10
Steels and Ferrous Metals	MAX-D300	10
Universal Thin Resin/rubber Alumina Blade	MAX-A300	10



MAXCUT-A Blade



MAXCUT-C Blade



MAXCUT-D Blade



MAXCUT-D-RT Blade



MAXCUT-E Blade



MAXCUT-I Blade

INSTRUCTION MANUAL

3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

4.5 Serial Cutting

1. Serial cutting is used to make 1 or more uniform thickness sections. The cuts typically have +/- 0.2 mm tolerances, however, cutting force, blade type and material can affect these tolerances.
2. To use this function the following formula is used to determine the distance that the table will move in the x-direction:

Distance (x-axis) = Thickness required + Thickness of the Blade

3. The approximate thickness of the abrasive blades is listed in the table below and are listed in the SETTING table. For increased accuracy is recommended that the blade thickness be measured with a micrometer and the value added to the CUSTOM thickness screen.

Serial Cut Settings	
X-Distance	10.0 mm
Cut Number	1 0

Set number and width of serial cuts

Blade Selection	
Thickness	1.52 mm
	0.060 inch
Blade	MAX-A <input checked="" type="checkbox"/>
Diameter	12 INCH <input checked="" type="checkbox"/>
Custom	0.0 mm

Select blade type and diameter or manually enter thickness with CUSTOM

Blade Thickness	12-inch (300 mm)
MAX-A	0.055-inch (1.40 mm)
MAX-C	0.078-inch (1.98 mm)
MAX-D	0.078-inch (1.98 mm)
MAX-E	0.078-inch (1.98 mm)
MAX-VHS	0.078-inch (1.98 mm)
MAX-I	0.078-inch (1.98 mm)
DMAX	0.060-inch (1.52 mm)

4.6 Fixturing sample

For proper clamping, use the appropriate clamping vises to securely hold the sample in place. It is recommended that both sides of the part be clamped to avoid pinching of the blade (possibly breaking the blade) and to minimize burning of the workpiece during cutting.

! For proper fixturing, take into account the initial stress on the samples.

Fixture examples



Fast clamping

Left lever (Cat. No. QCL-1100)
Right lever (Cat. No. QCR-1100)


Note: vises have a slot in them so that the face does not rotate



UNIVERSAL Clamping vise
(Cat. No. MG-01) 45mm (1.8inch)



MEGA Saddle
Clamping vise
(Cat. No. M-702)

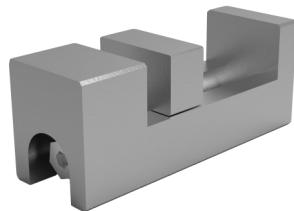

MEGA Fastener Vise
(Cat. No. MEGA-711)



5.0 Maintenance

5.1 Introduction

The recommended maintenance for the **MEGA-T300A** includes the following:

Important: Clean the table feed shaft after using the cutter to keep the abrasive from building up in the bushing. If the table feed begins to become more difficult use WD-40 to clean the feed shaft until the feed improves.

1. Use an anti-corrosion cutting fluid such as the MAXCUT 2 cutting fluid. It is recommended that the cutting fluid be changed every 50 cuts or weekly for maximum effectiveness.
2. Rinse and dry the table and vises after use.
3. When the cutter is not in use keep the hood OPEN to prevent humidity and corrosion build-up inside the cutting chamber.
4. If cutter will not be used for an extended period of time, coat the table and vises with a water displacement oil such as WD-40.
5. Use only recommended abrasive blades to minimize overheating the bearings and motor. Do not exceed the recommend amp rating of the motor.
6. Grease the appropriate fittings periodically (minimally after every 200 cuts or every other month, whichever occurs first).
7. Clean and apply a film of WD-40 on the table feed screw and bearing guides.
8. Periodically spray WD-40 between the inner arbor flange and the bearing housing to keep debris from building up.
9. Immediately replace any broken or damaged parts.

Note: corrosion is considered a maintenance issue and not a warranty issue

5.2 Cleaning outside cabinet

The cabinet and front shield should be cleaned occasionally with a moistened cloth. Do not use any chemicals or cleaning abrasives. We recommend an anti-condensation gel be used on the front panel.

5.0 Maintenance

5.3 Instructions for Maintenance X-Axis

The recommended X-Axis maintenance for the **MEGA-T300A** includes the following:

Important: Clean the table feed shaft after using the cutter to keep the abrasive from building up in the bushing. If the table feed begins to become more difficult use WD-40 to clean the feed shaft until the feed improves.

1. Clean periodically the bottom of the workbench to keep debris out of the X-Axis motor and transmission.
2. Clean the top of the table base after a couple of cuts to avoid debris going into the X-Axis screw.
3. Make sure to clean rails and guides to not over force the X-Axis motor.
4. Make sure nothing interferes with the X-Axis movement when is going back to origin to not stuck the motor.
5. Clean the detection surface of the sensors from the top to avoid interferences. (Applies for both axis Y/X)

In case of disassembly of the workbench follow these directions:

1. Make sure the M5 Hex screw that holds the X-Axis screw with the transmission is tight.
2. Align properly the transmission on the top plate.
3. Make sure the wiring in the bottom is not stressed and it's covered properly (in case replace the motor is needed)
4. Make sure the X-Axis bracket was properly installed on the bottom of the X-Axis plate and the bolts are properly tightened.

Note: corrosion is considered a maintenance issue and not a warranty issue

▲ INSTRUCTION MANUAL

3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

6.0 Trouble Shooting

Problem	Cause	Solution
No power or function	a. Unit is disconnected from main electrical power supply. b. Main power switch is off. c. Emergency stop button engaged.	a. Verify electrical source and connection. b. Turn on main power switch. c. Release by turning clockwise.
Main motor does not operate	a. Hood not closed. b. Overload relays activated	a. Close hood b. Restart after resetting relay.
Pump motor does not operate	a. Plug is disconnected	a. Check plug.
Excessive vibration during cutting	a. Not cutting with enough force b. Specimen locked too far back on cutting table c. Incorrect blade selection	a. Increase cutting force. b. Move specimen forward. c. Use correct blade.
Corrosion of table of vises	a. Improper cleaning and drying b. Incorrect or old cutting fluid. c. Cutting chamber has become a humidity / corrosion chamber d. Long term storage oxidation	a. Rinse and dry table and vises. b. Use an anti-corrosion cutting fluid such as the MAXCUT 2. It is recommended that the cutting fluid be changed every 50 cuts or weekly for maximum effectiveness. c. Open hood when not being used d. If cutter will not be used for an extended period of time, coat the table and vises with a water displacement oil such as WD-40.
ERROR 01*	Sudden stopping of motor after red light flashes a. Feed rate too high b. Amperage draw is too high c. Improper blade selection Sudden stopping of motor with no red light warning a. Feed rate too high b. Amperage draw is too high c. Improper blade	Turn off main power for several seconds to reset and correct problem: a. Reduce feed rate b. Decrease cutting force c. Use correct blade. Turn off main power for several seconds to reset and correct problem: a. Reduce feed rate b. Decrease cutting force c. Use correct blade.
Hood does not open when motor is off	a. Trying to open the hood before the safety switch releases (approximately 2-3 seconds) b. Safety switch still does not disengage	a. Wait until the safety lock out clicks b. Try pushing down on the hood to release the safety pin in the safety lockout switch
Motor stops immediately after starting	a. X-axis table feed is activating the sensor	a. Move the X-axis table to the right so the x-coordinate is less than the maximum value (approx. 22 mm)

* If the motor stalls and an ERROR 01 is displayed, power off the cutter and wait until you here a small click from the driver for it to reset, typically <15 seconds.



6.1 Trouble Shooting Error Codes

Error Code	Code Description	Solution
ERR 01	IPM Error	Disassemble and return to supplier
ERR 02	Over current	Check the setting parameter
ERR 03	EEPROM Error	Rewrite the default value
ERR 04	Current sampling	Rewrite parameters
ERR 05	Bus under or over voltage	Check the voltage and parameters
ERR 06	Communication anomaly	Check the wire connection
ERR 07	Locate failed	Check the encoder signal; Or the current feedback is wrong reposition
ERR 08	Lost step when running	Normally, the reverse angle is 2 degrees and the reverse angle is 20 degrees when positioning will be wrong reposition
ERR 09	Overload	Cutting too fast, too much load, reduce the speed.
ERR 10	Module temperature is too high	Check whether the fan works properly
ERR 12	Power error	Disassemble and return to supplier
ERR 13	Encode Failure	Check wire connection or the encoder was damaged
ERR 14	Motor get stuck	Check the motor or cutting speed is too high

7.0 Metallographic Cutting Basics

Philosophy:

1. Cutting is the most important metallographic sample preparation step
2. Minimize damage during sectioning
3. Less sectioning damage = reduced number of grinding and polishing steps


Note: In some cases cutting damage cannot be removed by subsequent polishing operations (especially for brittle samples)

Time Well Spent:

1. Spend more time on sectioning in order to decrease overall sample preparation time
2. Protect integrity of the specimen with proper sectioning

Wheel feed


- High speed (3500 rpm)
- High load (15 amps)

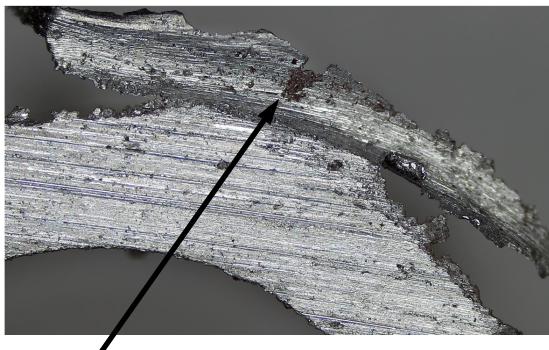

Coarse finish with
large cutting burr

Table feed

- Low speed (1500 rpm)
- Low load (2-5 amps)

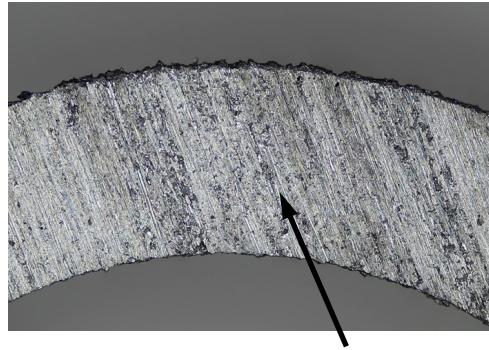


Finer finish with
minimal cutting burr

7.0 Metallographic Cutting Pentagon

Five variables for successful sectioning:

- **Abrasive Type**
 - Silicon Carbide
 - Alumina
- **Abrasive Bond**
 - Resin
 - Rubber
- **Cutting Force**
 - Control by motor current
 - Cutter Speed
Low (1500-2000 rpm)
 - Medium (2000-2500 rpm)
 - High (2500-3500 rpm)
- **Coolant**
 - Water-based
 - Oil Emulsion-based

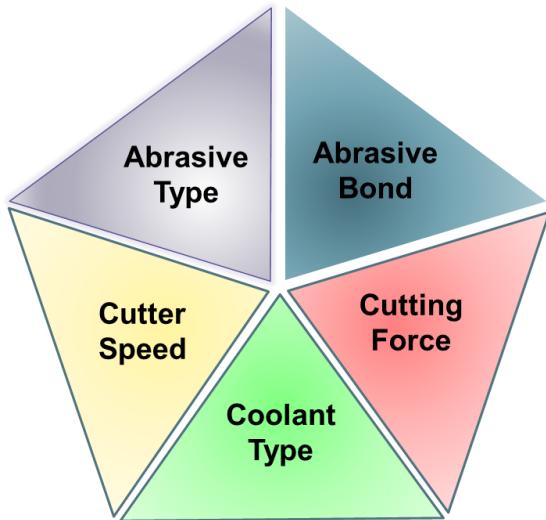

Metallographic Abrasive Blades:

- **Resin Bonded Blades**
 - Breakdown easily
 - Less smell and burning
 - More versatile
 - Alumina abrasives
- **Rubber Bonded Blades**
 - Longer life
 - Produces a burnt rubber smell
 - More likely to burn the sample
 - More difficult to match the blade to the material being cut
 - SiC and Alumina abrasives




Machine Parameters:

- **Wheel Speed**
 - Slower speed for softer metals (1500-2000 rpm)
 - Medium speed for most metals (2500-3000 rpm)
 - Higher speed for harder samples (3500 rpm)
- **Cutting Force**
 - Controlled by limiting current to the motor
 - Softer metals lower load (2-5 amps)
 - Most metals medium load (6-10 amps)
 - Harder samples higher loads (11-15 amps)


Recommended Cutting Procedures:

Material	Abrasive / Bond	Blade	Cutting Speed	Cutting Force
Soft non-ferrous metals (aluminum, brass, zinc, etc.)	Alumina / resin	MAX-E	1500 rpm	2-5 amps
Hard non-ferrous metals (titanium, zirconium, etc.)	SiC / rubber	MAX-C	2500 rpm	5-10 amps
Soft steels	Alumina / resin	MAX-E	1500 rpm	5-10 amps
Hard and case hardened steels	Alumina / resin	MAX-VHS	2500 rpm	10-15 amps
General purpose steel and ferrous metals	Alumina / resin	MAX-D	2500 rpm	10-15 amps
Universal thin blade	Alumina / rubber	MAX-A	2500 rpm	10-15 amps
Industrial longer life blades	Alumina / resin	MAX-I	3500 rpm	10-15 amps
Ceramics and glass	Diamond / metal	DMAX	3500 rpm	10-15 amps

8.0 Drawing

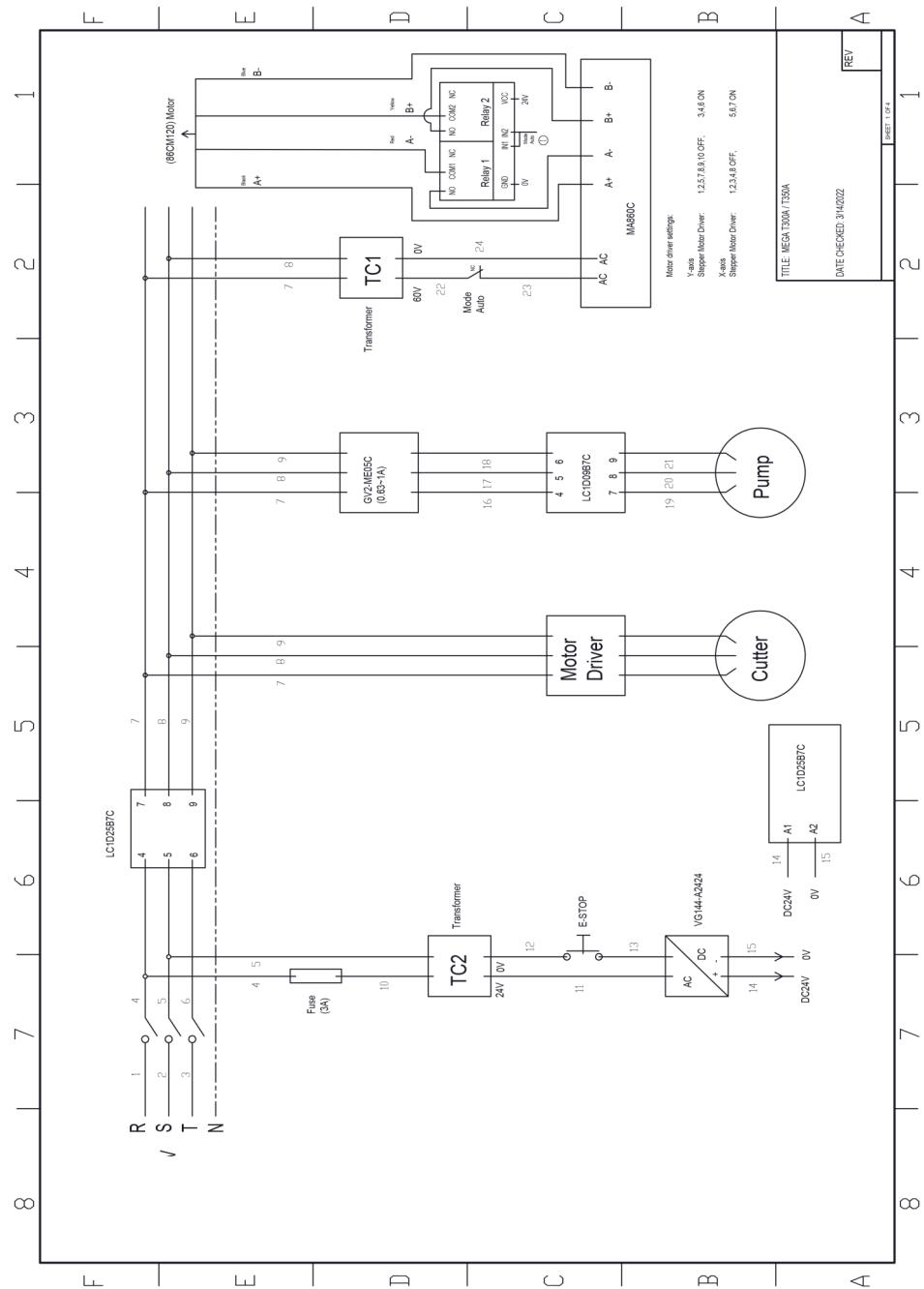

Please read this instruction manual carefully and follow all installation, operating and safety guidelines. 30

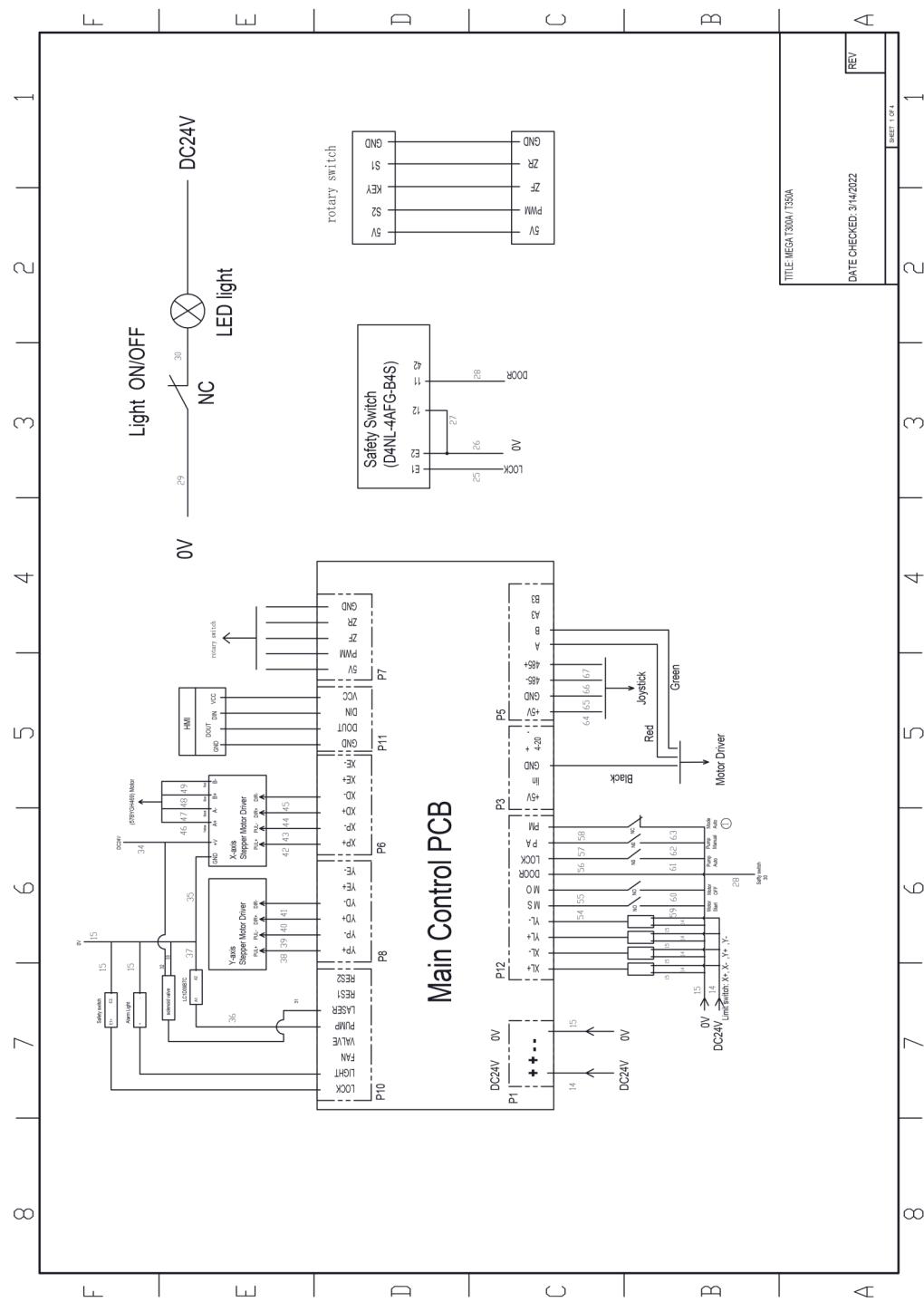


Please read this instruction manual carefully and follow all installation, operating and safety guidelines. 31





3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>



9.0 Exploded Drawings

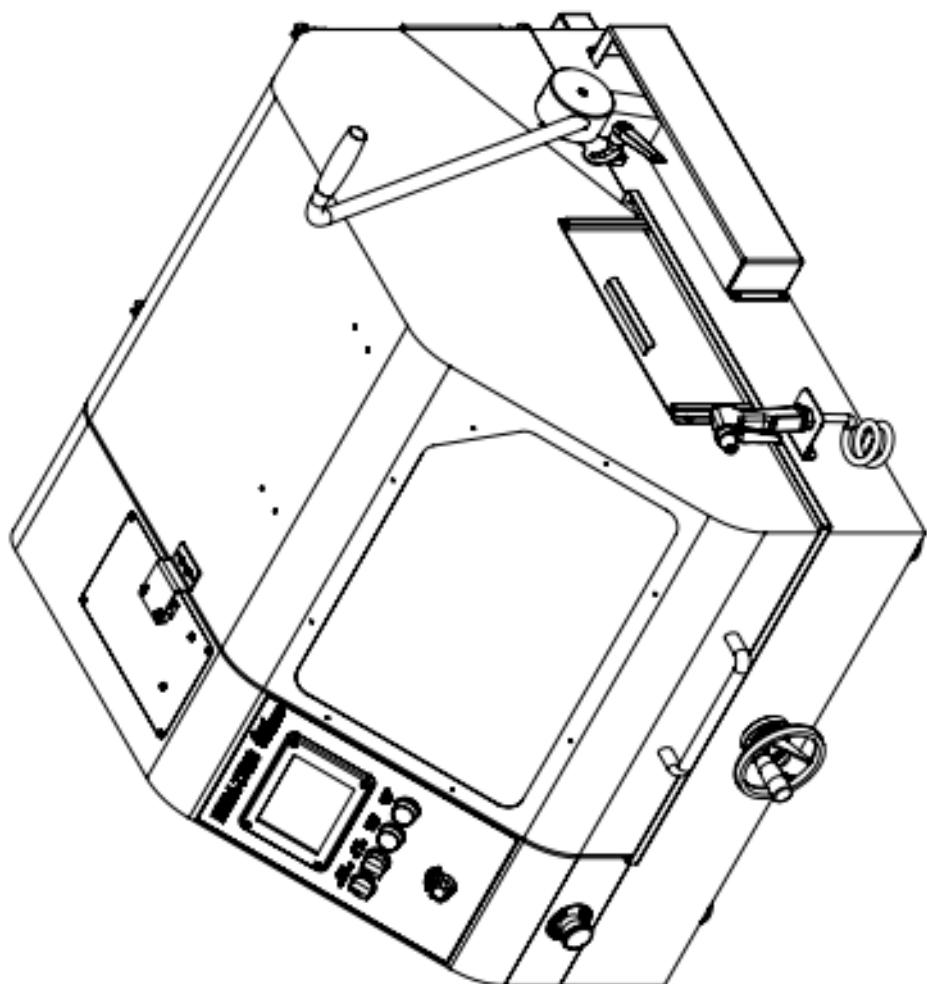






3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

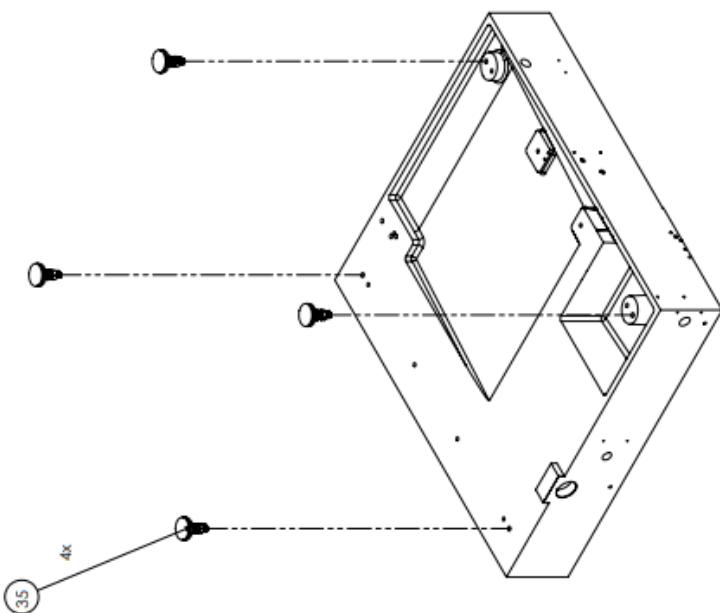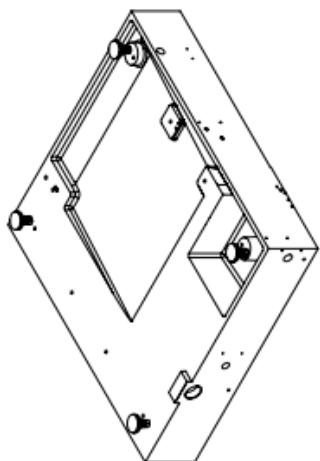

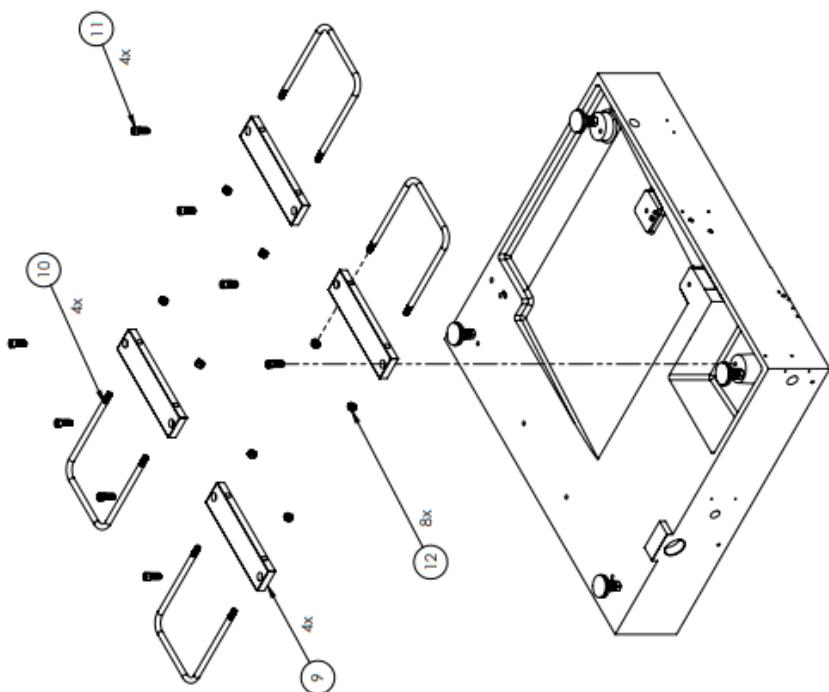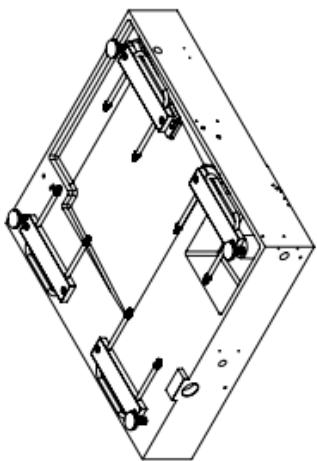





3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

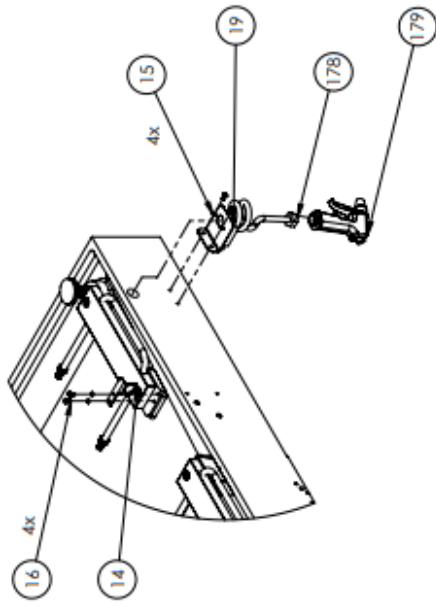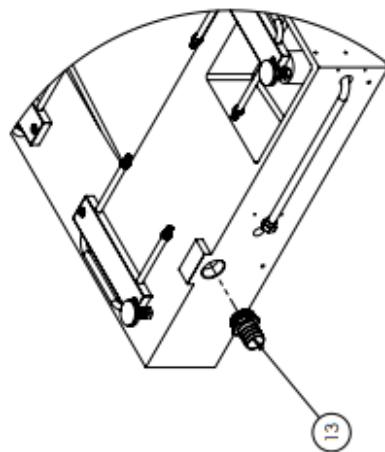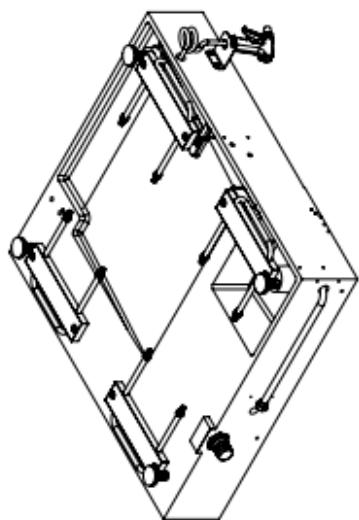







3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

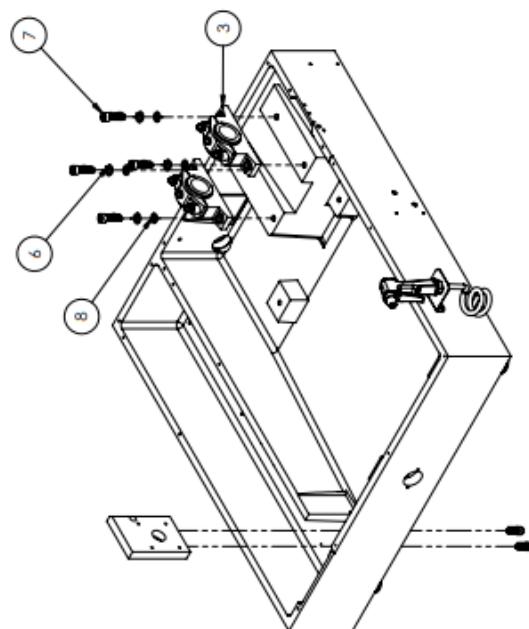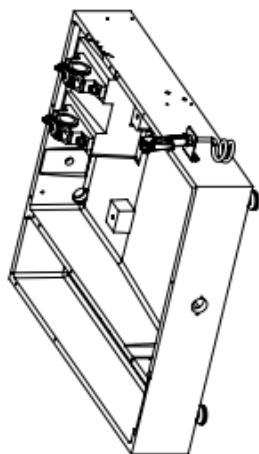





3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

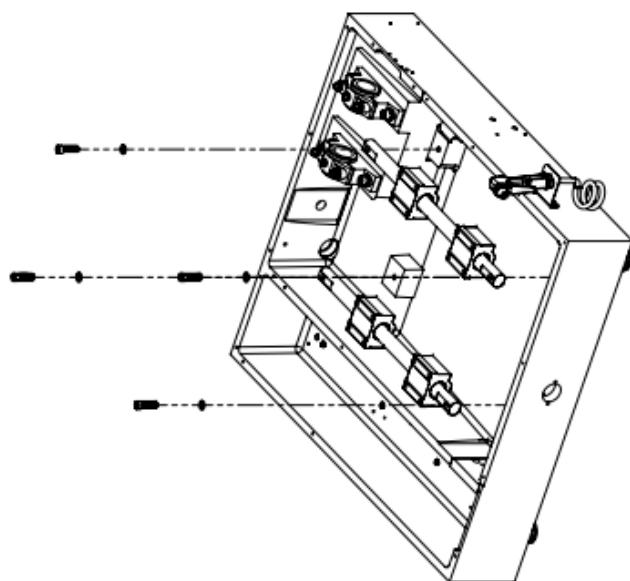
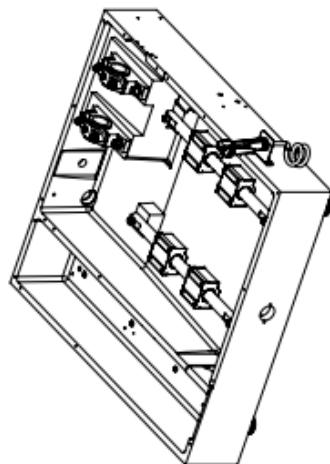


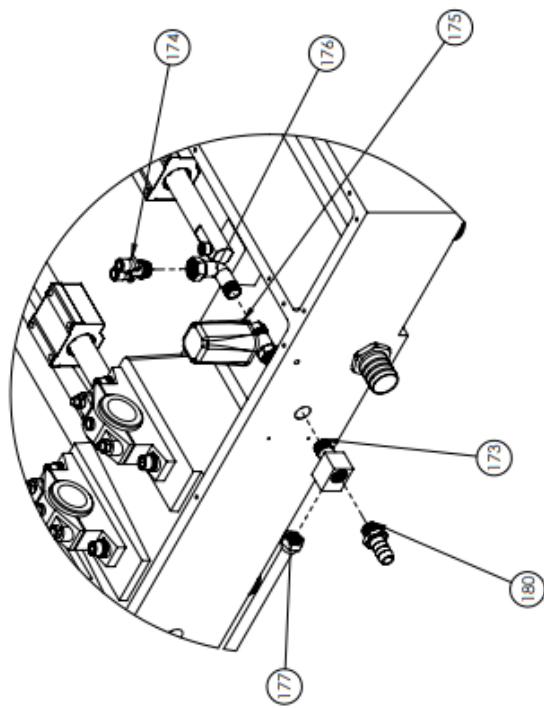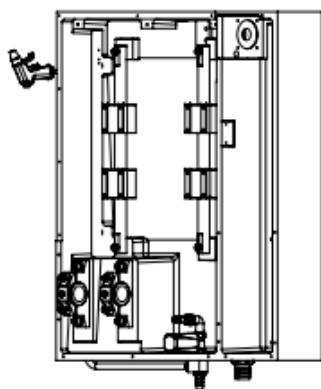




3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

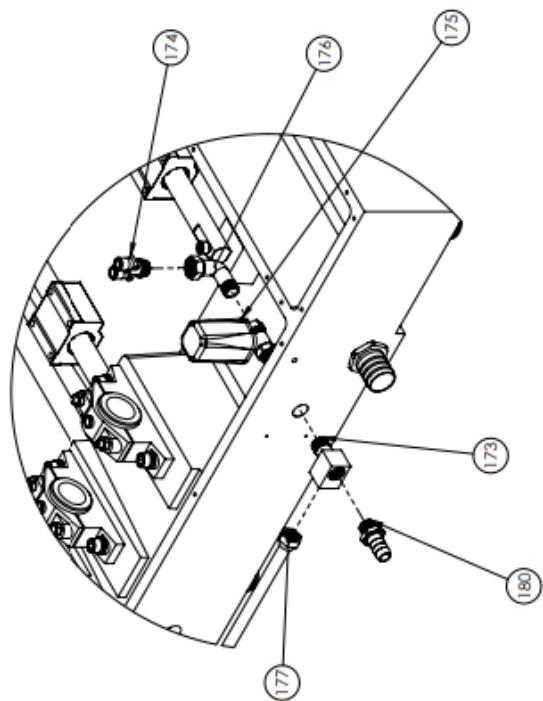
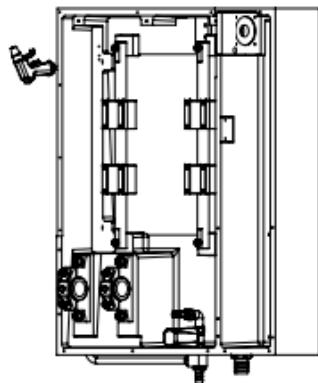


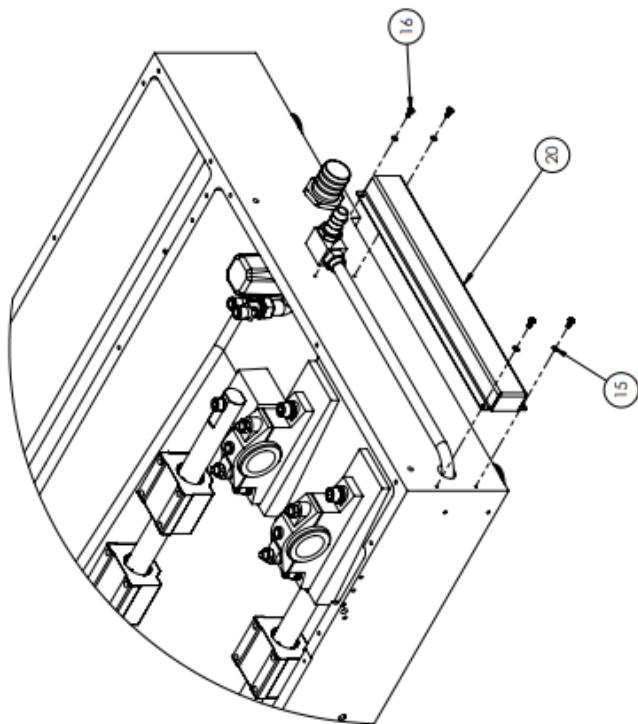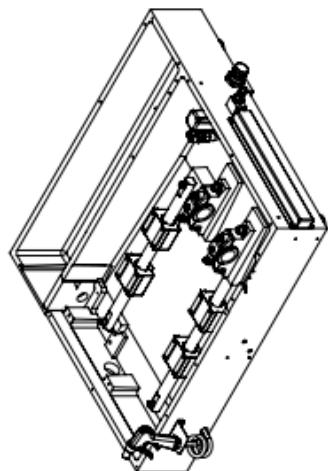

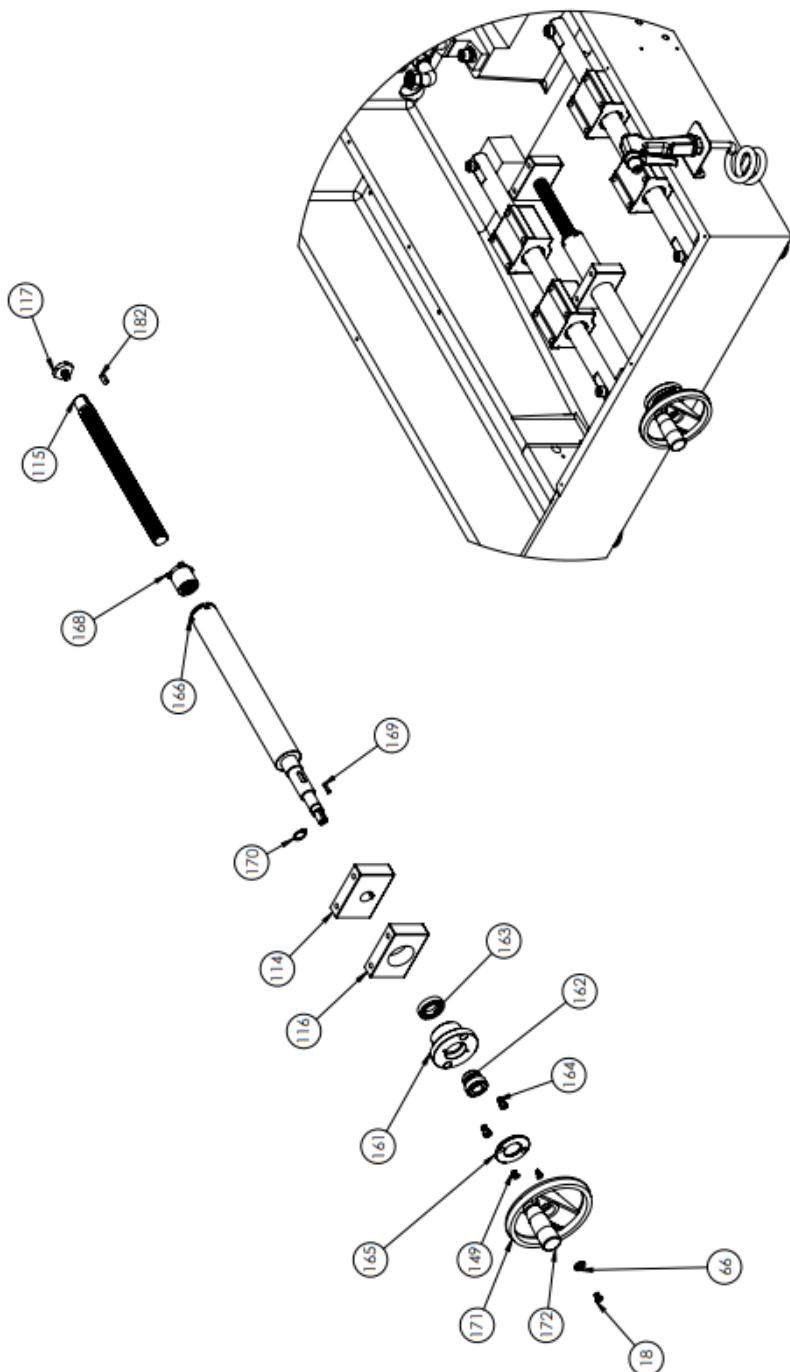






3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

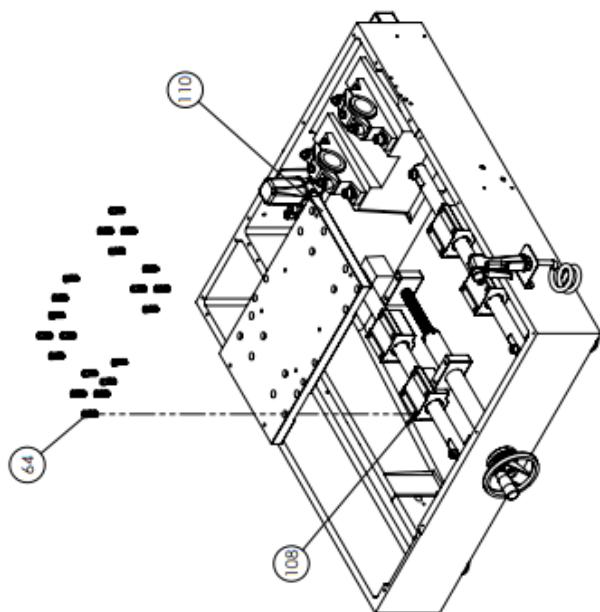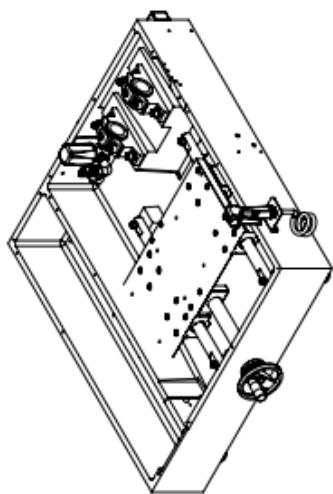





3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

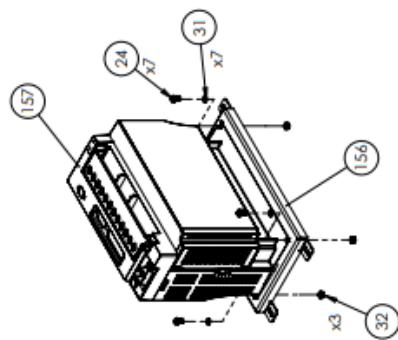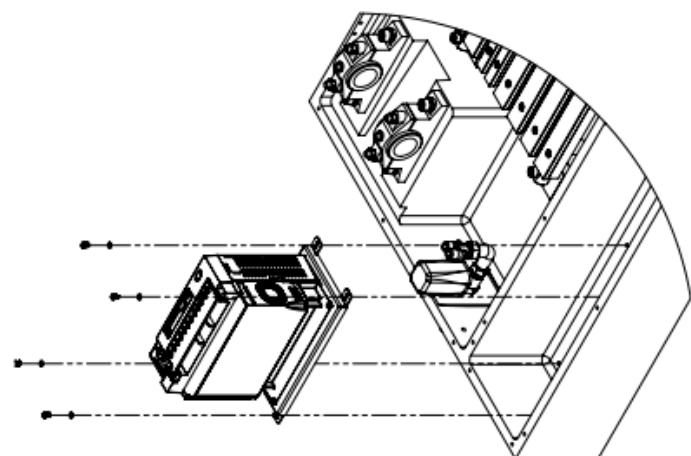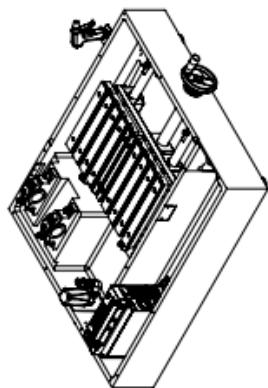






3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

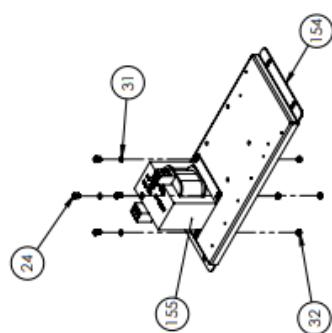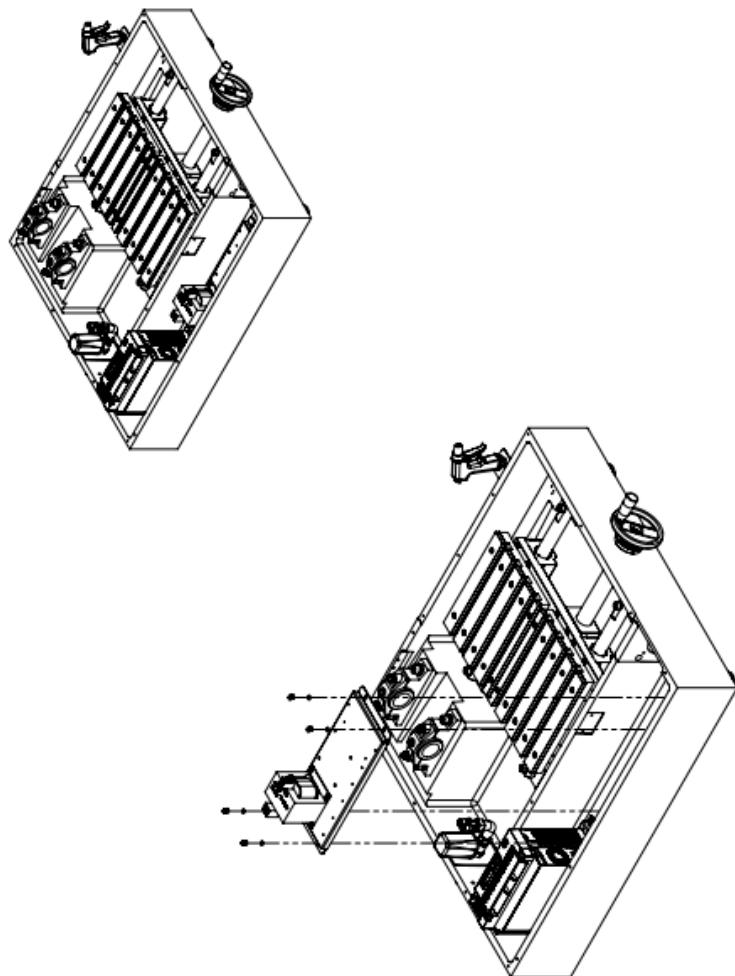





3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

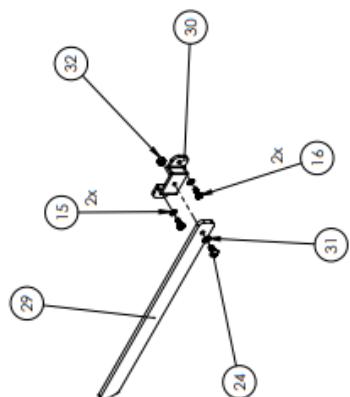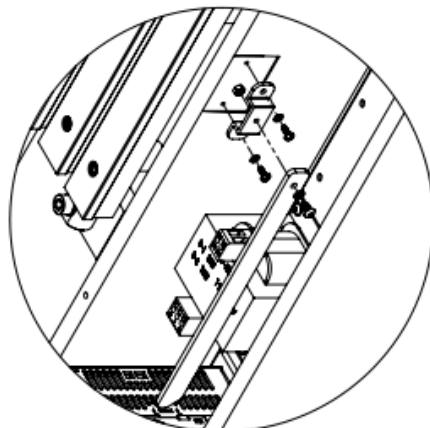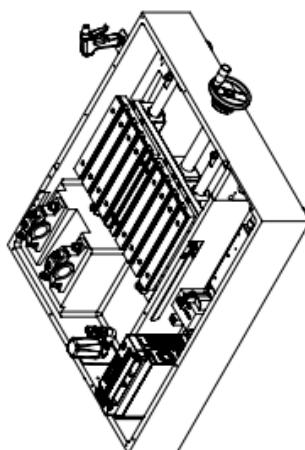





3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

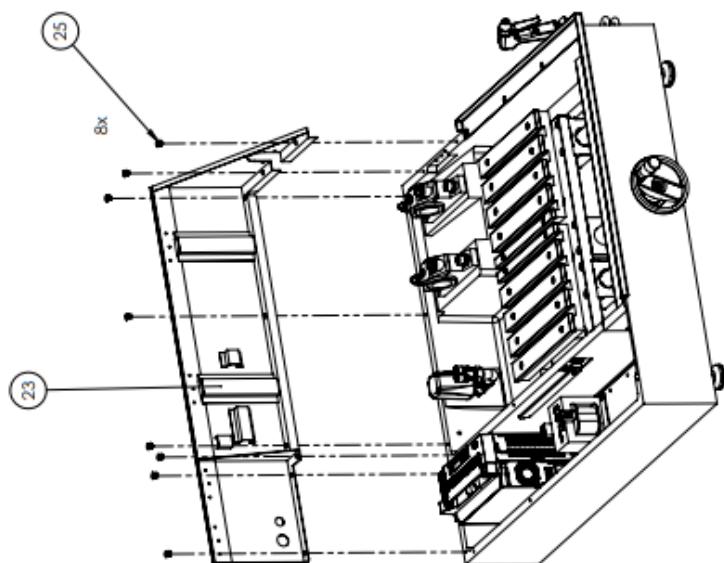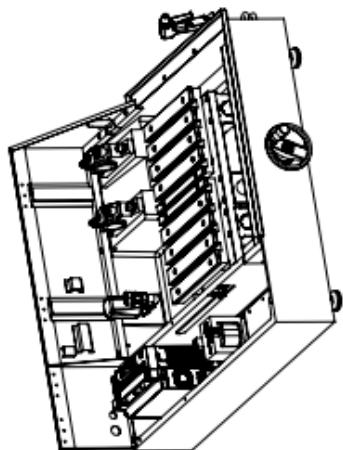




3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

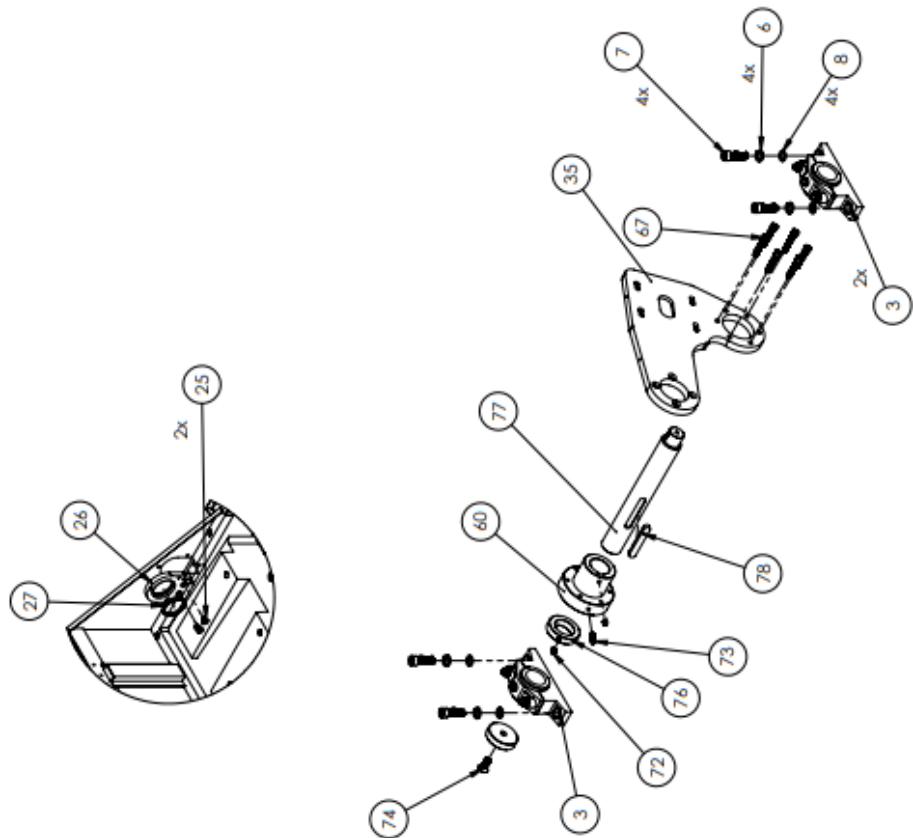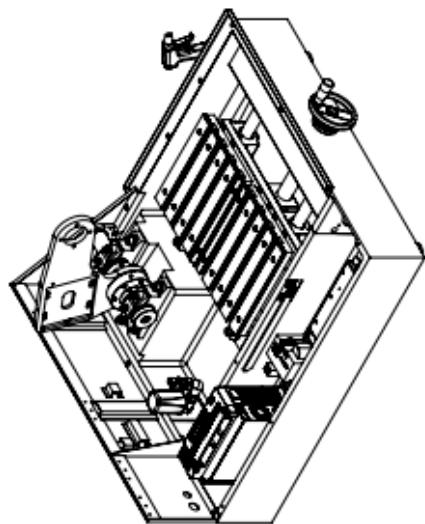





3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

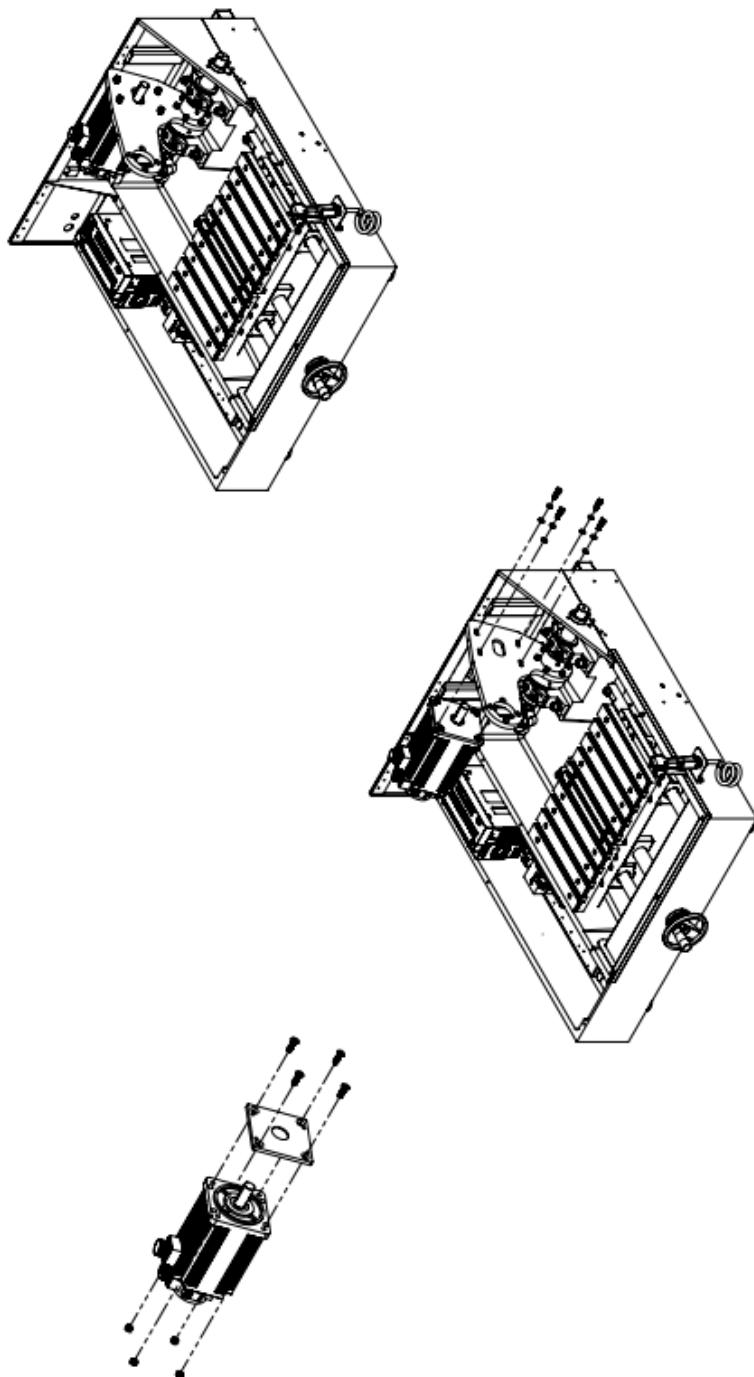




3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

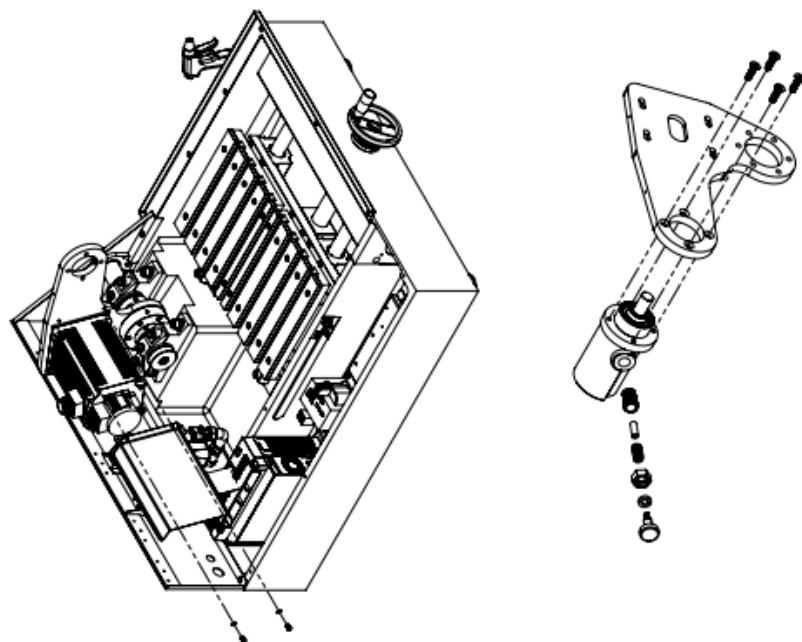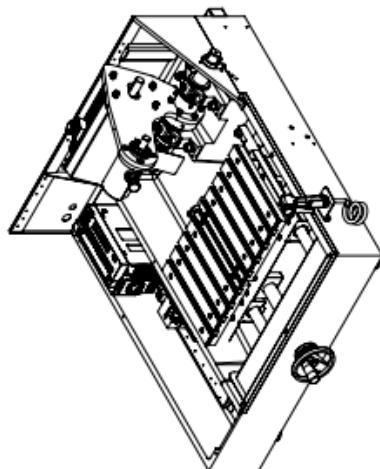




3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

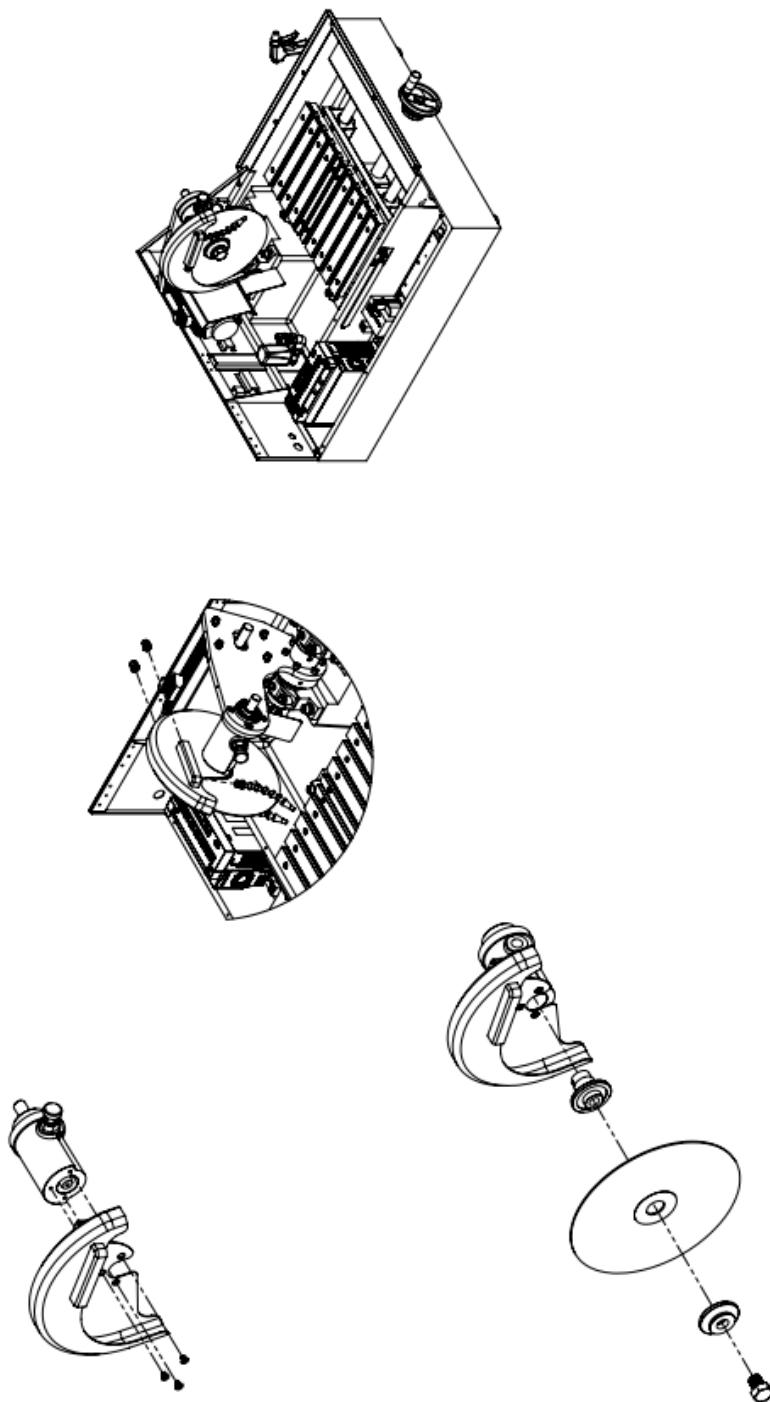




3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

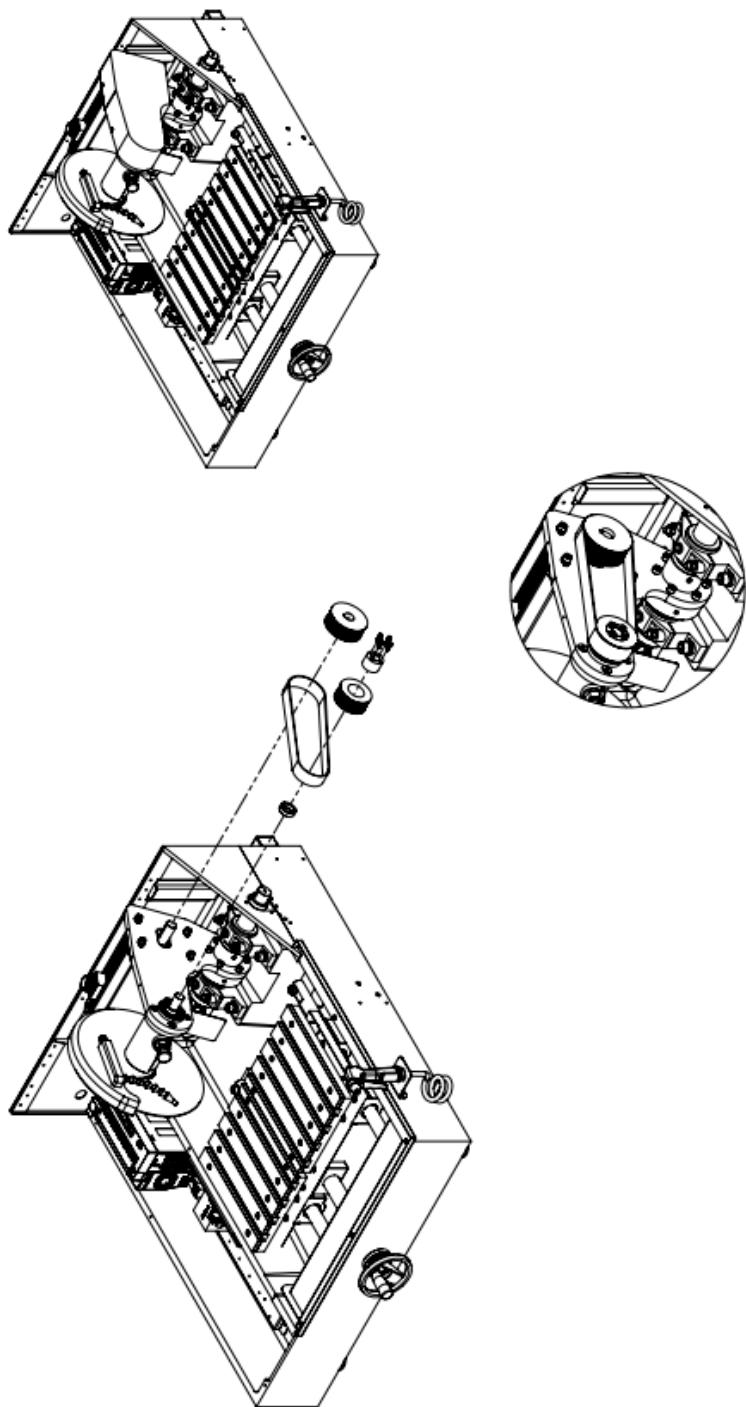




3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

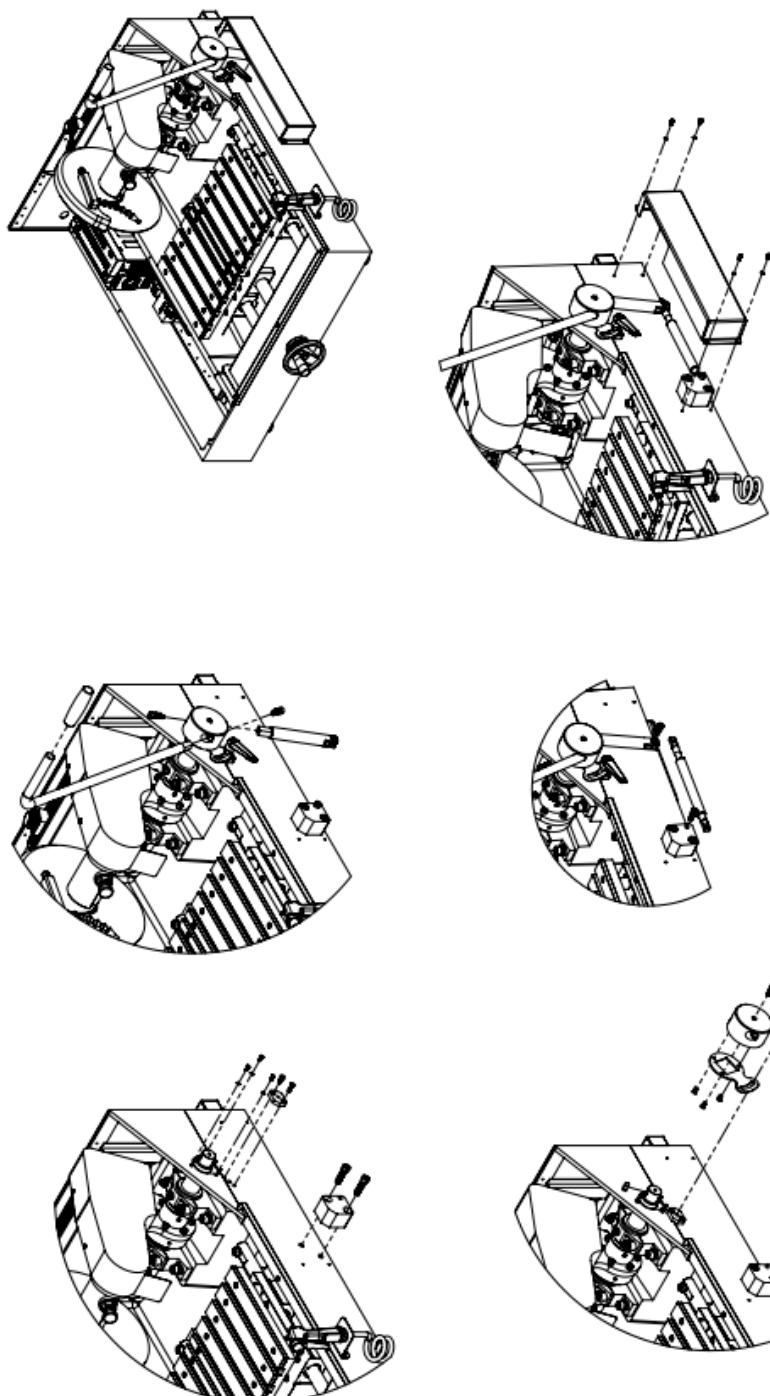






3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

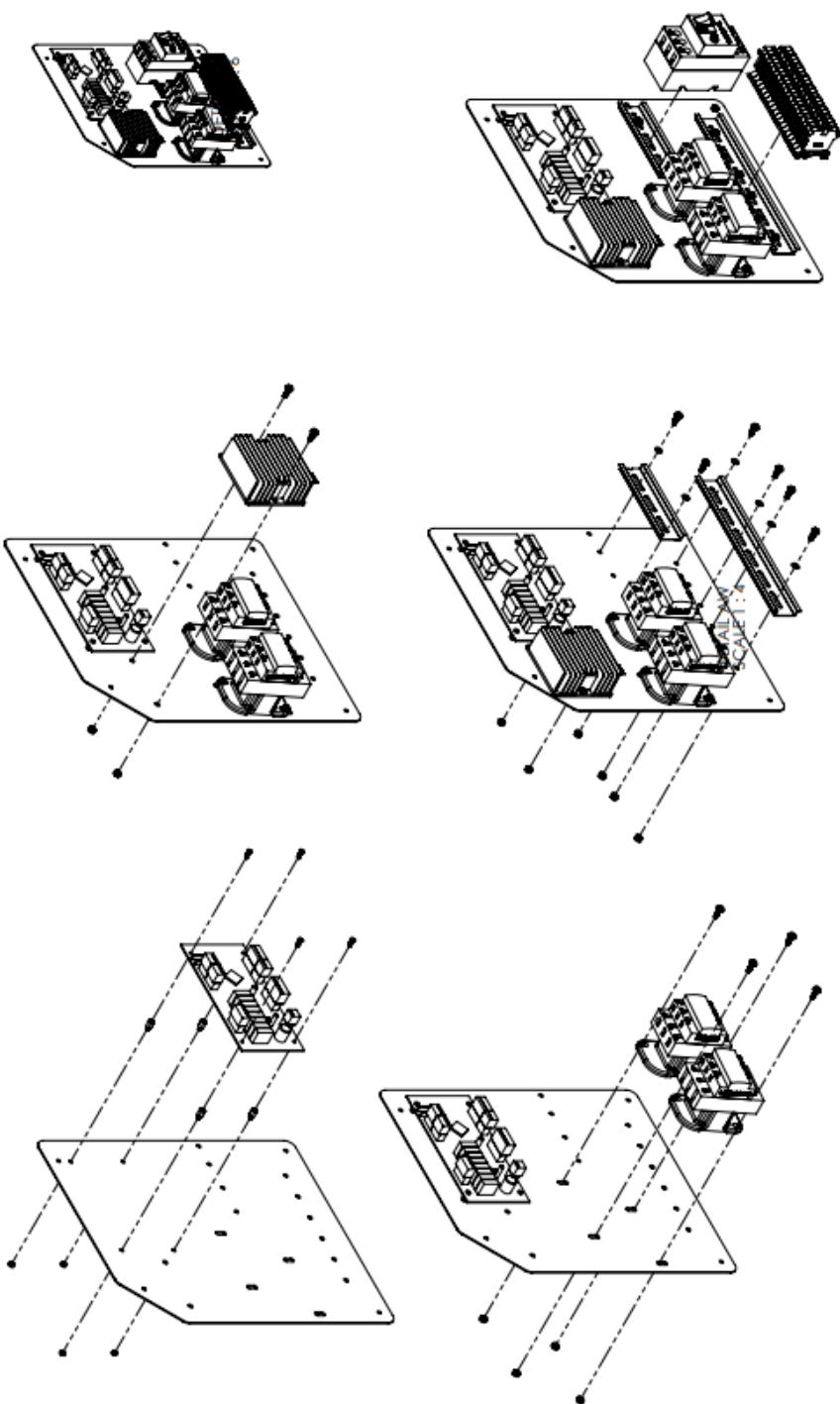








3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

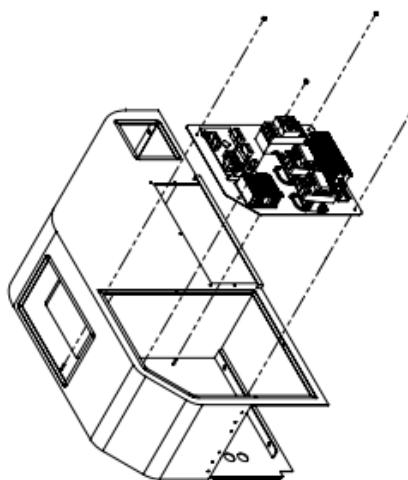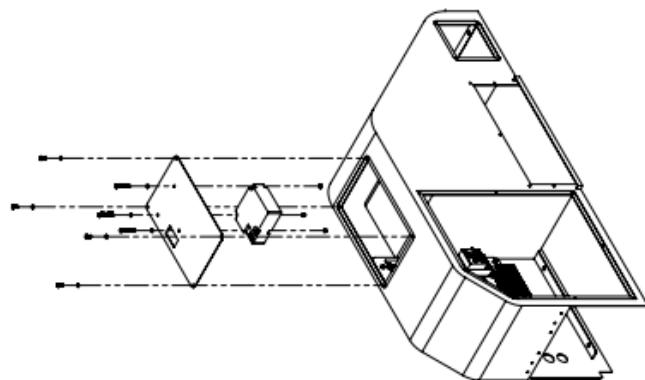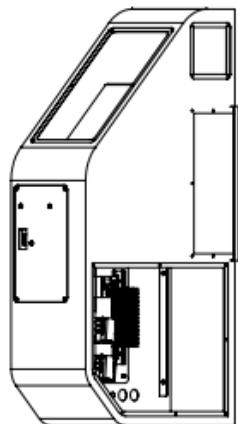








3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

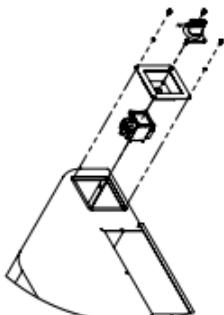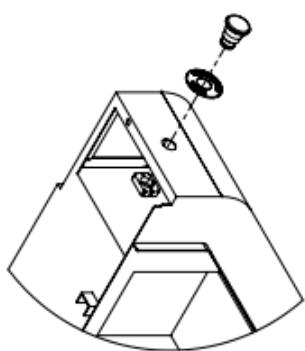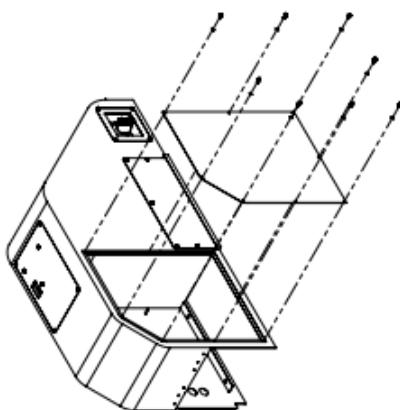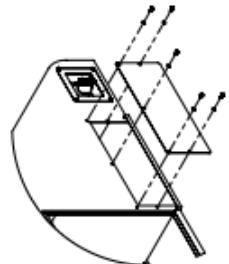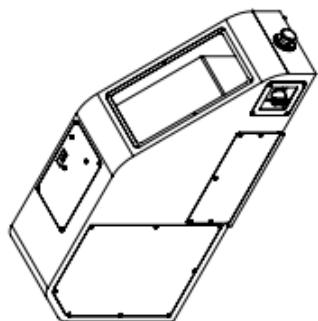



3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

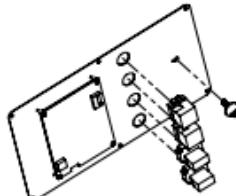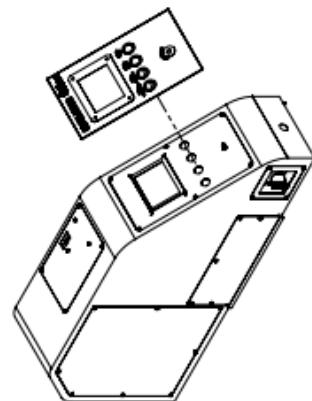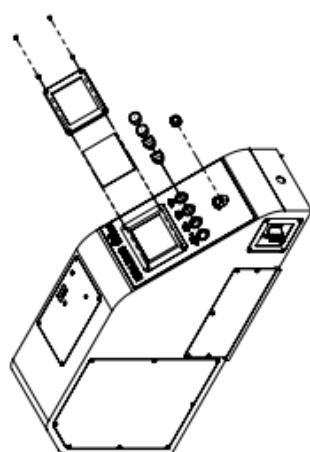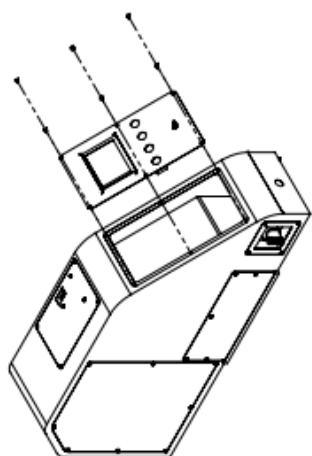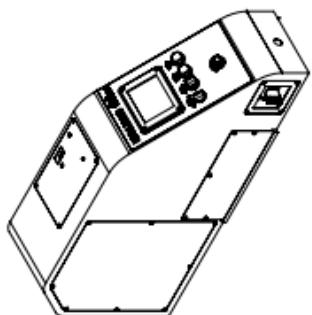



3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1-520-882-6598 Fax +1-520-882-6599 email: pace@metallographic.com Web: <https://www.metallographic.com>

