



# HV-1000Z Microhardness Tester



## INSTRUCTION MANUAL

3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: http://www.metallographic.com



Equipment Type:	Microhardness Tester (Vickers and Knoop)
Model:	<b>HV-1000Z</b>
Electrical Requirements:	110 Volts (single-phase)
Frequency:	50/60 Hz
Manual Revision Date:	April 2014

Please read this instruction manual carefully and follow all installation, operating and safety guidelines.



# HV-1000Z Microhardness Tester

## INSTRUCTION MANUAL

3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: http://www.metallographic.com

### Contents

	PAGE
<b>Warranty</b>	ii
<b>1.0 Product Description</b>	1
<b>2.0 Unpacking, Shipping and Installation</b>	3
<b>3.0 Safety Guidelines</b>	14
<b>4.0 Operation</b>	15
<b>5.0 Maintenance</b>	19
<b>6.0 Trouble Shooting</b>	20
<b>7.0 Spare Parts</b>	21
<b>8.0 Microhardness Testing Basics</b>	27
<b>9.0 Appendix RS232 Connector</b>	28
<b>9.1 Appendix MHT measurement software (optional)</b>	35

Please read this instruction manual carefully and follow all installation, operating and safety guidelines.



# HV-1000Z Microhardness Tester



## INSTRUCTION MANUAL

3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: <http://www.metallographic.com>

### WARRANTY

#### Terms and Conditions applying to all PACE Technologies Products

##### **1. LIMITED WARRANTY AND DISCLAIMER:**

PACE Technologies Products are warranted for one year from the purchase date to be free from defects in material and workmanship under correct use, normal operating conditions, and proper application. PACE Technologies obligation under this warranty shall be limited to the repair or exchange, at PACE Technologies option, of any PACE Technologies Product or part which proves to be defective as provided herein. PACE Technologies reserves the right to either inspect the product at Buyer's location or require it to be returned to the factory for inspection. Buyer is responsible for freight to and from factory on all warranty claims. The above warranty does not extend to goods damaged or subjected to accident, abuse or misuse after release from PACE Technologies warehouse, nor goods altered or repaired by anyone other than specifically authorized PACE Technologies representatives. PACE Technologies shall not in any way be responsible for the consequences of any alteration, modification or misuse unless previously approved in writing by an officer of PACE Technologies.

PACE TECHNOLOGIES MAKES NO EXPRESS WARRANTIES OTHER THAN THOSE WHICH ARE SPECIFICALLY DESCRIBED HEREIN. Any description of the goods sold hereunder, including any reference to Buyer's specifications and any description in catalogs, circulars and other written material published by PACE Technologies, is the sole purpose of identifying such goods and shall not create an express warranty that the goods shall conform to such description.

THIS WARRANTY IS EXPRESSLY IN LIEU OF ALL OTHER WARRANTIES, EXPRESSED OR IMPLIED. THERE ARE NO IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR PARTICULAR PURPOSE. THIS WARRANTY STATES PACE TECHNOLOGIES ENTIRE AND EXCLUSIVE LIABILITY AND BUYER'S EXCLUSIVE REMEDY FOR ANY CLAIM FOR DAMAGES IN CONNECTIONS WITH PACE TECHNOLOGIES PRODUCTS. PACE TECHNOLOGIES WILL IN NO EVENT BE LIABLE FOR INCIDENTAL OR CONSEQUENTIAL DAMAGES WHATSOEVER, NOR FOR ANY SUM IN EXCESS OF THE PURCHASE PRICE.

##### **2. LIABILITY CAP:**

PACE Technologies maximum aggregate liability for loss and damage arising under, resulting from or in connection with the supply or use of the Equipment and Consumables provided under this purchase, or from the performance or breach of any obligation (s) imposed hereunder, whether such liability arises from any one or more claims or actions for breach of contract, tort, (including negligence), delayed completion, warranty, indemnity, strict liability or otherwise, unless otherwise limited by the terms hereof, shall be limited to one hundred percent (100%) of the purchase price.

##### **3. DELIVERY:**

Customer assumes and shall bear the risk of all loss or damage to the Products from every cause whatsoever, whether or not insured, and title to such Products shall pass to Customer upon PACE Technologies delivery of the Products to the common carrier of Pace Technologies choice, or the carrier specified in writing by Customer, for shipment to Customer. Any claims for breakage, loss, delay, or damage shall be made to the carrier by the Customer and Pace Technologies will render customer reasonable assistance in prosecuting such claims.



## HV-1000Z Microhardness Tester

### INSTRUCTION MANUAL

3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: <http://www.metallographic.com>

#### **4. ACCEPTANCE:**

Customer shall inspect the Products promptly upon receipt of delivery. Unless customer objects in writing within thirty (30) business days thereafter, customer shall be deemed to have accepted the Products. All claims for damages, errors, or shortage in Products delivered shall be made by Customer in writing within such five (5) business day period. Failure to make any claim timely shall constitute acceptance of the Products.

#### **5. PAYMENT:**

Customer agrees to provide timely payment for the Products in accordance with the terms of payment set forth on the reverse side hereof or in any proposal submitted herewith. If any payment is not paid on or before its due date, Customer shall pay interest on such late payment from the due date until paid at the lesser of 12% per annum or the maximum rate allowed by law.

#### **6. DEFAULT:**

If Buyer is in default (including, but not limited to, the failure by Buyer to pay all amounts due and payable to Seller) under the work or purchase order or any other agreement between Buyer and Seller, Buyer's rights under the warranty shall be suspended during any period of such default and the original warranty period will not be extended beyond its original expiration date despite such suspension of warranty rights.

#### **7. MISCELLANEOUS PROVISIONS:**

This agreement has been made in and shall be governed by the laws of the State of Arizona. These terms and conditions and the description of the Products on the reverse side hereof or in any proposal submitted herewith constitute the entire agreement and understanding of the parties with respect to this sale and supersede all prior and contemporaneous agreements or understandings, inducements or representations, expressed or implied, written or oral, between the parties with respect hereto. Any term or provision of this Agreement may be amended, and any observance of any term of this Agreement may be waived, only by a writing signed by the party to be bounds. The waiver by a party of any breach shall not be deemed to constitute a waiver of any other breach. Should suit be brought on this Agreement, the prevailing party shall be entitled to recover its reasonable attorneys' fees and other costs of suit including costs and attorneys' fees incurred on appeal or in collection of any judgment.

## 1.0 Product Description

### 1.1 General Description



The HV-1000Z is a Microhardness tester for testing metallographic specimens.

The HV-1000Z is a variable load Microhardness tester with a load range of 10 grams to 1000 grams and has a auto turret for switching between the indenter and the objectives.



# HV-1000Z Microhardness Tester



## INSTRUCTION MANUAL

3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: <http://www.metallographic.com>

### 1.2 Technical Specifications

Electrical specifications:	220V single-phase (50/60 Hz)
Test forces:	0.098N(10gf) 0.245N(25gf) 0.49N(50gf) 0.9807N(100gf) 1.961N(200gf) 2.942N(300gf) 4.903N(500gf) 9.807N(1000gf)
Magnification:	100X, 400X
Dwell time of test force:	1-60 seconds
Max. height of specimen:	65 mm (2.6-inches)
Loading:	Automatic loading and unloading
Weight:	Approx. 55 lbs (25 kg)
Dimensions (WxHxD):	Approx. 16.5" x 7" x 18.5" (420 mm x 180 mm x 470 mm)
Working temperature:	70° - 85°F (23 - 28°C)



3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: http://www.metallographic.com

## 2.0 Unpacking, Shipping and Installation

### 2.1 Unpacking

Unit is delivered in a box. Unpack and check for completeness of parts.

Measures WxHxD: 21"x17"x21" (530x430x530 mm)

Weight: Approximately 75 lbs (35 kg)

### 2.2 Shipping

When moving box, lift with handles.



**! Caution:** Heavy equipment. Take care to avoid bodily injury.


## 2.3 Installation

! Install unit carefully! Improper installation voids warranty.

The **HV-1000Z** should be placed on a flat stable vibration free surface.

### 2.3.1 Remove shipping securing string and tape

Remove shipping ties

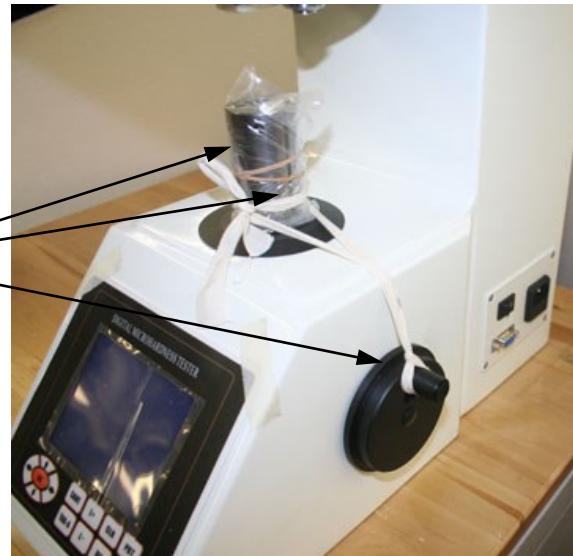


Remove indenter securing tape








3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: http://www.metallographic.com

### 2.3.2 Install leveling feet

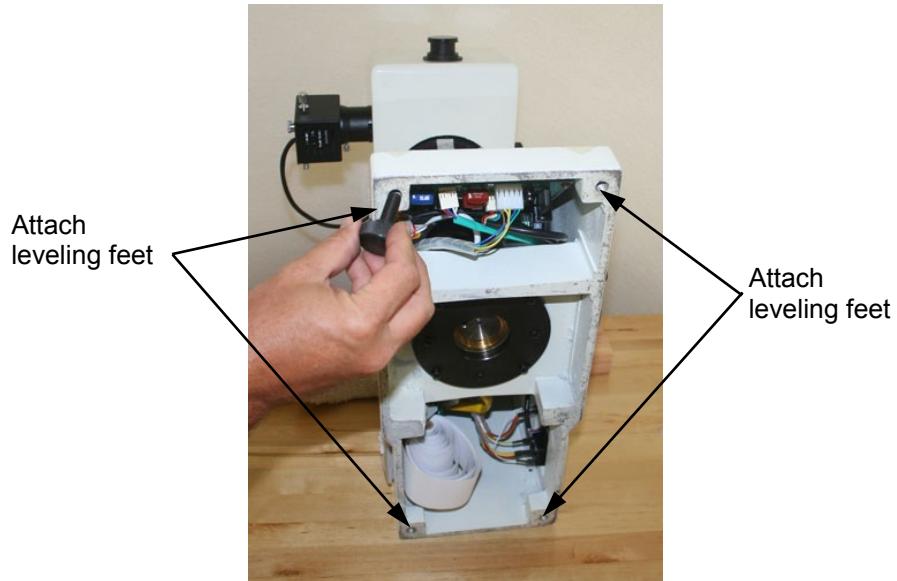


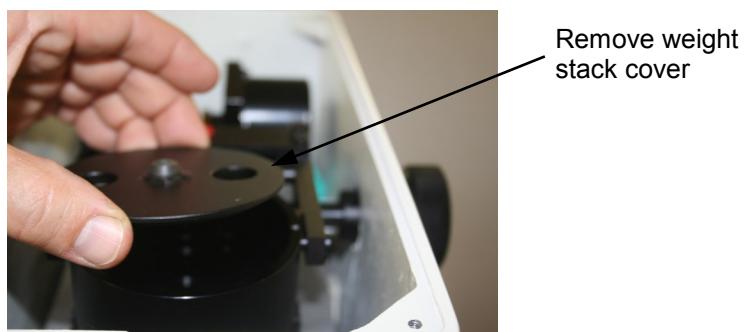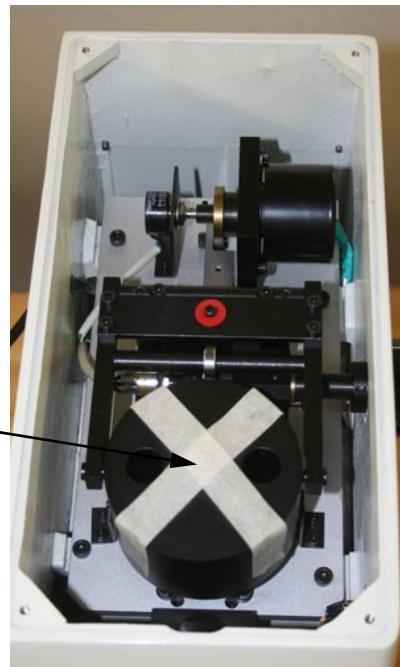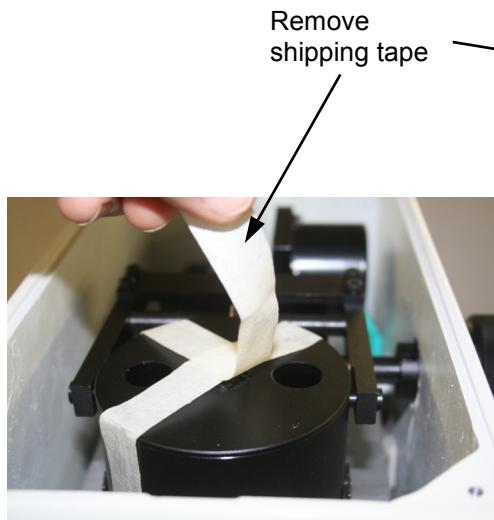**2.3.3 Remove shipping securing tape from weight holder**




3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: http://www.metallographic.com


### 2.3.4 Remove shipping securing screws for Loading lever

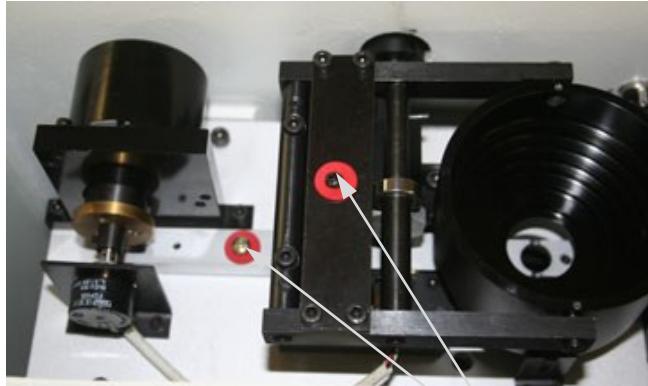


Shipping securing screws

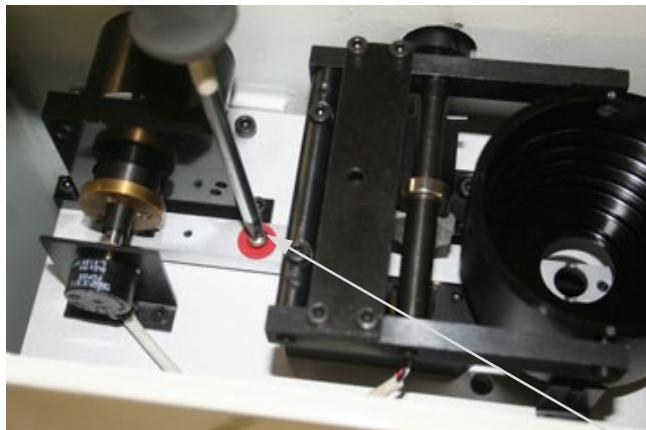


Remove shipping securing screws

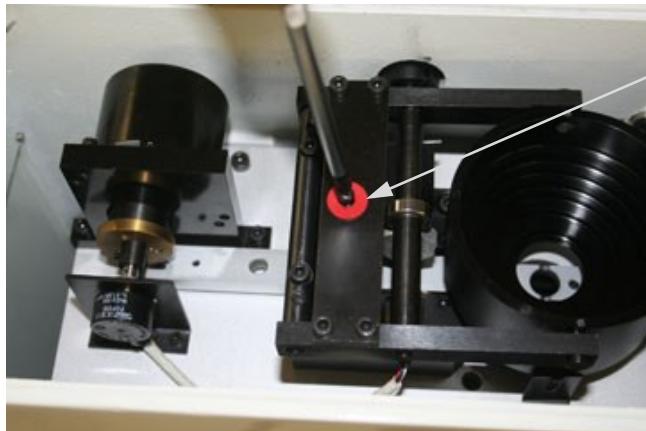

### 2.3.5 Weight load stack

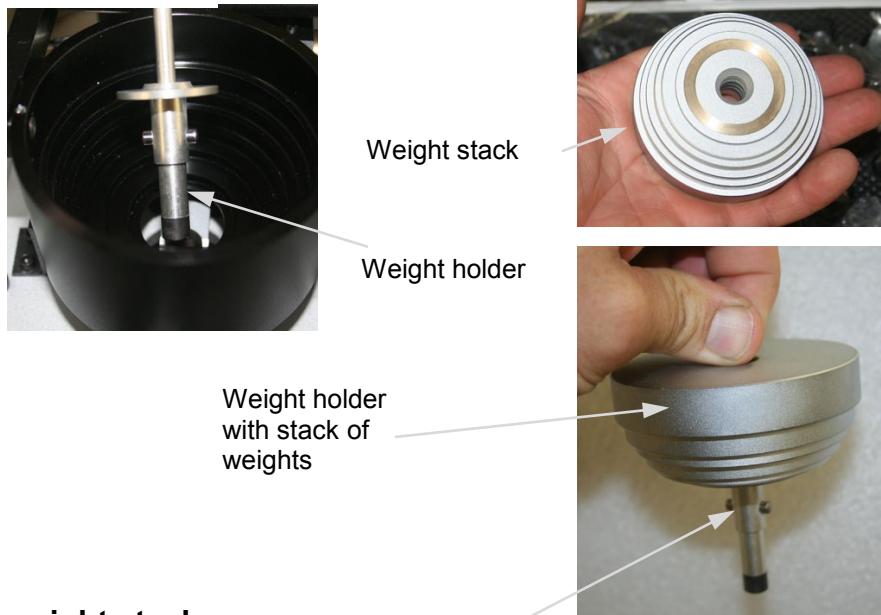


#### 2.3.5.1 Installing weight stack

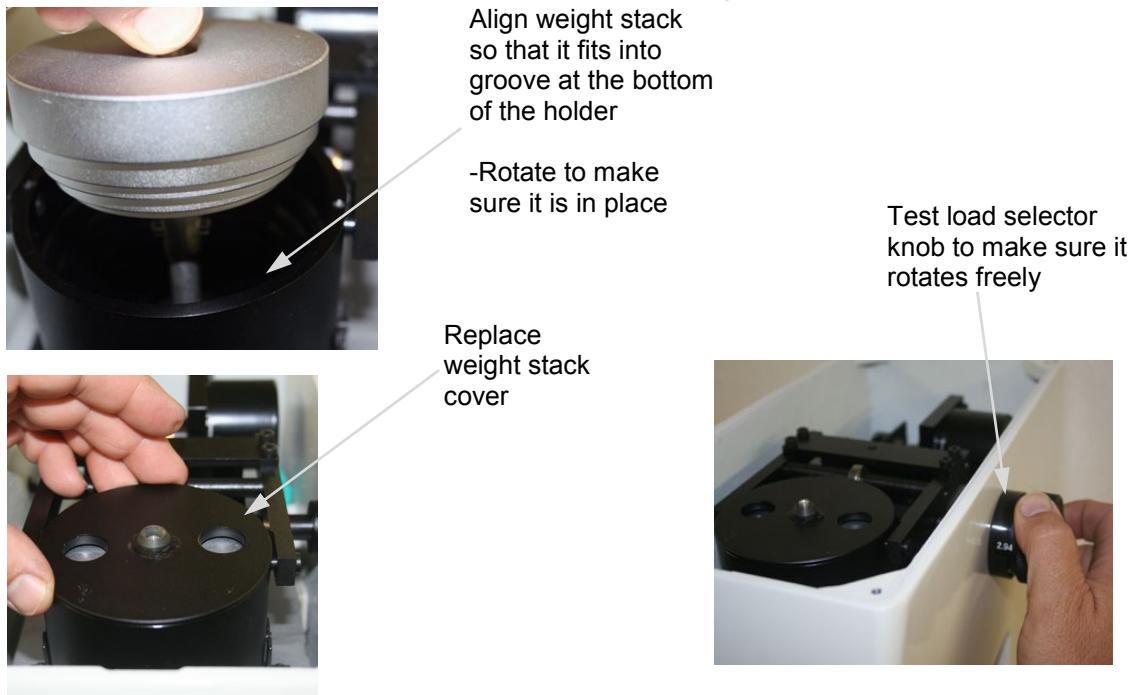



3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: http://www.metallographic.com

### 2.3.6 Installing filar eyepiece

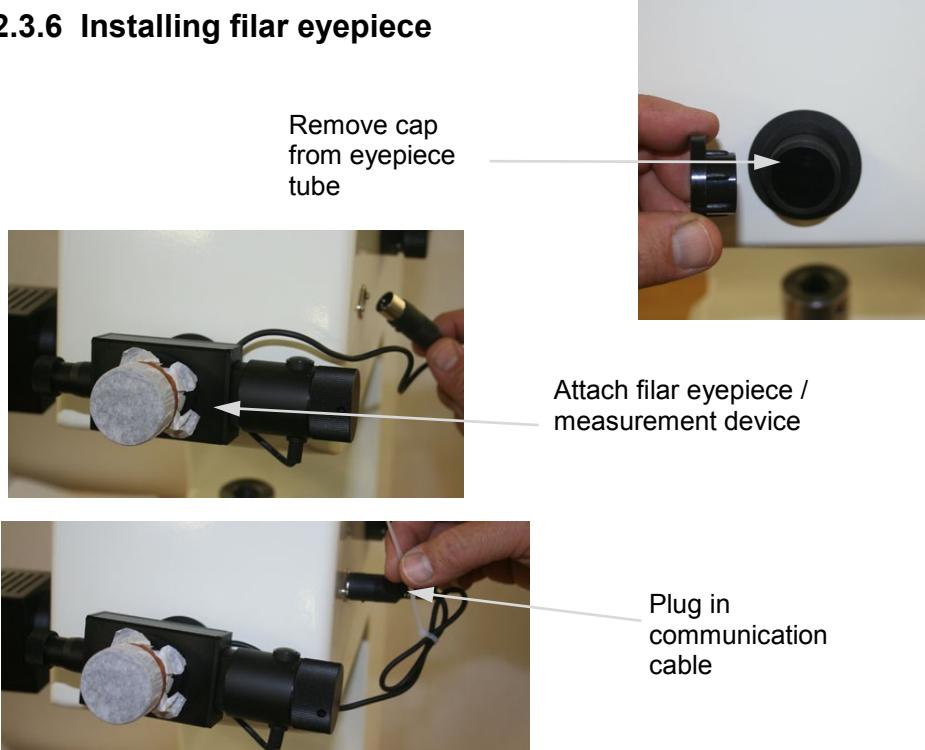

### 2.3.7 Installing camera



-Insert c-mount adapter into port and attach camera.

To center image in camera the following adjustment may need to be made by moving lever shown. First remove filar and unscrew lens, lever is inside.

**2.3.8 Installing stage**

Lower stage height and attach stage



Tighten set screw

**2.3.9 Leveling unit**


Place bubble level on stage and adjust feet height to level the unit



3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: http://www.metallographic.com

### 2.3.10 Installing digital camera (optional)


- Remove MHT top panel.
- Remove plastic cover on lens splitter
- Screw in camera adapter
- Replace MHT cover plate
- Screw on camera

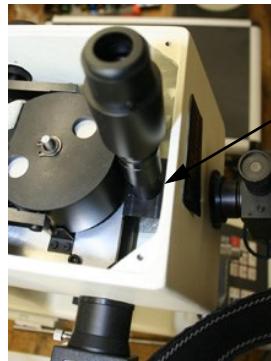

Remove plastic cover and screw in camera adapter



Remove plastic cover and screw in camera adapter

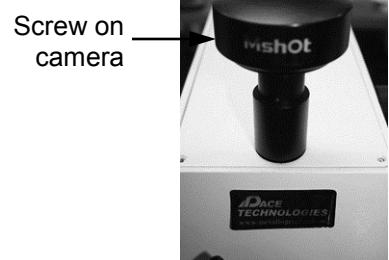

Screw on camera

#### 2.3.10.1 Aligning light path for camera

- If the light path for the camera is off center it can be adjusted by adjusting the angle of the prism
- Remove filar
- Unscrew eyepiece lens
- With a screwdriver move lever left or right to adjust



Remove filar and unscrew eyepiece lens



Camera / eyepiece splitter prism adjustment (prism shown out of unit above)

#### 2.3.10.2 Focusing Digital camera to filar eyepieces


- Focus through filar eyepieces
- Adjust focal length of camera adapter to match focus



Loosen locking nut and adjust focal length, tighten to lock in place

### 2.3.11 Replacing Indenter

#### 2.3.11.1 Remove Indenter

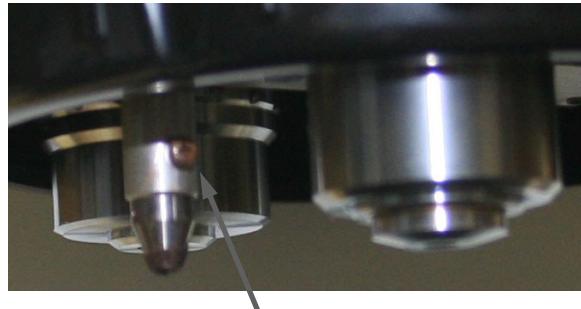


Loosen screw and  
remove indenter

#### 2.3.11.2 Install Indenter

With turret in the indenting position  
(indenter forward)

-Align red dot on indenter so it is facing forward for proper alignment  
(note for Knoop indenters the some additional adjustment may be required to align the indenter)

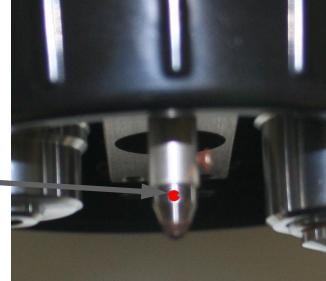

After loading indenter - run a sample indent to properly set the indenter

### 2.3.11.3 Aligning the indenter with the objectives

To align the indenter with the objective there are three set screws on the objective.

-remove the turret plastic cover (three screws)  
-Adjust the set screws (note if the indenter is too far off it may need to be rotated in 90 degree increments for Vickers or 180 degrees for Knoop)



Remove three screws to remove plastic cover



Adjustment set screws

Adjustment set screws



3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: http://www.metallographic.com

### 2.3.12 Calibrating Measurements

The unit has been calibrated before shipment, however, if it is necessary to recalibrate the unit the magnification of the image (optics) will need to be changed.

Prior to adjusting the optics it is recommended that the following be checked first:

1. Verify that the filar is pushed all the way into the light tube
2. Verify that the light tube is not loose (turn clockwise to tighten)
3. Verify that the unit level (use bubble leveler)
4. Verify that the load stack is hanging free and weight stack arm crossbar is positioned in the loading arm slot
5. Check that the indenter is not binding in its holder.
6. Focus the eyepiece so the filar lines are in sharp focus.

#### 2.3.12.1 Calibrating Procedure

Step 1. remove filar from light tube



Step 2. remove light tube cover plate



Step 3. Unscrew light tube

Step 4. Adjust sleeve on light tube.

Example:

To increase the measured hardness turn the locking ring clockwise so the tube threads deeper into the prism (shorter tube length). As an approximation 1/2 turn of the ring changes the measured value by about 2% (1000 gram load, 1/2 turn increased hardness from 405-415)



Step 5. reverse procedure to reinstall and retest hardness.  
Repeat procedure until calibrated.



# HV-1000Z Microhardness Tester

## INSTRUCTION MANUAL

3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: http://www.metallographic.com

### 3.0 Safety Guidelines

#### 3.1 Warning Sign

- ! This sign points to special safety features on the machine.

#### 3.2 Safety Precautions

- ! Careful attention to this instruction manual and the recommended safety guidelines is essential for the safe operation of the **HV-1000Z**.
- ! Proper operator training is required for operation of the **HV-1000Z**. Any unauthorized mechanical and electrical change, as well as improper operation, voids all warranty claims. All service issues need to be reported to the manufacturer / supplier.

- ! Operate unit as specified in this manual.
- ! Disconnect power before opening unit.
- ! Lower stage to avoid damaging indenter or lens when not in use.
- ! Cover unit with dust cover when not in use to eliminate dust contamination.

#### 3.3 Emergency Statement

Always follow proper operational guidelines and avoid contact with lubricants and abrasives.



3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: http://www.metallographic.com

### 4.0 Operation

#### STEPS FOR HV-1000Z OPERATION IN BRIEF

(Please see explanation of important points on next pages before proceeding)

1. Mount the sample so it is flat on the stage. For non-parallel mounts a leveling vice is recommended (optional)
2. Turn on the HV-1000Z (allow the auto turret to calibrate itself before beginning any operation). Focus on an edge of the sample and then move to the location for the indentation. Adjust intensity of light (L- and L+), set Filar to ZERO (press CLR to zero), select the Load and enter the Dwell Time
3. Press Start to make Indentation. Note: The auto turret will automatically rotate the indenter into place, apply the load and then rotate back to the 40X objective.
4. Measure D1 / D2 (for Vickers) (PRESS BUTTON ON FILAR TO ACCEPT), rotate the Filar micrometer counter-clockwise, measure D2 (PRESS BUTTON ON FILAR TO ACCEPT), and then rotate to original position clockwise). For Knoop measure only the major axis, D1 is measured (PRESS BUTTON ON FILAR TO ACCEPT). NOTE THE FIRST INDENT IS NOT RECORDED INTO MEMORY AS IT IS COMMON PRACTICE FOR THE FIRST INDENT TO BE USED TO SET THE INDENTOR.
5. Note the HV or HK value
6. If the measurement is bad, press (DEL) and re-measure indent.
7. If more impressions need to be taken, move stage to relevant location and go to step 3
8. PRT to print results or to send to Hyper Terminal
9. Lower stage and remove the sample and/or leveling vice. Switch the instrument off.



Measurement button





# HV-1000Z Microhardness Tester

## INSTRUCTION MANUAL

3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: http://www.metallographic.com

### PRECAUTIONS

Clamp the sample properly on the self-leveling vice. This will make the sample exactly perpendicular to the indenter.

**Important :** This is very important as a tilted sample can damage the indenter. It is also very important to have a flat specimen. A rounded specimen will have varying height as it moves and can touch and damage the indenter.

The stage should be at its lowest position before placing the leveling vice on the stage. The stage can be lowered to its lowest position by rotating the focusing handle in the counter-clockwise direction.

**Important :** Take utmost care to insure that the sample or the leveling vice does not touch the indenter or the optics as this may damage the indenter

Switch on the HV-1000Z and focus the specimen under 10x objective ( Care should be taken to see that the indenter is not touching any part of the sample or the leveling vice when you index the indenter / objective )


**Important :** When raising the stage to get the sample in focus, always look at the sample as it is raised from the side ( NOT from the eyepiece ) and stop raising it when there is a gap of approximately 0.1 mm between the sample and the objective

Index the indenter over the sample (Care should be taken to insure that the indenter is not touching any part of the sample or the self leveling vice when you index the indenter / objective)

Next swing the 40x objective into place and make sure that the indenter again does not touch anything when being indexed. Focus if necessary

**Important :** When raising the stage to get the sample in focus, always look at the sample as it is raised from the side ( NOT from the eyepiece ) and stop raising it when there is a gap of approximately 0.1 mm between the sample and the objective tip. Then view through the eyepiece and lower the stage to do the fine focusing

#### 4.1 Control Panel

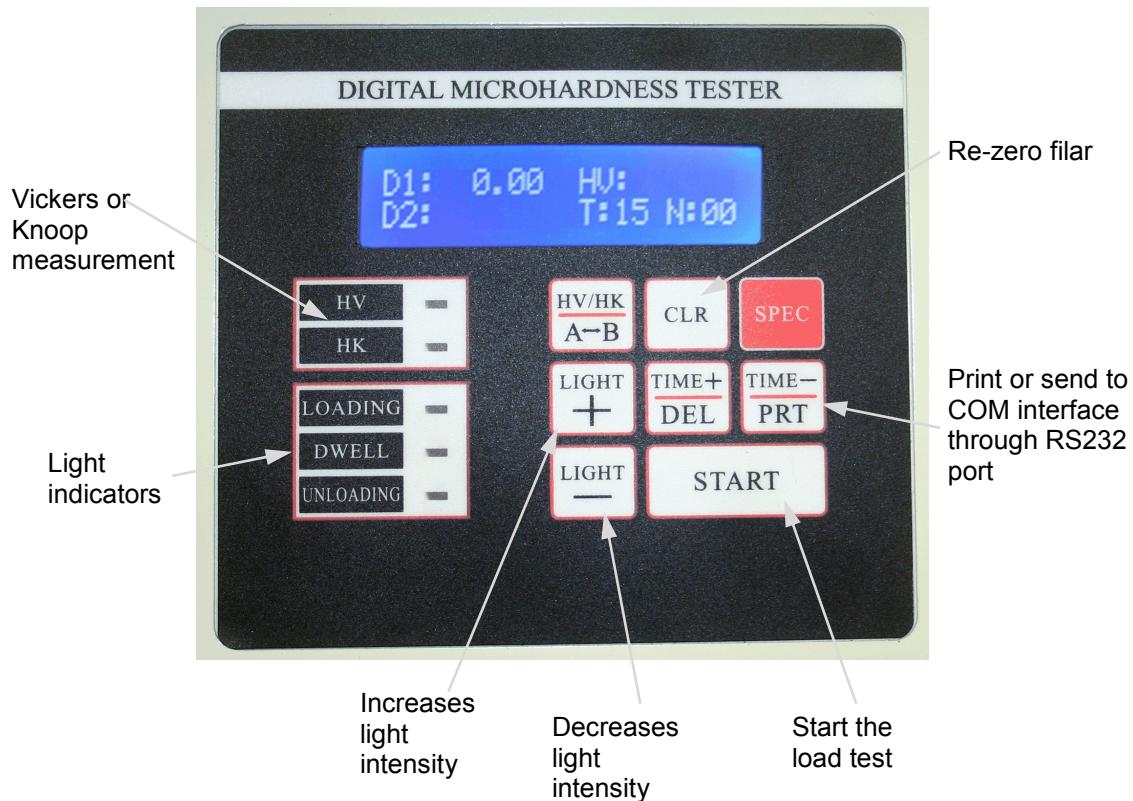

**Setting DATE:** To set the date, immediately after the date shows on the display after the WELCOME, press the TIME+ or TIME- key. Set the year and press SPEC TIME+ , set the month, press SPEC TIME+ to adjust the day. To exit date field press SPEC.



## HV-1000Z Microhardness Tester

### INSTRUCTION MANUAL

3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: http://www.metallographic.com

CLR	The CLR key is used to reset the filar to zero
PRT	<p>The PRT key is used to print or connect the interface; the key serves to print the data of the present test or send them to the COM interface of a PC through RS232.</p> <p>Note: see Appendix for setting up Hyper Terminal to PC</p>
L+	The L+ key is used to increase the intensity of the light. A continuous tone indicates the maximum intensity has been reached.
L-	The L- key is used to reduce the intensity of the light. A continuous tone indicates that the intensity of the light is at a minimum.
START	The START key is used to begin the test
DEL	Delete's entry after the measurement has been made. Note you must enter a value with the button on the filar before you DEL. Otherwise you will not be able to continue.



# HV-1000Z Microhardness Tester



## INSTRUCTION MANUAL

3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: <http://www.metallographic.com>

### 5.0 Maintenance

#### 5.1 Introduction

The **HV-1000Z** requires very minimal maintenance. However, to increase the life of the Microhardness tester, it is suggested that the unit be covered when not in use.

#### 5.2 Cleaning outside cabinet

The cabinet should be cleaned occasionally with a moistened cloth. Do not use any chemicals or cleaning abrasives.


## 6.0 Trouble Shooting

**More extensive trouble shooting, repair guides, video's, parts list are provided online at [www.metalgraphic.com](http://www.metalgraphic.com) or**

<http://www.metalgraphic.com/PACE-service/MHT-service.html>

Problem	Cause	Solution
No power or function	a. Unit is disconnected from main electrical power supply b. Main power switch is off c. Blown fuse	a. Verify electrical source and connection. b. Turn on main power switch. c. Replace fuse
Indenter does not move	a. Unit needs to be reset b. Internal damage	a. Turn off power and restart b. Call service center.
Chipped indent	Broken indenter	Replace indenter
Indenter not in alignment with objective	Objective not centered on indenter	-Align objective -If grossly out of alignment the indenter may need to be rotated in 90 degree increments for Vickers and 180 degrees for Knoop indenters
Indent not symmetric	Indenter may need to be rotated	-Rotate indenter 90 degree increments for Vickers and 180 degrees for Knoop indenters
Indenter not indenting properly	Indenter balancing pin out of position (see image below)	-Place dislodged part back into position (see below)

Align this into position

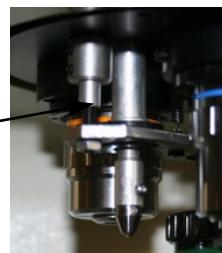



3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: http://www.metallographic.com

## 7.0 Spare Parts


Part no.	Description	Image
	<b>Electrical components</b>	
<b>MNT-1050</b>	HV-1000Z control board	
<b>MHV2-E-CB</b>	MHV-2000 control board	
<b>MHT-010</b>	MHT driver control board	
<b>MHT-011</b>	MHT driver control input board	
<b>MHT-012</b>	MHT indenter motor	
<b>MHT-013</b>	MHT turret motor	
<b>MHT-014</b>	MHT 110/220 volt transformer	



## HV-1000Z Microhardness Tester

### INSTRUCTION MANUAL

3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: http://www.metallographic.com

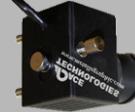Part no.	Description	Image
	<b>Electrical components</b>	
<b>MHT-E-LH</b>	MHT lamp housing	
<b>MHT-Bulb</b>	MHT light bulbs (each)	
<b>811-881</b>	PC-RS232 software cable	
<b>811-882</b>	PC RS-232 to USB software cable	
<b>MHT-E-PRI</b>	MHT Printer-1	
<b>MHT-PRI2</b>	MHT Printer-2	
<b>CORD-110</b>	110V USA power cord	
<b>CORD-220R</b>	220V round prong power cord	
<b>CORD-220F</b>	220V flat prong power cord	



3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: http://www.metallographic.com

Part no.	Description	Image
<b>Optical components</b>		
<b>MHT-001</b>	MHT filar eyepieces	
<b>823-311</b>	MHT-10X objective	
<b>823-331</b>	MHT-40X objective	

Part no.	Description	Image
<b>Camera components</b>		
<b>823-361</b>	MHT 1X camera adapter (c-mount)	
<b>MP50</b>	5 MP digital camera	
<b>Microhardness-Pro</b>	Microhardness Pro software for measuring indent size	



## HV-1000Z Microhardness Tester

### INSTRUCTION MANUAL

3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: http://www.metallographic.com

Part no.	Description	Image
	<b>Camera components</b>	
MHT-003	MHT X-Y table w/o micrometers	
MHT-ENG	MHT micrometers (English units) (2/set)	
MHT-MET	MHT micrometers (Metric units) (2/set)	
MHT-M-VI	Vickers Indenter	
MHT-M-KI	Knoop Indenter	
MHT-M-IRS	Indenter retention screw	
MHV-S	Precision screw feed height adjustment unit for microhardness testers	
MHT-M-HAW	Height adjustment wheel	
MHT-FEET	MHT feet (4/pkg)	
823-902	Dust cover for MHT	

Please read this instruction manual carefully and follow all installation, operating and safety guidelines. 24



3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: http://www.metallographic.com

Part no.	Description	Image
	<b>Camera components</b>	
<b>811-100-1</b>	Printer paper rolls for MHT-1 printer (2.25-inch width)	
<b>811-100-2</b>	Printer paper rolls for MHT-2 printer (1.75-inch width)	
<b>MHT-P-TPC</b>	Turret plastic cover	
<b>MHT-2000-W</b>	MHV-2000 additional 1000 gram weight	
<b>MHT-W</b>	10 to 1000 gram weight stack for MHT testers	
<b>MHT-ROD</b>	Weight stack holder for MHT testers	



## HV-1000Z Microhardness Tester

### INSTRUCTION MANUAL

3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: http://www.metallographic.com

Part no.	Description	Image
	<b>Fixturing vises</b>	
<b>823-741</b>	MHT sample leveler	
<b>MHT-BUB</b>	Bubble Level	
<b>MHT-004</b>	MHT Fork-shaped test table	
<b>MHT-005</b>	MHT Thin specimen test table	
<b>MHT-006</b>	MHT Fine wire test table	




3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: <http://www.metallographic.com>

## 8.0 Microhardness Testing Basic

Microhardness testers are both mechanical and optical measuring tools. The indent is produced by applying a known load to the specimen and then measuring the size of the appropriate diagonals either optically or with image analysis software.


Microhardness is primarily determined with either a Knoop or Vickers indenter under test loads in the range of 1 to 2000 gram-force. Microhardness is used to measure the hardness of specific phases, small particles, and for brittle materials. Figure 6-1 shows the Knoop microhardness indents for a heat treated steel specimen.

Knoop hardness (HK) number is based on the size of the indent that a rhombic-based, pyramidal diamond indenter produces under a known applied load. The HK number is calculated by dividing the applied load (kilogram-force) by the projected area of the indentation (square millimeters).

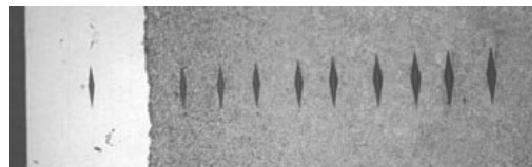

**Figure 6-1** Knoop hardness indents in a hardened steel.

The Vickers hardness (HV) number is obtained by dividing the applied load in kilogram-force by the surface area of the indentation. The area of the indentation produced from the Vickers square-based pyramidal diamond is determined by the mean distance between the two diagonals of the indentation (Figure 6-2).

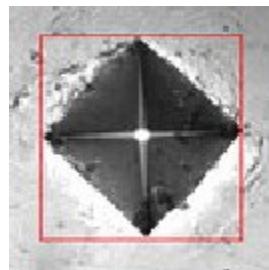

**Figure 6-2** Vicker indent.

## 9.0 Appendix RS232 Connector


For Window XP HyperTerminal is already installed. If you are using WIN7 the program needs to be installed (contact us for CD with program).

Also if your computer has an RS232 port use this port instead of the USB port. If you do not have an RS232 port you will need to install a conversion driver (optional)

First determine which scenario fits your computer.

- A. Window XP with RS232 port (Step 4 only)
- B. Windows XP without RS232 port (need to install conversion programs - PL203\_Prolific\_DriverInstaller\_v10518) (Step 1,2 and 4)
- C. Window WIN7 without RS232 port (need to install conversion programs - PL203\_Prolific\_DriverInstaller\_v10518 and Hyperterm software) (Step 1,2,3)
- D. Window WIN7 with RS232 port (need to install Hyperterm) (Step 1 and 3).

STEP 1. Load Install Files to Desktop or other easy to find location and connect cable to tester and turn on

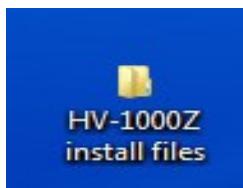


3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: http://www.metallographic.com


**STEP 2.** Computer does not have an RS232 port (need to adapt to USB port)

Use the following procedure:


- Open HV-1000Z Install File Folder

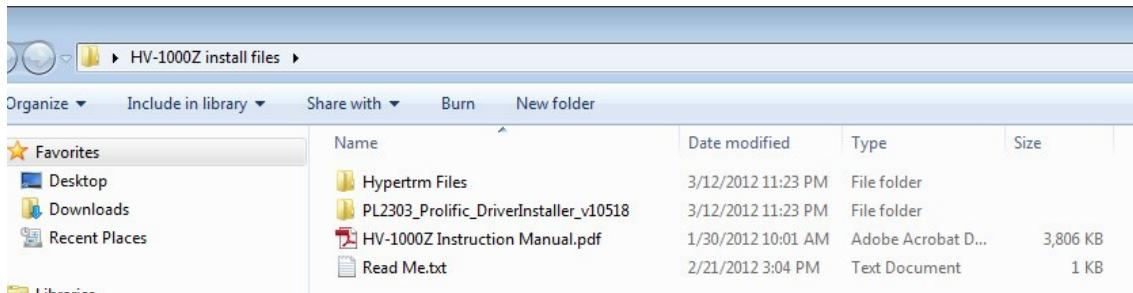


- Double click or open **PL203\_Prolific\_DriverInstaller\_v10518** file

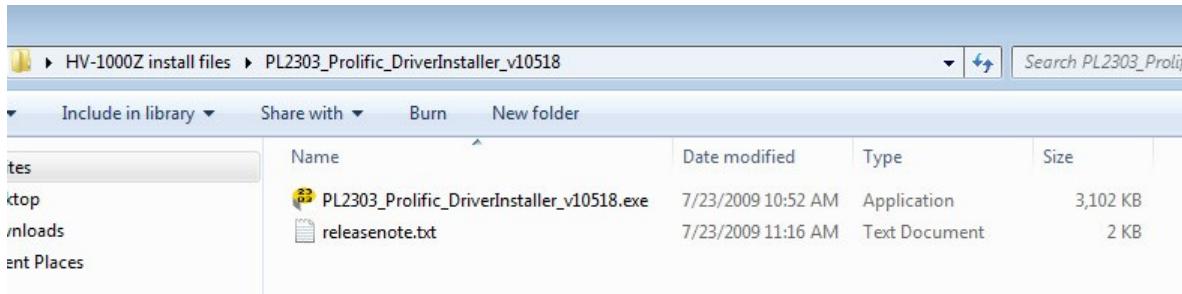

- Double click or run **PL203\_Prolific\_DriverInstaller\_v10518.exe**  
-Follow installation procedure (Click on NEXT)

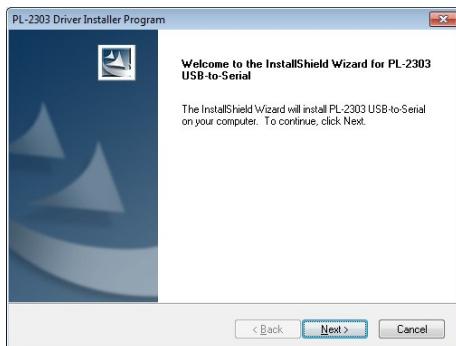


-After installation is complete (Click on FINISH) and restart computer.

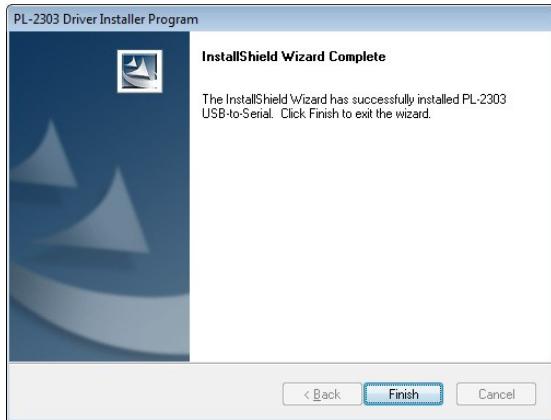

STEP 3. Computer uses WIN 7 and does not have Hyperterm installed.

- Double click to open HV-1000Z installation file

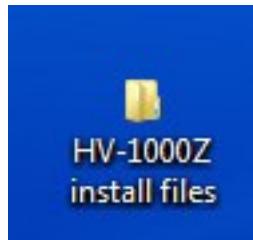


Name	Date modified	Type	Size
Hyperterm Files	3/12/2012 11:23 PM	File folder	
PL-2303	Date created: 3/12/2012 11:23 PM	File folder	
Hyperterm	2012-03-12 10:01 AM	Adobe Acrobat D...	3,806 KB
Readme.txt	Files: hyperterm.chm, hyperterm.dll, hyperterm.exe, ... 3/12/2012 3:04 PM	Text Document	1 KB




3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: http://www.metallographic.com


- Double click or open Hyperterm Files folder
- Double click or Run hypertrm.exe

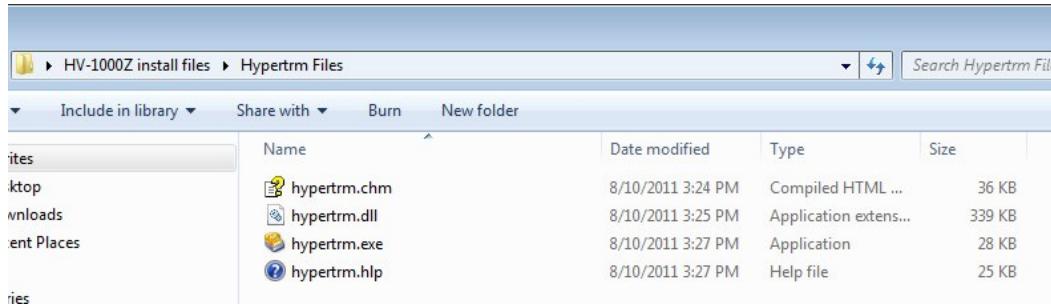

- Answer NO to question

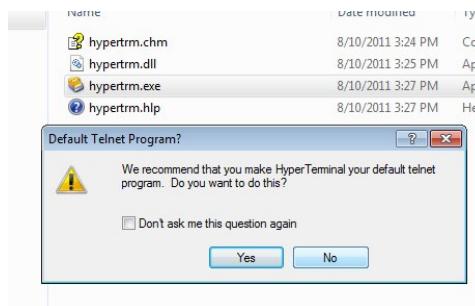


-Name: Input RS232. Click OK

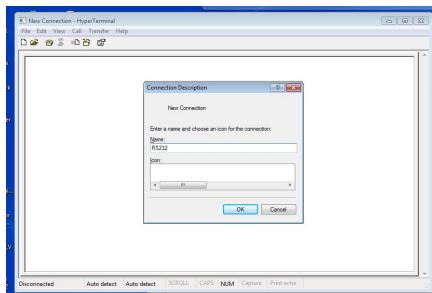


- If hardness tester was on and connect the correct COM port should be selected. Click OK. Note if connect using TCP/IP (Winsock) is displayed the computer was not rebooted after adding RS232 to USB conversion and needs to be rebooted.

-Change the Bits per seconds to 9600. Click OK

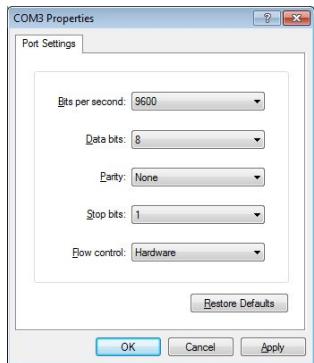

-HyperTerminal will now be available to accept the data by pressing PRT

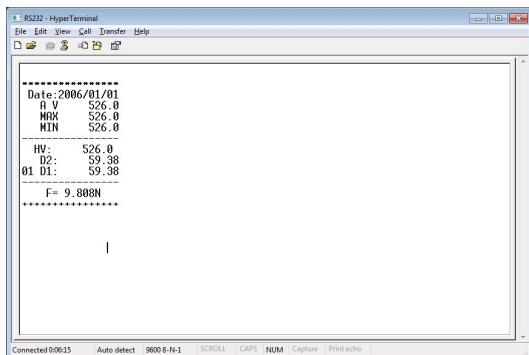

**-NOTE: FOR WIN7** if you turn off computer or close the HyperTerminal connection you will need to do the following to access the program the next time you use it (A bit inconvenient, however this is because Microsoft does not support this 3rd party software anymore).

1. Open Hyperterm Folder (on Desktop if saved there during installation)
2. Double click on Hyperterm.exe
  - Choose NO
  - Cancel or “X” out of window.
3. File > Open > Select RS232.ht file (if using same USB port as it was installed on)


Please read this instruction manual carefully and follow all installation, operating and safety guidelines.




3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: http://www.metallographic.com

#### STEP 4. To set-up with Windows XP

-Click → Start, ALL Programs → Accessories → Communications → Hypertrm

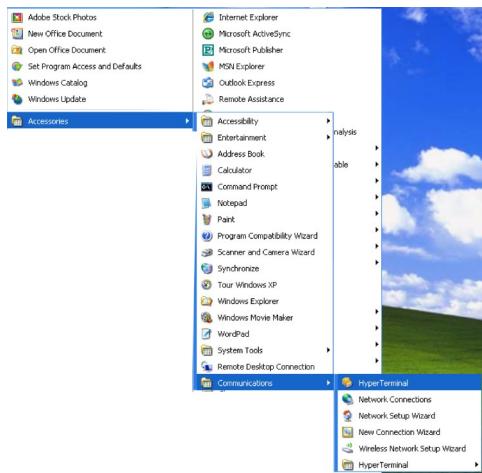


-In the new connection box, enter “RS232” in the Name field and click OK

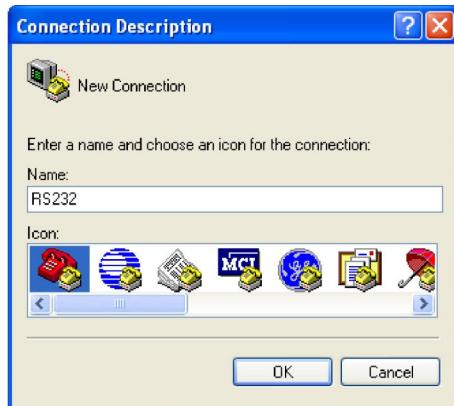

-In the Connect To dialogue box choose the appropriate COM port from the Connect using drop down box, Click OK



-In the Port Properties dialogue box, select 9600 for Baud rate, and leave the other parameters unchanged, Click Apply, OK and HyperTerminal will be connected to the tester.

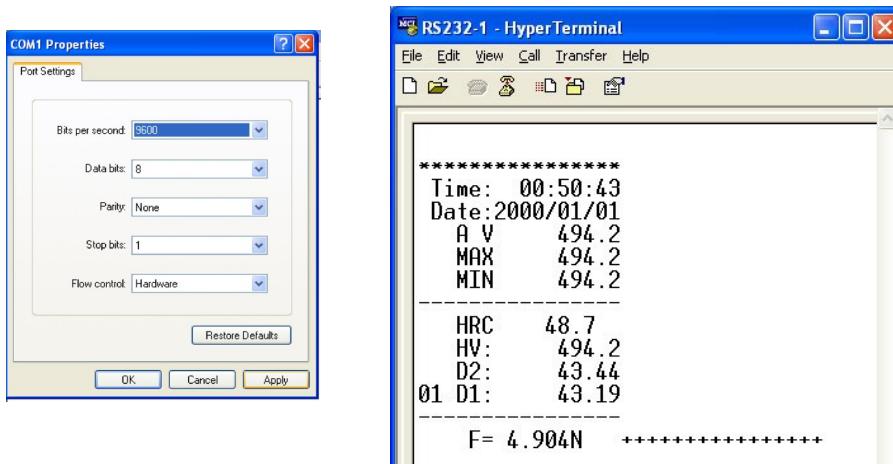


### 9.1 Appendix MHT measurement software (optional)

#### 9.1.1 General Description


**Microhardness Pro 4.2** is an easy to use image analysis software for measuring microhardness and for generating easy to distribute Excel reports. It is a single screen Windows based system.

The **Weld Check** software can handle both gray monochrome (8 bit) and color (24 bit) images. Multiple images of any size can be opened and displayed on the screen for analysis or comparison. The software supports the most common formats like BMP, JPEG, TIFF, PNG, GIF and PSD. Live images also can be observed and captured on the same platform. Since the system is in a Windows environment, graphs and charts displayed on the monitor can be quickly transferred into other Windows based programs such as, MS Word, MS Excel or any other commercial Windows based software for the purpose producing reports and presentations.


Main screen

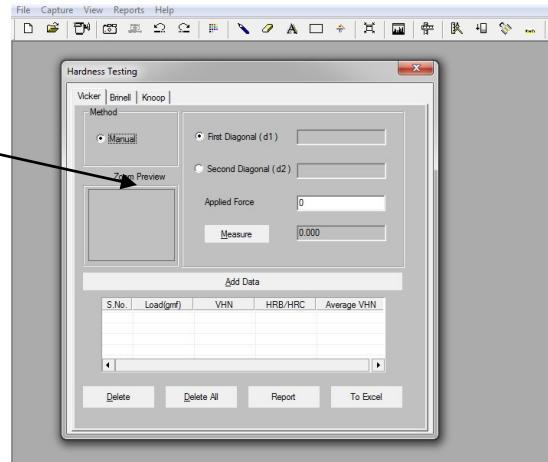

Capture screen

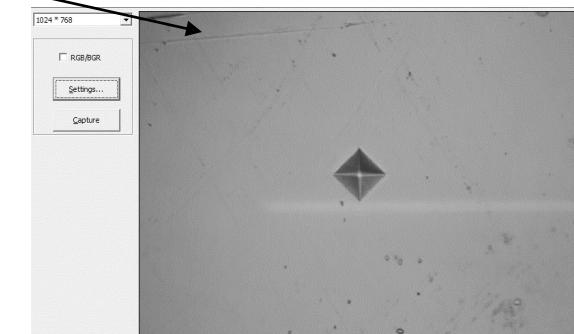

Vickers Indent

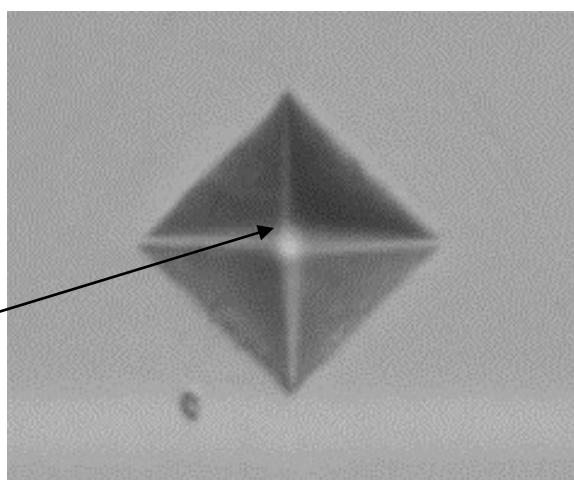


### 9.1.2 Installing Software

#### 9.1.2.1 Insert software CD and click on setup

Files Currently on the Disc (10)			
 data1	4/7/2010 4:12 AM	Cabinet File	753 KB
 data2	4/7/2010 4:12 AM	Cabinet File	79,362 KB
 engine32	8/4/2002 10:59 PM	Cabinet File	387 KB
 back final	4/21/2009 10:55 PM	BMP File	2,305 KB
 data1	4/7/2010 4:12 AM	HDR File	92 KB
 layout.bin	4/7/2010 4:12 AM	BIN File	1 KB
 setup.boot	4/7/2010 4:11 AM	BOOT File	327 KB
 setup	8/4/2002 10:59 PM	Application	96 KB
 setup	4/7/2010 4:11 AM	Configuration sett...	1 KB
 setup	10/20/2009 11:26 ...	InDesign XML Inte...	158 KB

#### 9.1.2.2 Follow prompts

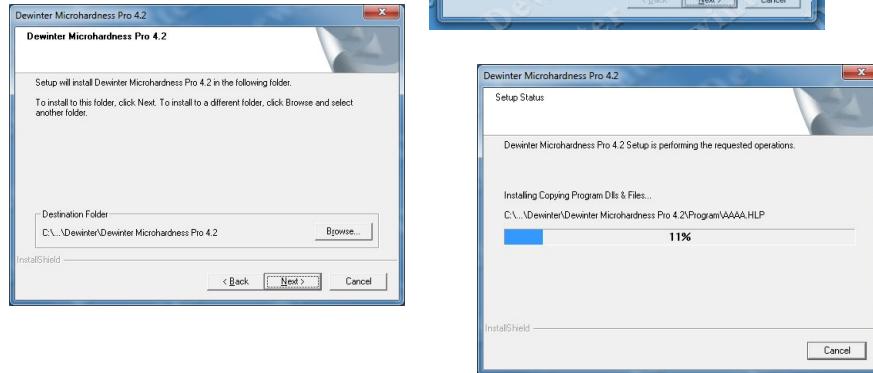



3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metallographic.com Web: http://www.metallographic.com

### 9.1.2.3 Load driver for USB dongle key

Click on Install Driver and then Exit when completed





### 9.1.2.4 Click on Finish and then Reboot the computer

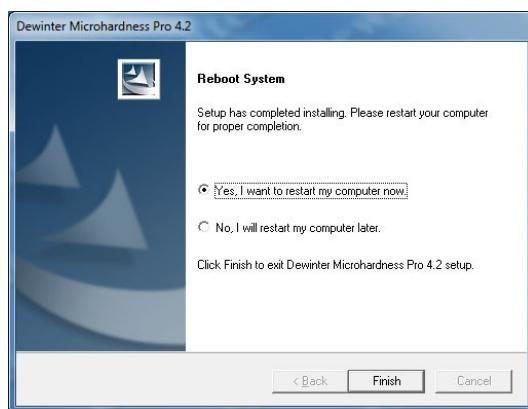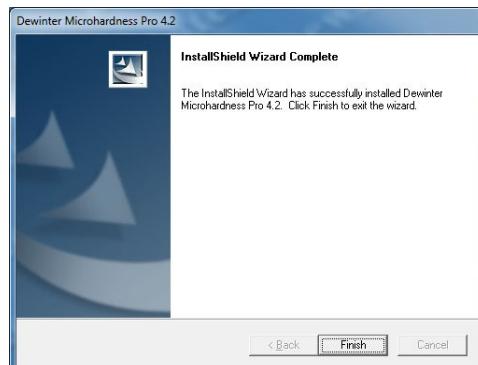



### 9.1.3 Camera set-up

The initial step after loading the software is to set-up the camera

Select the **Capture** drop down menu. Then click **Select Source** and choose from the available sources.

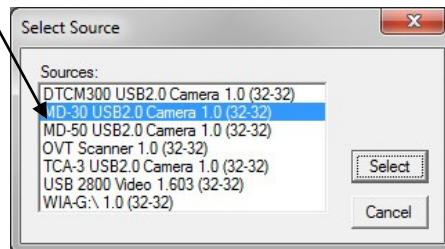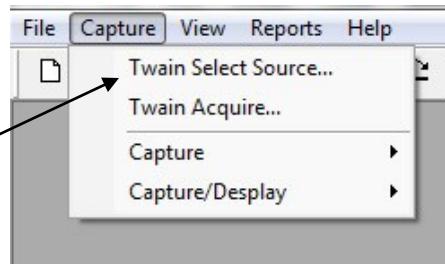


Select the **Capture** drop down menu. Then click **Select Twain Acquire**

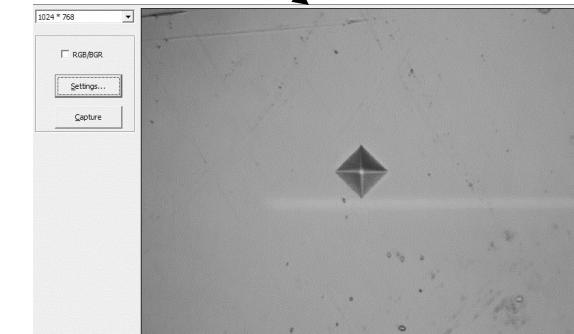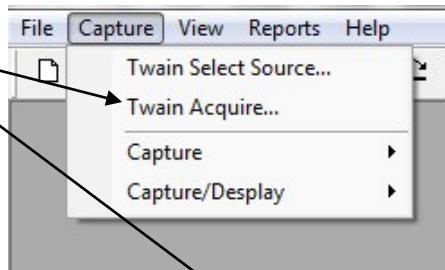

### 9.1.4 Calibration



After setting up the camera the next step is to calibrate the optics of the microhardness tester with the computer being used (note if the computer is changed or if any of the optics are changed in the future, this recalibration process needs to be repeated).

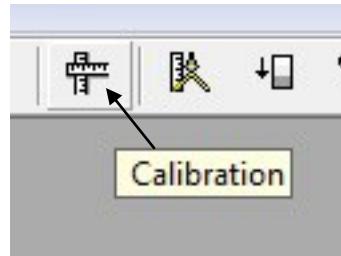

#### 9.1.4.1 Calibration slide

Using a calibration standard capture the image for the 40X objective and save the image in the x-direction and in the y-direction

The images shown are for a calibration slide having 0.1 mm increments (100 microns).

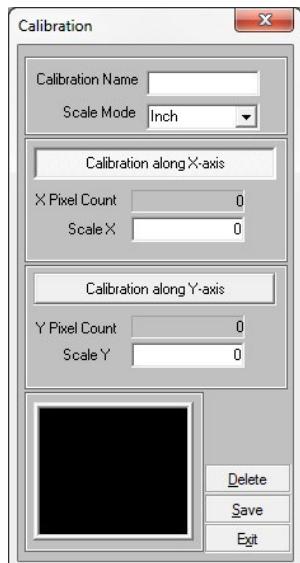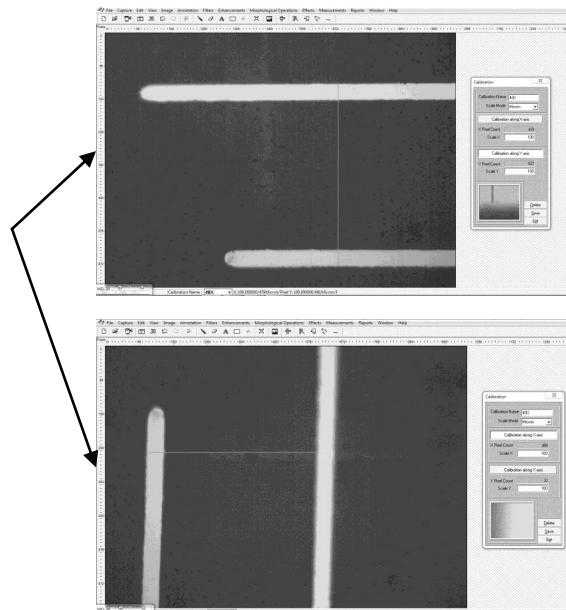

Enter the required information and click on Calibration along X or Y-axis.

Click and drag at the front edge of the starting scale line to the front end edge of the ending line.

Save the calibration (to delete an old calibration click on delete and choose the calibration to be deleted).


### 9.1.5 Measurement

Select hardness test under the measurement icon


#### 9.1.5.1 Diagonal Measurement

Click on D1 diagonal radio button and click and drag the mouse to perform the measurement.


Repeat for D2

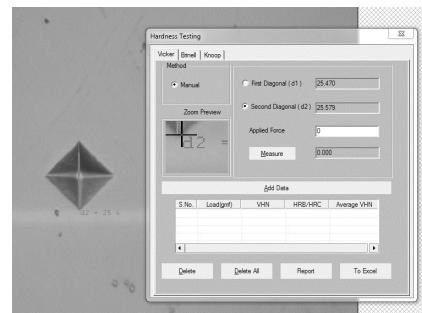


#### 9.1.5.2 Enter load (gmf)

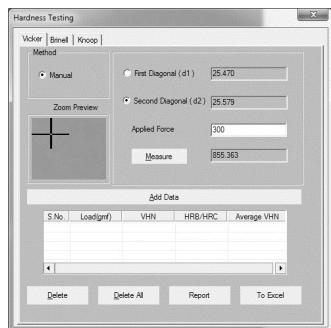

#### 9.1.5.3 Click on Measure


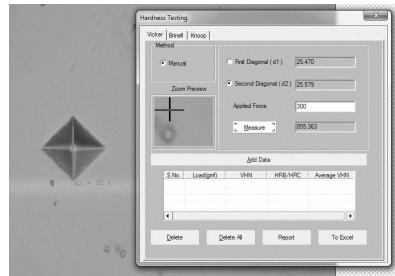


#### 9.1.5.4 Click on Add Date

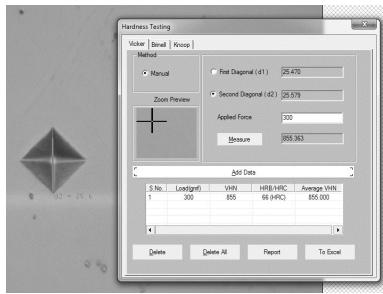

3601 E. 34th St. Tucson, AZ 85713 USA Tel. +1 520-882-6598 Fax +1 520-882-6599 email: pace@metalgraphic.com Web: http://www.metalgraphic.com

### 9.1.5.5 Click on Excel

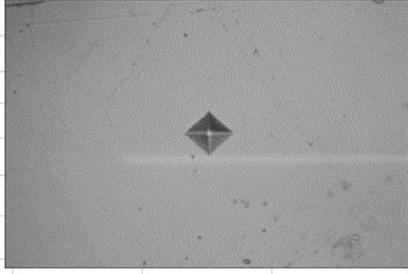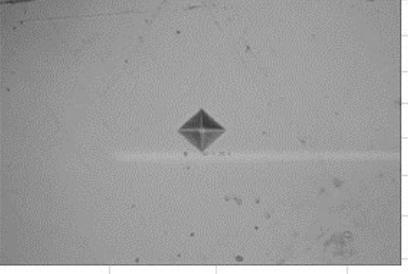MICROMEASUREMENT TEST REPORT								
7	NAME :				APPLICATION :			
8	EVALUATION DATE :				OPERATOR :			
9	SAMPLE INFO ID :				MICROSCOPE OBJ :			
Sr.No Load ( gmf ) VHN HRC/HRB Average VHN								
12	1	300	855	66 (HRC)	855			
13								
14								
15								
16								
17								
18								
19								
20								
21								
22								