



<b>Equipment Type:</b>	400 mm (~16-inch) Manual Abrasive Cut-off Machine
<b>Model:</b>	<b>MEGA-T400S</b>
<b>Electrical Requirements:</b>	480 Volts (3-phase)
<b>Frequency:</b>	50/60 Hz
<b>Motor Horsepower:</b>	Variable speed Servo Motor: 12 KW (16 hp) @480V, 3-phase
<b>Manual Revision Date:</b>	March 16, 2023

This instruction manual is provided with each piece of delivered equipment.

Please read this instruction manual carefully and follow all installation, operating and safety guidelines.

## WARRANTY

### Terms and Conditions applying to all PACE Technologies Products

#### 1. LIMITED WARRANTY AND DISCLAIMER:

PACE Technologies Products are warranted for two years from the purchase date to be free from defects in material and workmanship under correct use, normal operating conditions, and proper application. PACE Technologies obligation under this warranty shall be limited to the repair or exchange, at PACE Technologies option, of any PACE Technologies Product or part which proves to be defective as provided herein. PACE Technologies reserves the right to either inspect the product at Buyer's location or require it to be returned to the factory for inspection. Buyer is responsible for freight to and from factory on all warranty claims. The above warranty does not extend to goods damaged or subjected to accident, abuse or misuse after release from PACE Technologies warehouse, nor goods altered or repaired by anyone other than specifically authorized PACE Technologies representatives. PACE Technologies shall not in any way be responsible for the consequences of any alteration, modification or misuse unless previously approved in writing by an officer of PACE Technologies. **Note: Corrosion is considered a maintenance issue and not a warranty issue.**

PACE TECHNOLOGIES MAKES NO EXPRESS WARRANTIES OTHER THAN THOSE WHICH ARE SPECIFICALLY DESCRIBED HEREIN. Any description of the goods sold hereunder, including any reference to Buyer's specifications and any description in catalogs, circulars and other written material published by PACE Technologies, is the sole purpose of identifying such goods and shall not create an express warranty that the goods shall conform to such description.

THIS WARRANTY IS EXPRESSLY IN LIEU OF ALL OTHER WARRANTIES, EXPRESSED OR IMPLIED. THERE ARE NO IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR PARTICULAR PURPOSE. THIS WARRANTY STATES PACE TECHNOLOGIES ENTIRE AND EXCLUSIVE LIABILITY AND BUYER'S EXCLUSIVE REMEDY FOR ANY CLAIM FOR DAMAGES IN CONNECTIONS WITH PACE TECHNOLOGIES PRODUCTS. PACE TECHNOLOGIES WILL IN NO EVENT BE LIABLE FOR INCIDENTAL OR CONSEQUENTIAL DAMAGES WHATSOEVER, NOR FOR ANY SUM IN EXCESS OF THE PURCHASE PRICE.

#### 2. LIABILITY CAP:

PACE Technologies maximum aggregate liability for loss and damage arising under, resulting from or in connection with the supply or use of the Equipment and Consumables provided under this purchase, or from the performance or breach of any obligation (s) imposed hereunder, whether such liability arises from any one or more claims or actions for breach of contract, tort, (including negligence), delayed completion, warranty, indemnity, strict liability or otherwise, unless otherwise limited by the terms hereof, shall be limited to one hundred percent (100%) of the purchase price.

#### 3. DELIVERY:

Customer assumes and shall bear the risk of all loss or damage to the Products from every cause whatsoever, whether or not insured, and title to such Products shall pass to Customer

upon PACE Technologies delivery of the Products to the common carrier of Pace Technologies choice, or the carrier specified in writing by Customer, for shipment to Customer. Any claims for breakage, loss, delay, or damage shall be made to the carrier by the Customer and Pace Technologies will render customer reasonable assistance in prosecuting such claims.

**4. ACCEPTANCE:**

Customer shall inspect the Products promptly upon receipt of delivery. Unless customer objects in writing within thirty (30) business days thereafter, customer shall be deemed to have accepted the Products. All claims for damages, errors, or shortage in Products delivered shall be made by Customer in writing within such five (5) business day period. Failure to make any claim timely shall constitute acceptance of the Products.

**5. PAYMENT:**

Customer agrees to provide timely payment for the Products in accordance with the terms of payment set forth on the reverse side hereof or in any proposal submitted herewith. If any payment is not paid on or before its due date, Customer shall pay interest on such late payment from the due date until paid at the lesser of 12% per annum or the maximum rate allowed by law.

**6. DEFAULT:**

If Buyer is in default (including, but not limited to, the failure by Buyer to pay all amounts due and payable to Seller) under the work or purchase order or any other agreement between Buyer and Seller, Buyer's rights under the warranty shall be suspended during any period of such default and the original warranty period will not be extended beyond its original expiration date despite such suspension of warranty rights.

**7. MISCELLANEOUS PROVISIONS:**

This agreement has been made in and shall be governed by the laws of the State of Arizona. All disputes arising under or relating to the purchase of the equipment shall be brought and resolved solely and exclusively in the State of Arizona, Pima County. These terms and conditions and the description of the Products on the reverse side hereof or in any proposal submitted herewith constitute the entire agreement and understanding of the parties with respect to this sale and supersede all prior and contemporaneous agreements or understandings, inducements or representations, expressed or implied, written or oral, between the parties with respect hereto. Any term or provision of this Agreement may be amended, and any observance of any term of this Agreement may be waived, only by a writing signed by the party to be bound. The waiver by a party of any breach shall not be deemed to constitute a waiver of any other breach. Should suit be brought on this Agreement, the prevailing party shall be entitled to recover its reasonable attorneys' fees and other costs of suit including costs and attorneys' fees incurred on appeal or in collection of any judgment., errors, or shortage in Products delivered shall be made by Customer in writing within such five (5) business day period. Failure to make any claim timely shall constitute acceptance of the Products.

**8. RESTOCKING FEE:**

All Returns are subject to a restocking charge equal to 15% (fifteen percent) of the Invoice, unless the Goods are proved to be non-conformed by PACE Technologies.



**Contents**

	PAGE
<b>Nameplate</b>	i
<b>Warranty</b>	ii
<b>1.0 Product Description</b>	1
<b>2.0 Shipping, Unpacking and Installation</b>	6
<b>3.0 Safety Guidelines</b>	9
<b>4.0 Start-Up and Operation</b>	12
<b>5.0 Maintenance</b>	18
<b>6.0 Trouble Shooting</b>	19
<b>7.0 Cutting Basics</b>	21
<b>8.0 Electrical and Mechanical Drawings</b>	24

## 1.0 Product Description

The **MEGA-T400S** is a manual wheel feed or table feed abrasive cut-off machine for cutting materials ranging from soft aluminum metals to hardened tool steels. It is ideal for the metallographic laboratory, as well as for small industrial or production applications.

The **MEGA-T400S** is very robust and durable with its cast aluminum alloy and stainless steel construction. Featuring a corrosion-free T-slot table, the **MEGA-T400S** is a very versatile wheel/table feed floor model metallographic cutter.


### Powerful motor

The **MEGA-T400S** is equipped with a powerful belt driven servo motor. In addition, the motor has an inductive brake for a quicker response of halting the blade (approx. 2-3 seconds).

### Manual feed and table feed cutting

The **MEGA-T400S** is an easy-to-use manual cutter with simple controls. Cutting can be accomplished with either a wheel feed or table feed.


The **MEGA-T400S** features a removable side port panel for sectioning longer samples.

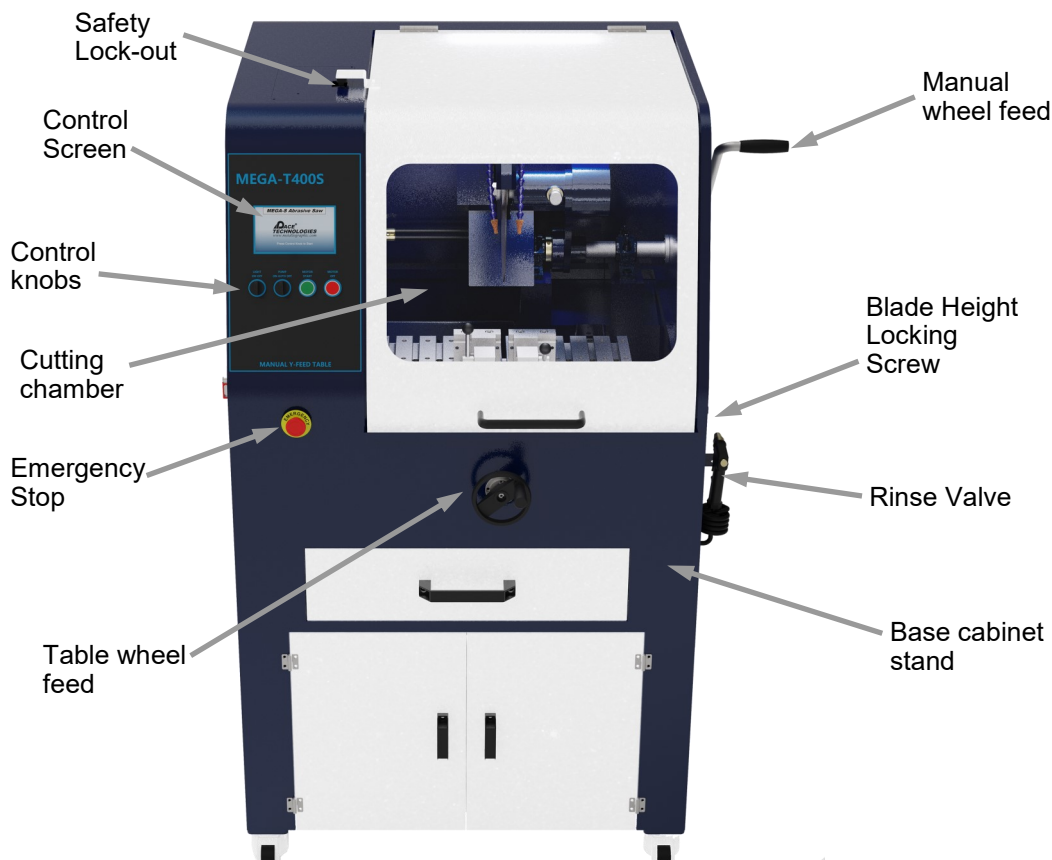




Coolant applied through the shroud directly to the blade.

Flexible coolant hose for application to cutting area (optional).

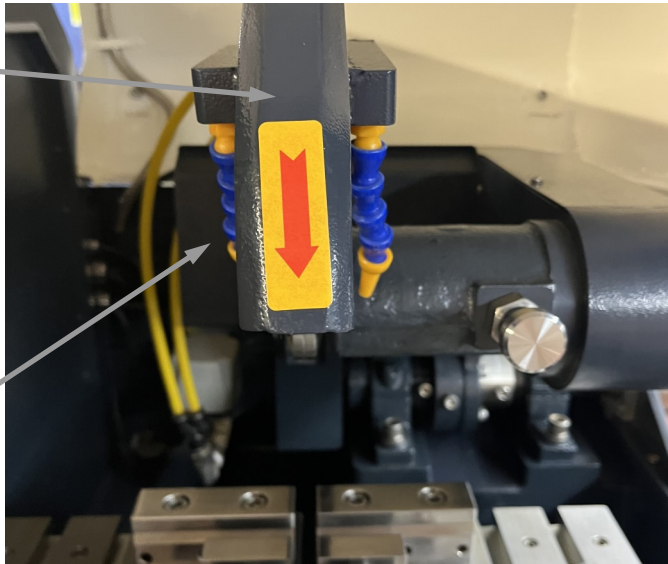

Easy to use vises.

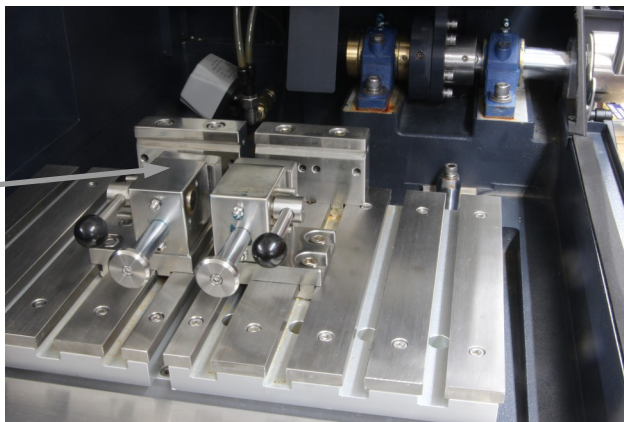




## 1.2 Technical Specifications


<b>Electrical specifications:</b>	480V, 3-phase
<b>Motor power:</b>	Constant torque variable speed Servo Motor: -14.4 KW (19.2 hp) 480V
<b>Speed (rpm):</b>	Variable 1000-3000 rpm
<b>Adjustable cutting load:</b>	2-15 amps (max. 20 amps)
<b>Electronic motor brake:</b>	Stops in ~3-5 seconds
<b>Cut-off wheels:</b>	16 inch mm (~400 mm) diameter
<b>Wheel arbor:</b>	32 mm (~1.25 inch) diameter
<b>Maximum sample diameter:</b>	6-inch (150 mm) diameter 3 x 6-inches (75 x 150 mm) (HxD)
<b>Wheel feed:</b>	Manual
<b>Table feed:</b>	Manual
<b>Vertical movement:</b>	5 inch (125 mm)
<b>Weight:</b>	Approx. 880 lbs (400 kg)
<b>Dimensions (WxDxH):</b>	<b>Hood Closed:</b> 48 x 51 x 63 inch (1220 cm x 1300 cm x 1600 cm) <b>Hood Open:</b> 48 x 51 x 66 inch (1220 cm x 1300 cm x 1665 cm)
<b>Table dimensions (WxD):</b>	11 x 18.25 inch t-slot (280 x 470 mm)
<b>Cutting chamber:</b>	Cast aluminum block construction
<b>Hood and body:</b>	Fabricated steel with shock resistant plastic window
<b>Working temperature:</b>	32° - 100°F (0 - 40°C)
<b>Shipping temperature:</b>	32° - 100°F (0 - 40°C)
<b>Storage temperature:</b>	32° - 100°F (0 - 40°C)
<b>Recirculation system (included):</b>	14 gallons (53 liters)

<b>EU Directives:</b>	Machinery directive 2006/42/EC RoHS Directive 2011/65/EU
-----------------------	---

## 1.3 Mechanical Schematic

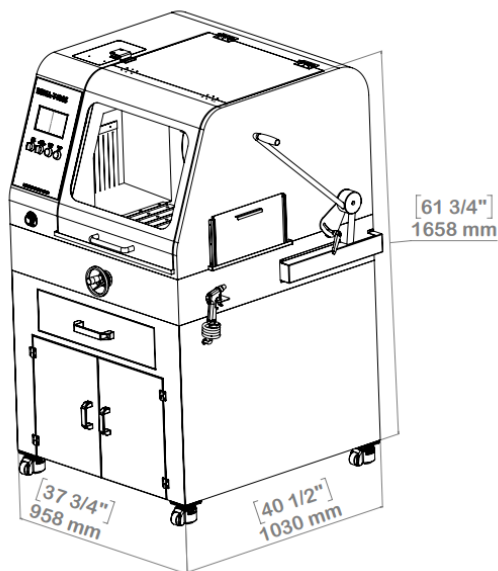

**MEGA-T400S** hood closed dimensions

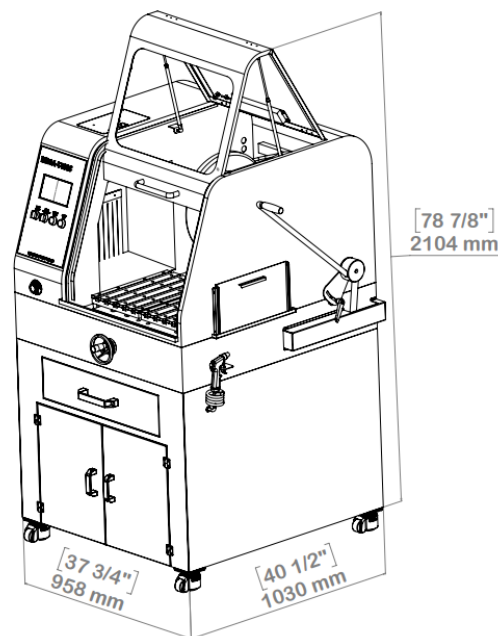

**MEGA-T400S** hood open dimensions

**Note:** Installation of the **MEGA-T400S** should be on a flat sturdy surface, with easy access to electrical connections.

## 2.0 Unpacking, Shipping and Installation

### 2.1 Unpacking

Unit is delivered in a crate. Unpack and check for completeness of parts.

**Measures WxHxD:** Approx. 51x51x75 inches

**Weight:** Approx. 1000 pounds


### 2.2 Shipping

When moving crate, lift from bottom. To move **MEGA-T400S** use support bars provided with the machine.

The **MEGA-T400S** contains sensitive electronic and mechanical components. **Do not drop.**




**! CAUTION:** Heavy equipment. Take care to avoid bodily injury.

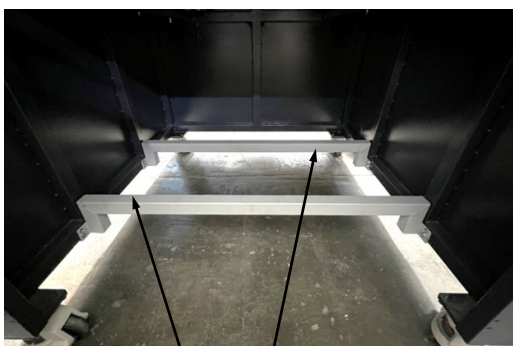

Machine includes support bars to lift machine and move with the aid of a forklift. The casters are adjusted so the saw can be rolled into place easier. After fixing the position, retract the wheels on the caster.

Once the cutter is in place, remove the forklift support bars so the tank will move inside the cabinet.

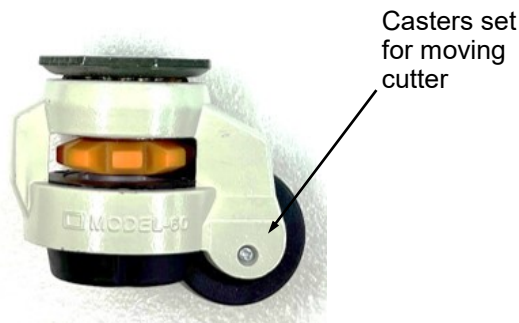

Casters set for stationary support



(Installation continued on next page)


## 2.3 Installation

**!** Install unit carefully! Improper installation voids warranty.

The **MEGA-T400S** should be placed on a flat stable surface.

**!** Connect the electrical connections to a 3-phase power line. The motor driver max amp draw is 20 amps. However, the operating power is set to 15 amps.

Not  
used

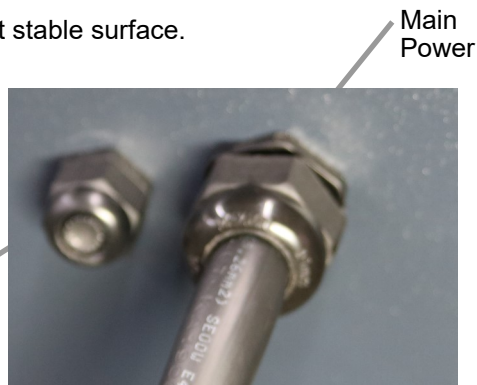

**!** Electrical Connections - The standard saw is designed to operate at 208-220V 3-phase (120V-120V-120V-ground) (optional 380V)

Verify voltage on nameplate



Main  
Power  
switch

**!** Connect pump connector, supply line from cutter to the pump and insert the drain hose into the basket filter

Coolant  
supply  
line to  
cutter




Electrical  
connector  
for pump



**!** Verify the direction of rotation of the cut-off wheel. The wheel should turn from top to bottom as viewed from the front of the machine. If not reverse any two of the 3-phase power wires.

(Installation continued on next page)

## Recirculating tank connections

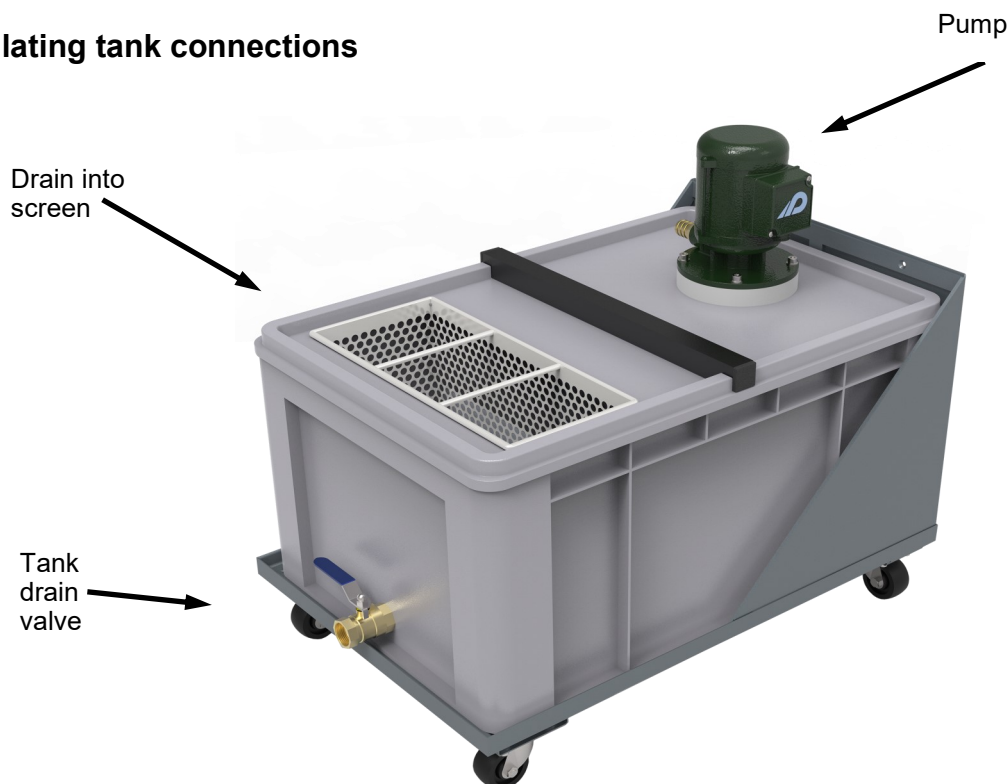

**External coolant supply:** Attach 3/4 inch tube between pump and cutter.

**Drain:** 1.5-inch (38 mm) tube.

**Electrical connections:** Connect electrical power cable to source.

**NOTE:** Inspect the operating voltage on the name plate.

**Electrical connection for external coolant supply:** Power for recirculation system comes from the **MEGA-T400S**,

**Recommended Coolant:** To help minimize corrosion apply a coolant with a anti-corrosion additive or an oil emulsion coolant is recommended. **Note:** Recommended that the coolant be changed every 50 cuts or once a week, whichever occurs first.

(Installation continued on next page)

### **3.0 Safety Guidelines**

#### **3.1 Warning Sign**

**!** This sign points to special safety features on the machine.

#### **3.2 Safety Precautions**

Careful attention to this instruction manual and the recommended safety guidelines is essential for the safe operation of the **MEGA-T400S**.

Proper operator training is required for the safe operation of the **MEGA-T400S**. Any unauthorized mechanical and electrical change, as well as improper operation, voids all warranty claims. All service issues need to be reported to the manufacturer / supplier.

**!** Before operating, cutting chamber hood must be closed. After cutting, the safety latch will not open for approximately 2-8 seconds after pressing the emergency or red stop (red) button.

**!** Use only certified cut-off wheels from a professional supplier. Improper blades selection voids warranty. (For appropriate blade selection, refer to the Abrasive Blade Selection Guidelines Chart in Section 4.4)

**!** Disconnect power before opening the main unit.

**!** Replacement parts should be installed only by qualified personnel.

**!** Securely clamp the part / sample to the working table. During cutting, consider that the part may pinch and cause jamming of the cut-off wheel. Use the appropriate clamping devices to avoid this occurrence.

**!** Never start a cut under load.

**!** Make sure that the cut-off blade is rotating down and into the sample.


#### **3.3 Emergency Statement**

The **MEGA-T400S** abrasive cutter has been designed for cutting metallographic specimens. DO NOT CUT oversize samples. Always follow proper operational guidelines and avoid contact with moving parts, lubricants and abrasives. Seek appropriate medical care for cutting injuries.

## Recirculation Tank

The **MEGA-T400S** is equipped with a coarse filter to remove large particulate and broken blades

**NOTE:** Dispose of the old coolant in accordance with federal, state and local regulations.

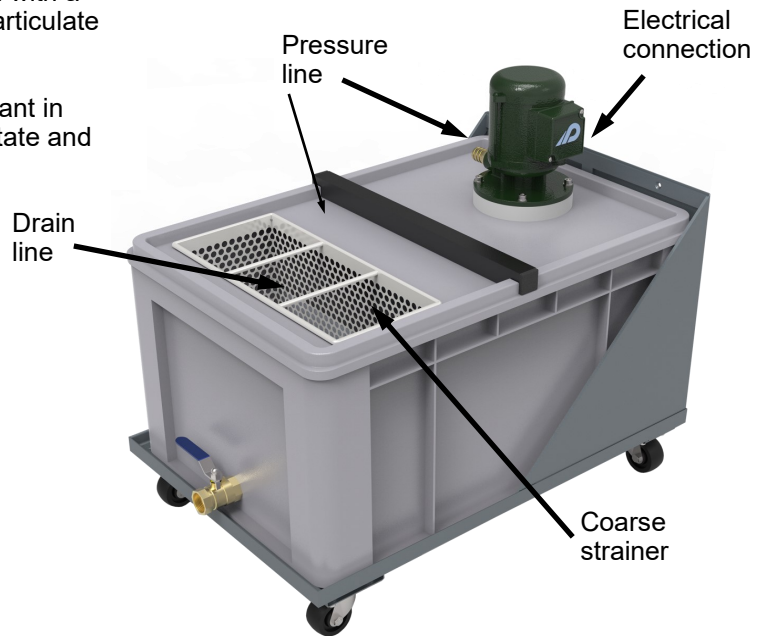

## Rinse line

After use, water off the cutting chamber to prevent debris and corrosion build-up.

**NOTE:** Leave hood open to dry after cleaning.





## 3.4 Safety Tests

**!** Examine and verify that the **MEGA-T400S** safety devices and operating performance are in good working condition before use. The following safety checks are considered significant:



Emergency stop switch

### Emergency stop switch

<b>Test:</b>	Activate the main switch and the close hood; depress emergency stop switch.
<b>Proper Response:</b>	Machine powers down and motor stops in 5-10 seconds
<b>Malfunction:</b>	Machine does not lose power.
<b>Corrective measure:</b>	If system does not power down, disconnect power supply cord and call service technician.



Magnetic Safety switch

### Magnetic safety lock out switch

<b>Test:</b>	Activate the main switch and close the hood; turn cut-off wheel ON then OFF; Try to open the hood.
<b>Proper Response:</b>	The hood does not open for approximately 3-5 seconds.
<b>Malfunction:</b>	If hood opens immediately.
<b>Corrective measure:</b>	If the cut-off wheel does not power down, disconnect the power supply cord and call a service technician.

## HOOD (Closing / Opening)

1. When closing the hood, press down on the hood to ensure that the safety switch is engaged.
2. If the hood lock does not unlock **DO NOT FORCE OPEN**, or try pushing down on the hood to release the safety pin.




## 4.0 Start-up and Operation

### 4.1 General

The **MEGA-T400S** is a manual wheel and manual table feed cutter.

Note: When first powering up the **MEGA-T400S** the initialization time for the saw is approximately 15 seconds.  
 The servo motor is also designed to increase speed slowly and typical takes 4-5 seconds to reach full speed.

### 4.2 Control Panel

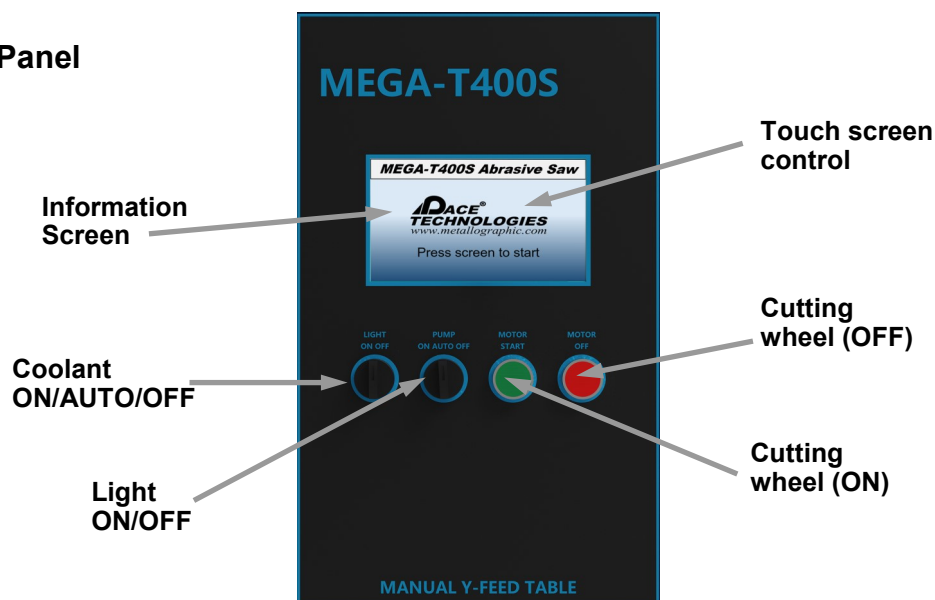

<b>Cutting wheel ON:</b>	Starts the cutting wheel.
<b>Cutting wheel OFF:</b>	Stops the cutting wheel (inductive brake engages to slowing wheel). The hood will not open until safety lock releases in approximately 5-8 seconds after pressing the stop button.
<b>Coolant ON/AUTO/OFF:</b>	Operates the coolant pump in Auto or Manual mode. <b>Must be selected before starting the motor.</b>
<b>Light ON/OFF:</b>	Turns inside chamber light On or Off.


**! Discontinued use if safety lock fail**

## START UP SCREEN



**Start up screen** - Press control panel knob to progress to information screen


## INFORMATION READY SCREEN

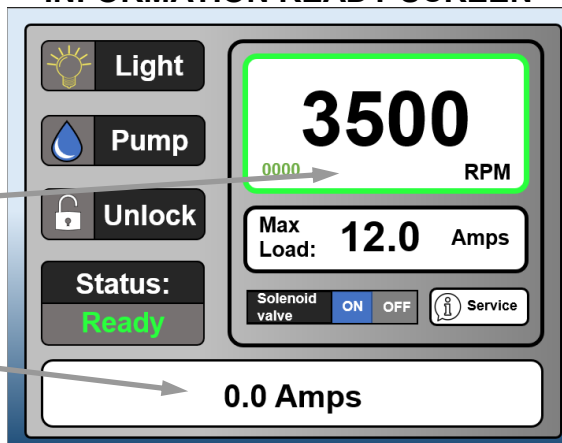

Set Speed

Cutting force (amps) Setting

## OPERATING SCREEN

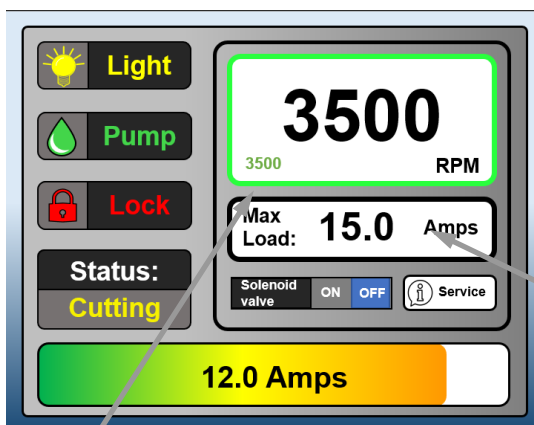


Actual Operating Speed

Percentage or fraction of amps to set amount

**Information screen** - set speed and force (amps) by selecting parameter on the touch screen and adjusting the parameter in the sub-screen. Repeat for other parameter.

## 4.3 Wheel and Table Feed

### 4.3.1 Wheel feed cutting

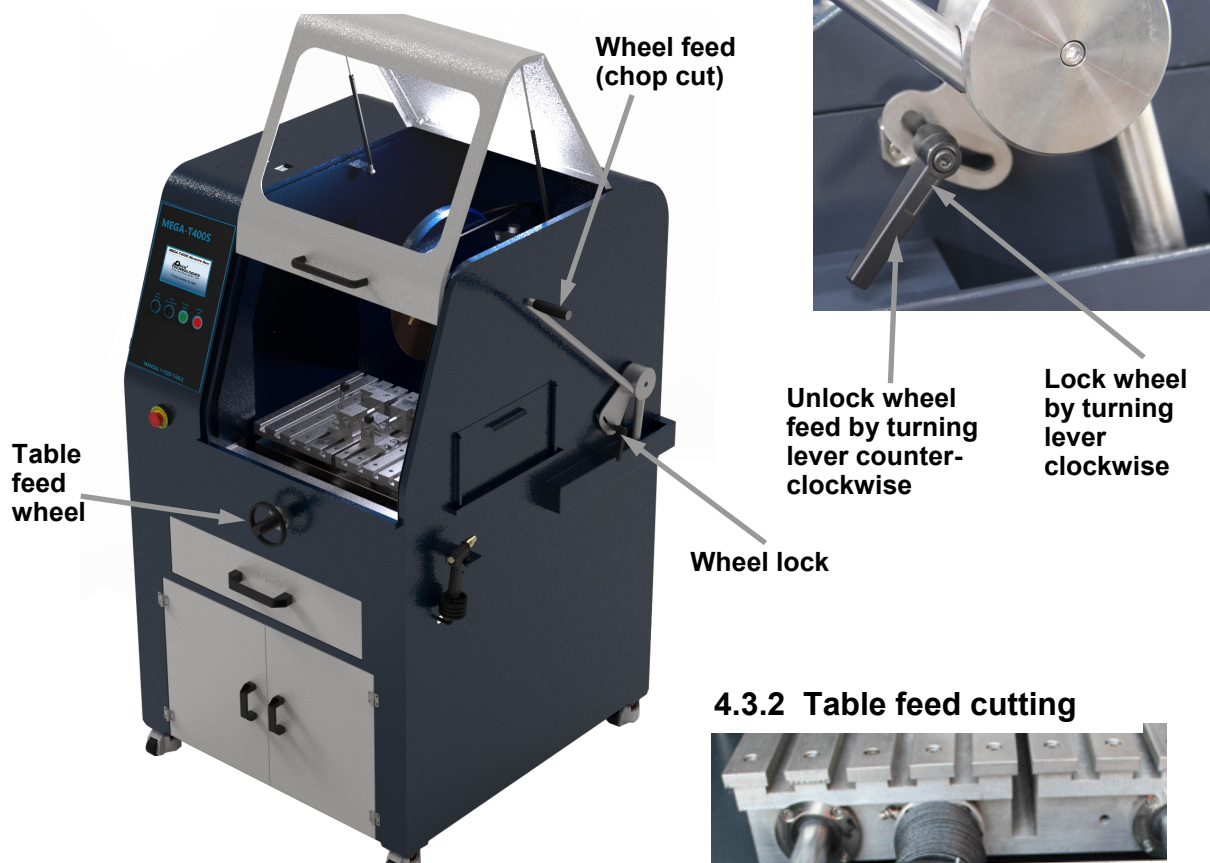

### 4.3.2 Table feed cutting



## 4.4 Changing abrasive cut-off wheels

1. Remove blade locking bolt (reverse threaded).
2. Position new 350 mm abrasive cut-off wheel into position.


**! Use only certified abrasive cut-off wheels.**

4. Gently tighten blade locking bolt (Note: locking bolt is reverse threaded).

### Changing Blade

**Remove:** Loosen by turning clockwise

**Tighten:** Turn counter-clockwise

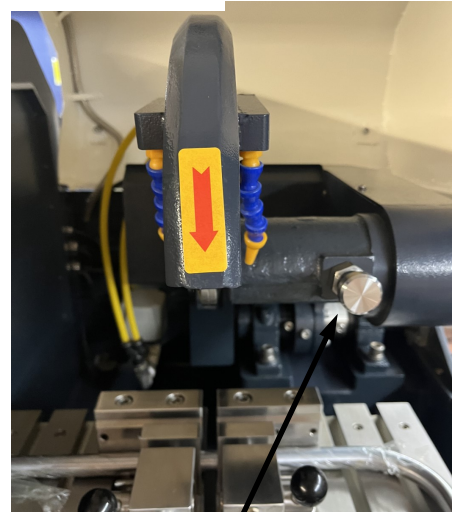

**To prevent shaft from turning, push in the locking rod to hold in**

## ABRASIVE Cutting CONSUMABLES SELECTION GUIDELINES CHART

### Abrasive Cutting Fluids

<a href="#">Pace Product Name</a>	<a href="#">Catalog Number</a>	<a href="#">Packaging</a>
MAXCUT™ Abrasive Cutting Fluid (32 oz)	MAXCUT 1000-32	32 oz
MAXCUT™ Abrasive Cutting Fluid (128 oz)	MAXCUT 1000-128	1 gallon
MAXCUT™ 2 Cutting Fluid with Anticorrosion Additive (32 oz)	MAXCUT2 1000-32	32 oz
MAXCUT™ 2 Cutting Fluid with Anticorrosion Additive (128 oz)	MAXCUT2 1000-128	1 gallon

### 350 mm (~14/16-inch) MAXCUT™ Abrasive Blades 32mm (~1¼-inch) Arbor

<a href="#">Pace Product Name</a>	<a href="#">Catalog Number</a>	<a href="#">Packaging</a>
Soft Non-ferrous Materials (Aluminum, Brass, Zinc, Etc.)	MAX-E350/400	10
Hard Non-ferrous Materials	MAX-C350/400	10
Soft Steels	MAX-E350/400	10
Case Hardened Steels	MAX-VHS350/400	10
Steels and Ferrous Metals	MAX-D350/400	10
Universal Thin Resin/rubber Alumina Blade	MAX-A350/400	10



## 4.5 Manual cutting (step-by-step procedure)

1. Connect electrical power and turn on the master power switch on the side panel before opening the hood.
2. Position fixtures and sample near to cut-off wheel.
3. Close the hood.
4. Turn on the coolant and cutter motor.
5. Hand feed the cutting wheel or table feed into the sample.
6. Turn off the cutter and allow the blade to completely stop before opening cover. Note magnetic safety switch has an approximately 3-5 second delay. The cutter hood cannot be opened until the safety switch is deactivated.

## 4.6 Fixturing sample

For proper clamping, use the appropriate clamping vises to securely hold the sample in place. It is recommended that both sides of the part be clamped to avoid pinching of the blade (possibly breaking the blade) and to minimize burning of the workpiece during cutting.

**!** For proper fixturing, take into account the initial stress on the samples.

### Fixture examples



#### Fast clamping


Left lever (Cat. No. QCL-1000)  
 Right lever (Cat. No. QCR-1000)

**Note:** vises have a slot in them so that the face does not rotate




#### UNIVERSAL Clamping vise

UNIVERSAL (Cat. No. MG-01)  
 Clamping height 45mm (1.8 inch)

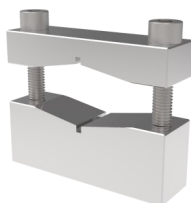

**MEGA Saddle  
Clamping vise  
(Cat No. MEGA  
702)**

**MEGA Fastener  
Vise (Cat No. MEGA-711)**

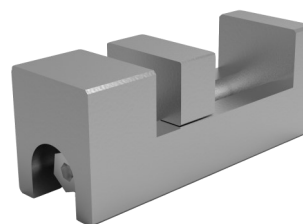



## 5.0 Maintenance

### 5.1 Introduction

The recommended maintenance for the **MEGA-T400S** includes the following:

1. **Important: Clean the table feed shaft after using the cutter to keep the abrasive from building up in the bushing. If the table feed begins to become more difficult use WD-40 to clean the feed shaft until the feed improves.**
2. Use an anti-corrosion cutting fluid such as the **MAXCUT 2** cutting fluid. It is recommended that the cutting fluid be changed every 50 cuts or weekly for maximum effectiveness.
3. Rinse and dry the table and vises after use.
4. When the cutter is not in use, OPEN the HOOD to prevent humidity and corrosion build-up inside the cutting chamber.
5. If cutter will not be used for an extended period, coat the table and vises with a water displacement oil such as WD-40.
6. Use only recommended abrasive blades to minimize overheating of the bearings and motor. Do not exceed the recommended amp rating of the motor.
7. Grease the appropriate fittings periodically (minimally after every 200 cuts or every other month, whichever occurs first).
8. Clean and apply a film of WD-40 on the table feed screw and bearing guides.
9. Periodically spray WD-40 between the inner arbor flange and the bearing housing to keep debris from building up.
10. Immediately replace any broken or damaged parts.

**Note: corrosion is considered a maintenance issue and not a warranty issue**

### 5.2 Cleaning outside cabinet

The cabinet and front shield should be cleaned occasionally with a moistened cloth. Do not use any chemicals or cleaning abrasives. We recommend an anti-condensation gel be used on the front panel.

## 6.0 Trouble Shooting

Problem	Cause	Solution
<b>No power or function</b>	a. Unit is disconnected from main electrical power supply. b. Main power switch is off. c. Emergency stop button engaged.	a. Verify electrical source and connection. b. Turn on main power switch. c. Release by turning clockwise.
<b>Main motor does not operate</b>	a. Hood not closed. b. Overload relays activated	a. Close hood b. Restart after resetting relay.
<b>Pump motor does not operate</b>	a. Plug is disconnected	a. Check plug.
<b>Excessive vibration during cutting</b>	a. Not cutting with enough force b. Specimen locked too far back on cutting table c. Incorrect blade selection	a. Increase cutting force. b. Move specimen forward. c. Use correct blade.
<b>Corrosion of table of vises</b>	a. Improper cleaning and drying b. Incorrect or old cutting fluid.  c. Cutting chamber has become a humidity / corrosion chamber d. Long term storage oxidation	a. Rinse and dry table and vises. b. Use an anti-corrosion cutting fluid such as the <b>MAXCUT 2</b> . It is recommended that the cutting fluid be changed every 50 cuts or weekly for maximum effectiveness. c. Open hood when not being used  d. If cutter will not be used for an extended period of time, coat the table and vises with a water displacement oil such as WD-40.
<b>ERROR 01*</b>	<b>Sudden stopping of motor after red light flashes</b> a. Feed rate too high b. Amperage draw is too high c. Improper blade selection  <b>Sudden stopping of motor with no red light warning</b> a. Feed rate too high b. Amperage draw is too high c. Improper blade	<b>Turn off main power for several seconds to reset and correct problem:</b> a. Reduce feed rate b. Decrease cutting force c. Use correct blade.  <b>Turn off main power for several seconds to reset and correct problem:</b> a. Reduce feed rate b. Decrease cutting force c. Use correct blade.
<b>Hood does not open when motor is off</b>	a. Trying to open the hood before the safety switch releases (approximately 2-3 seconds) b. Safety switch still does not disengage	a. Wait until the safety lock out clicks b. Try pushing down on the hood to release the safety pin in the safety lockout switch

\* If the motor stalls and an ERROR 01 is displayed, power off the cutter and wait until you hear a small click from the driver for it to reset, typically <15 seconds.



## 6.1 Trouble Shooting Error Codes

Error Code	Code Description	Solution
ERR 01	IPM error	Disassemble and return to supplier
ERR 02	Over current	Check the setting parameters
ERR 03	EEPROM error	Rewrite the default value
ERR 04	Current sampling	Rewrite parameters
ERR 05 (GF)	Bus under or over voltage	Check the voltage and parameters
ERR 06	Communication anomaly	Check the wire connection
ERR 07	Locate failed	Check the encoder signal; Or the current feedback wrong. Reposition and try again.
ERR 08	Lost step when running	Normally, the reverse angle is 2 degrees, and the reverse angle is 20 degrees when positioning will be wrong. Reposition and try again.
ERR 09	Overload	Cutting too fast, too much load, or reduce the speed.
ERR 10	Module temperature is too high	Check whether the fan works properly
ERR 12	Power failure	Disassemble and return to supplier
ERR 13	Encode failure	Check wire connection, or the encoder was damaged
ERR 14	Motor get stuck	Check the motor or cutting speed is too high

## 7.0 Metallographic Cutting Basics

### Philosophy:

1. Cutting is the most important metallographic sample preparation step
2. Minimize damage during sectioning
3. Less sectioning damage = reduced number of grinding and polishing steps


**Note:** In some cases cutting damage cannot be removed by subsequent polishing operations (especially for brittle samples)

### Time Well Spent:

1. Spend more time on sectioning in order to decrease overall sample preparation time
2. Protect integrity of the specimen with proper sectioning

#### Wheel feed


- High speed (3500 rpm)
- High load (10 amps)

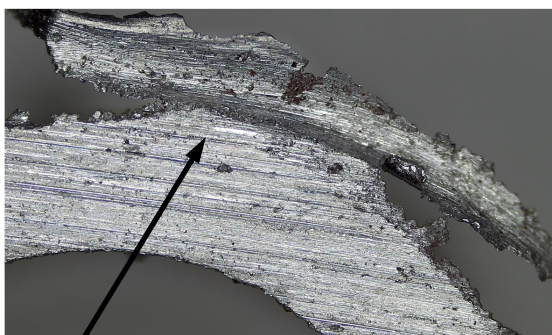

Coarse finish with  
large cutting burr

#### Table feed

- Low speed (1500 rpm)
- Low load (3 amps)

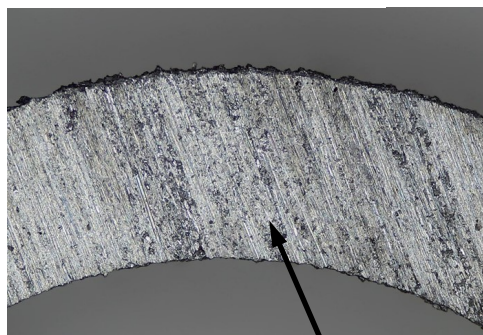


Finer finish with  
minimal cutting burr

## 7.1 Metallographic Cutting Pentagon

Five variables for successful sectioning:

- **Abrasive Type**
  - Silicon Carbide
  - Alumina
- **Abrasive Bond**
  - Resin
  - Rubber
- **Cutting Force**
  - Control by motor current
  - Cutter Speed
    - Low (1500-2000 rpm)
    - Medium (2000-2500 rpm)
    - High (2500-3500 rpm)
- **Coolant**
  - Water-based
  - Oil Emulsion-based

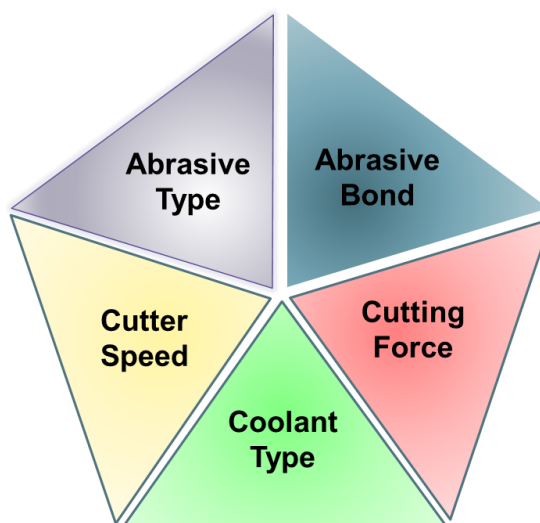

**Metallographic Abrasive Blades:**

- **Resin Bonded Blades**
  - Breakdown easily
  - Less smell and burning
  - More versatile
  - Alumina abrasives
- **Rubber Bonded Blades**
  - Longer life
  - Produces a burnt rubber smell
  - More likely to burn the sample
  - More difficult to match the blade to the material being cut
  - SiC and Alumina abrasives



## Machine Parameters:

- **Wheel Speed**
  - Slower speed for softer metals (1000-1500 rpm)
  - Medium speed for most metals (1500-2000 rpm)
  - Higher speed for harder samples (2500 rpm)
- **Cutting Force**
  - Controlled by limiting current to the motor
  - Softer metals lower load (3-5 amps)
  - Most metals medium load (6-10 amps)

## Recommended Cutting Procedures:

Material	Abrasive / Bond	Blade	Cutting Speed	Cutting Force
Soft non-ferrous metals (aluminum, brass, zinc, etc.)	Alumina / resin	MAX-E	1500 rpm	3-5 amps
Hard non-ferrous metals (titanium, zirconium, etc.)	SiC / rubber	MAX-C	2300 rpm	5-10 amps
Soft steels	Alumina / resin	MAX-E	1500 rpm	5-10 amps
Hard and case hardened steels	Alumina / resin	MAX-VHS	1900 rpm	10-15 amps
General purpose steel and ferrous metals	Alumina / resin	MAX-D	2500 rpm	10-15 amps
Universal thin blade	Alumina / rubber	MAX-A	2500 rpm	10-15 amps
Industrial longer life blades	Alumina / resin	MAX-I	3000 rpm	10-15 amps
Ceramics and glass	Diamond / metal	DMAX	3000 rpm	10-15 amps

The diagram illustrates the electrical control system for a machine, featuring a 3-phase 380V AC input and a 24V DC control circuit. The main power is supplied by a 3x380V LC1D25B7C circuit breaker. The power is then distributed to a Motor Driver and a Pump. The control circuit is powered by a 24V 0V transformer (TC2) and a 3A fuse. The control circuit includes a stop button (E-STOP) and a 24V DC power supply (VG144-A242) that provides DC24V 0V and DC24V 0V. The DC24V 0V line is connected to a 3x380V LC1D25B7C circuit breaker. The DC24V 0V line is also connected to a Motor Driver and a Pump. The Motor Driver is connected to a Motor. The Pump is connected to a 3x380V LC1D25B7C circuit breaker. The diagram includes a title block with fields for TITLE, DWG. No., DATE, NAME, and a revision table.

REV	DATE	NAME	DESIGN
1			

TITLE: \_\_\_\_\_  
 DWG. No.: \_\_\_\_\_

